

# COVID-19 Test Home Collection Kit



## Let's get started.

[WWW.PIXEL.LABCORP.COM](http://WWW.PIXEL.LABCORP.COM)

### Kit includes:




1 Saline  
collection tube



1 foam swab



1 Biohazard  
specimen bag

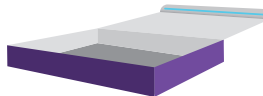

1 Pixel box



Registration  
Card



FedEx return pack  
UN 3373

**Intended use:** to collect a nasal swab for purpose of in vitro diagnostic testing in a CLIA laboratory.

IVD

Distributed by Pixel by LabCorp  
531 South Spring Street  
Burlington, NC 27215

This product has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA; This product has been authorized only for the detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens; and

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug, and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

PXC12820

# Self-Collection kit instructions

For In Vitro Diagnostic Use



1

**Please read all of the instructions before beginning collection. You must be 18+ to use this test kit.**

REGISTER ONLINE

2

**When you are ready to collect your sample, register your kit online at [www.pixel.labcorp.com/register](http://www.pixel.labcorp.com/register) and type in the 12-digit registration code that is located on your collection tube, Pixel box, and registration card.**

3




**Visit <https://bit.ly/fedexdropbox> to view FedEx drop box locations and pickup schedules. It's important to bring your sample to a drop box on the same day you collect it before the last Express pickup. Or refer to the flyer in your kit to schedule a FedEx pickup on the day of collection.**

4




**Wash and dry hands before opening the kit. Open your kit and place all the contents on a clean, dry surface.**

5

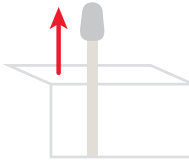


**Take the foam swab out of its package. Do not touch the tip of the swab with your hands.**

6

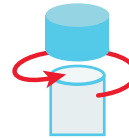


**Unscrew the top of the collection tube. Hold swab in one hand and collection tube in the other being careful not to spill the liquid. Do not drink the liquid.**

7

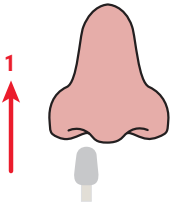


**Insert the tip of the foam swab into one nostril. The foam swab does not need to be inserted far – insert just until the foam tip of the swab is no longer visible. Rotate the swab in a circle around the entire inside edge of your nostril at least 3 times.**

8

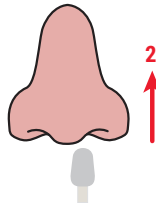

**Take the foam swab out of your nostril. Using the same end of the foam swab, repeat step 7 in your other nostril.**

9

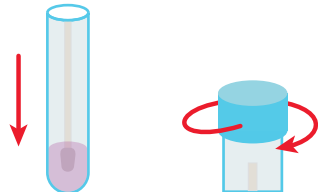


**Remove the foam swab from your second nostril and place in the collection tube. The end of the foam swab that went into your nose should be placed into the tube first so that it sits down in the liquid. Screw the top of the collection tube back on.**

10



**Wash and dry hands thoroughly again. Insert collection tube into the biohazard specimen bag. Seal the biohazard specimen bag by closing the zip lock seal.**

11

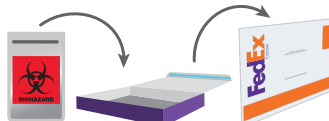

**Place the specimen bag back into your purple Pixel shipping box and close the lid. Place the shipping box into the FedEx return pack. Remove the adhesive cover strip and seal the return pack closed.**

12



**Deliver the postage paid, pre-addressed FedEx return pack to a FedEx drop box. It's important to bring your sample to a drop box on the same day you collect it before the last Express pickup. Or refer to the flyer in your kit to schedule a FedEx pickup on the day of collection.**