Phosphorus

COVID-19 RT-qPCR Test: at-home Ordering Guide

At-Home Collection Guide

For Prescription, In Vitro Diagnostic and Emergency Use Authorization ONLY

Ordering Note

A COVID-19 RT-qPCR Test order requires you to complete a detailed COVID-19 Medical Questionnaire that is reviewed by an independent provider. A portion of your fee goes to this independent physician for reviewing your test request.

Please follow collection instructions carefully to ensure an adequate specimen is collected for testing.

This product has not been FDA cleared or approved but has been authorized by FDA under an EUA for use by authorized laboratories. This product has been authorized only for the detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens; and, this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

1

Please visit our website ([COVID-19 RT-qPCR Test](#)) to request and pay for your test and related service fees. You must be 18 years or older to use this product.

2

You will be emailed a link to an online COVID-19 Medical Questionnaire to complete. A healthcare provider will authenticate the information, review and assess necessity for testing.

3

Once your test request has been approved, you will receive notification that your at-home saliva collection kit has shipped to your address.

4

Complete sample collection and return the specimen back to the Phosphorus Diagnostics laboratory via a prepaid FedEx Return Shipment Pack.

5

Once ready, you and your ordering healthcare provider will receive a notification via email with instructions for obtaining your results. You may choose to follow-up with your provider.

Please visit phosphorus.com/covid-19 to learn more about COVID-19 and what Phosphorus is doing to increase testing availability during this global pandemic.

Phosphorus

At-Home Collection Guide

For Prescription, In Vitro Diagnostic and Emergency Use Authorization ONLY

Collection Kit Contents

Phosphorus Kit Box

Prepaid FedEx Return Shipment Pack

Biohazard Bag

Absorbent Pouch

DNA Genotek OGD-510 Saliva Collection Device

For sample collection tube

Patient First & Last Name
Date of Birth (mm/dd/yy)
Collection Date and Time (EST):

Sticky Notes for Sample Labeling

Pre-populated Test Requisition Form

This At-Home Patient Collection Guide

Sample Collection Instructions


Collection Note

Do NOT eat, drink (including water), smoke, or chew gum for 30 minutes before giving a saliva sample.

To ensure timely receipt of intact specimen, be sure to drop your packet with FedEx the same day that you collect your sample. Your sample must be received within 56 hours of specimen collection.


Instructions for Collecting and Returning Saliva Specimen

1

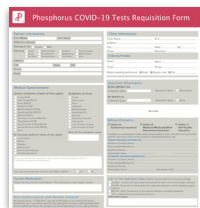Wash hands thoroughly for 20 seconds with soap upon receipt of package. Confirm that the information on your pre-populated test requisition form is correct.

2

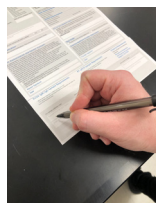Read and sign the test authorization and consent on the pre-populated requisition form to be returned.

3

Write your name, date of birth and date/time of specimen collection on two sticky notes. Place one on the sample collection tube; place the other in a space provided on the requisition form.

4


Wash your hands thoroughly for 20 seconds with soap. Follow the instructions provided with the saliva collection device. After collecting your sample, wash your hands again thoroughly for 20 seconds with soap.

5


Place the filled sample tube into absorbent pouch and then with signed requisition form into the biohazard bag and seal the bag. Press firmly along closure to securely zip. Wash your hands thoroughly for 20 seconds with soap after placing sample in the biohazard bag. Place the biohazard bag in the Phosphorus kit box.

6

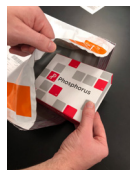Wash your hands thoroughly for 20 seconds with soap before placing the Phosphorus box within the prepaid UN3373 FedEx Return Shipment Pack. Drop the package at the nearest FedEx store or call 1.800.463.3339 to schedule a same day pickup.

7

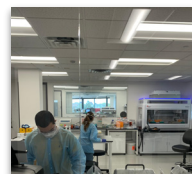Your sample will be shipped to the Phosphorus Diagnostics laboratory in Secaucus, NJ, and results will be available within 72 hours. Your results will be sent directly to you. A copy of your results will also be sent to your provider.

8






You will have the option to talk to a physician should you have any questions about your results.

Questions about collecting your specimen? Please get in touch.


To view video guided instructions on collection, please visit: <https://www.youtube.com/watch?v=ql2dc39nZ94&feature=youtu.be>
To view a copy of the collection instructions, please visit: <https://www.dnagenotek.com/US/support/collection-instructions/oragene-dx/OGD-510andOGD-610.html>

Phosphorus | support@phosphorus.com | phosphorus.com

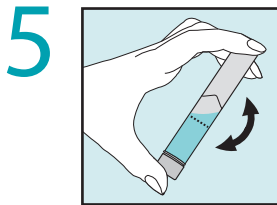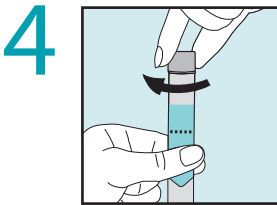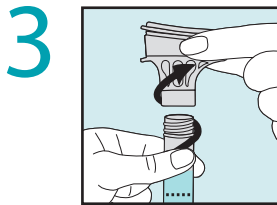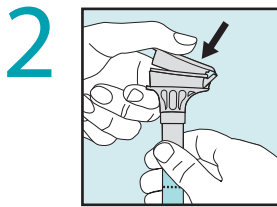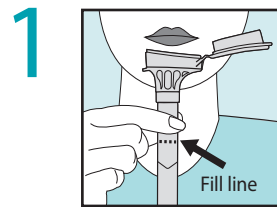DNAgenotek

Made in Canada
DNA Genotek Inc.
3000 - 500 Palladium Drive
Ottawa, ON, Canada K2V 1C2
Tel.: +1.613.723.5757
info@dnagenotek.com
www.dnagenotek.com

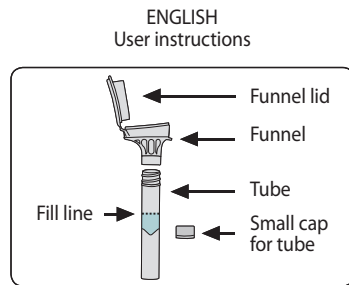Read all instructions
prior to collection.

Intended use: For the collection of human DNA from saliva samples.

Contents: Kit contains stabilizing liquid.

Warnings and precautions: Rx only. Wash with water if stabilizing liquid comes in contact with eyes or skin. Do NOT ingest. See MSDS at www.dnagenotek.com.

Small cap may pose a choking hazard.

Storage: 15°C / 30°C

Summary and explanation of the kit:

Oragene-Dx is a self-collection kit that provides the materials and instructions for collecting and stabilizing saliva specimens.

Collection precautions:

Do NOT eat, drink, smoke or chew gum for 30 minutes before giving your saliva sample.

Do NOT remove the plastic film from the funnel lid.

Procedure:

Most people take between 2 and 5 minutes to deliver a saliva sample following steps 1 to 5.

1 Spit into funnel until the amount of liquid saliva (not bubbles) reaches the fill line shown in picture #1.

2 Hold the tube upright with one hand. Close the funnel lid with the other hand (as shown) by firmly pushing the lid until you hear a loud click. The liquid in the lid will be released into the tube to mix with the saliva. Make sure that the lid is closed tightly.


3 Hold the tube upright. Unscrew the funnel from the tube.

4 Use the small cap to close the tube tightly.

5 Shake the capped tube for 5 seconds. Discard the funnel.

Label legend:

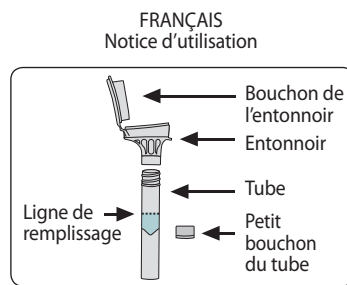- Consult package insert
- Collect saliva by (Use by)
- In vitro diagnostic medical device
- Catalog number
- CE Marking
- Caution, consult instructions for use
- Storage instructions
- Authorized Representative
- Manufacturer
- Lot number

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all rights reserved.