



# LabCorp At Home COVID-19 Test Home Collection Kit



1

**Read through the instructions completely before collecting your sample.** To see an instructional video or report a missing item in the kit, visit [www.LabCorp.com/AtHomeKit](http://www.LabCorp.com/AtHomeKit).

Las instrucciones en español están disponibles en línea en [www.LabCorp.com/AtHomeKit](http://www.LabCorp.com/AtHomeKit).



2

**Your specimen must be shipped Monday through Friday and on the same day that you collect it.** Do not collect your specimen until you confirm pickup times of the drop box or have confirmed your scheduled pickup.


Visit <https://bit.ly/fedexdropbox> to view FedEx drop box locations and pickup schedules. For more information, refer to the back of this instruction card.

3



**Wash and dry hands before opening the kit. Open your kit and place all the contents on a clean, dry surface.**

6

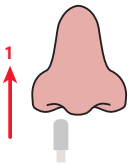


**Insert the tip of the swab into one nostril.** The swab does not need to be inserted far—insert just until the tip of the swab is no longer visible. Rotate the swab in a circle around the entire inside edge of your nostril at least 3 times.

4



**Take the swab out of its package.** Do not touch the tip of the swab with your hands.

7

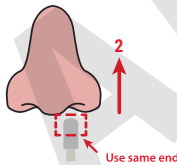


**Take the swab out of your nostril. Using the same end of the swab, repeat step 6 in your other nostril.**

5



**Screw off the top of the collection tube.** Hold swab in one hand and collection tube in the other, being careful not to spill the liquid. Do not drink the liquid.

8

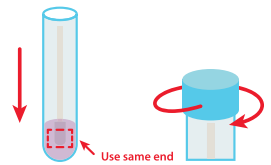

**Remove the swab from your second nostril and place in the collection tube.** The end of the swab that went into your nose should be placed into the tube first so that it sits down in the liquid. Screw the top of the collection tube back on.

9



**Wash and dry hands thoroughly again. Insert collection tube into the biohazard specimen bag.** Seal the biohazard specimen bag. Fold the specimen bag in half.

10



**Fill out the Specimen Confirmation Form.** Provide the date and time of collection, and complete the Personal Information section.

12



**Return the postage-paid, pre-addressed FedEx return pack to LabCorp.** You can deliver it to a FedEx drop box or schedule a pickup. It's important to ship your specimen Monday through Friday and on the same day you collect it.

13



**Results will be available in your LabCorp Patient™ account 1 to 2 days after your specimen arrives at the lab.** If you don't have an account, visit [patient.LabCorp.com](http://patient.LabCorp.com) to create one. Results will also be delivered to the ordering healthcare provider.

11



**Place the specimen bag and Specimen Confirmation Form into the white shipping box and close the lid.** Place the shipping box into the FedEx return pack. Remove the adhesive cover strip and seal the FedEx return pack closed.

## Ship Your Specimen to LabCorp via FedEx

**Your specimen must be shipped Monday through Friday on the same day that you collect it.**

You have 2 options for shipping your sample.

### FedEx Drop Box

### Schedule a FedEx Pickup

**Schedule** your FedEx pickup for a weekday (Monday through Friday).

Call FedEx at **800-463-3339** and answer some questions.

**FedEx:** "What are you calling about?" **You:** "FedEx Express" return pickup"

**FedEx:** "...a pickup request for a return. To arrange a pickup for today or

tomorrow, depending on availability, say "pickup." **You:** "Pickup"

**FedEx:** "Is the word "stamp" written anywhere on your return label?" **You:** "No"

**FedEx:** "What is the tracking number?" **You:** Provide the tracking number

(TRK#) from the return label

**FedEx:** "What kind of return are you calling about?" **You:** "FedEx Express"

**Provide additional information** as requested. You may be asked for the

pickup location address, number of packages to be picked up, day and time

window for the pickup, or if your shipment is currently ready.

**Note** the confirmation pickup number provided. If you are asked for a credit

card or account number, state that this is a FedEx Express return pickup and

provide the tracking number (TRK#) from the return label.

**Do not collect your specimen until you confirm pickup times of the drop box or have confirmed your scheduled pickup.** Visit <https://bit.ly/fedexdropbox> to view FedEx drop box locations and pickup schedules.



# Let's get started.

Las instrucciones en español están disponibles en línea en [www.LabCorp.com/AtHomeKit](http://www.LabCorp.com/AtHomeKit).

Kit includes:



**Intended use:** to collect a nasal swab for purpose of in vitro diagnostic testing in a CLIA laboratory. Rx only

IVD

Distributed by LabCorp  
531 South Spring Street  
Burlington, NC 27215

