



Simoa® SARS-CoV-2 N Protein Antigen Test

INSTRUCTIONS FOR USE

IFU 0002

Draft version 10.0

Quanterix Corporation
900 Middlesex Turnpike
Billerica, MA 01821

Customer Support
1-877-786-8749
Techsupport@quanterix.com

For use under an Emergency Use Authorization (EUA) Only

Prescription Use Only

For In Vitro Diagnostic Use Only

Read this package insert prior to use. Package insert instructions must be carefully followed. Reliability of assay results cannot be assured if there are deviations from the instructions in this package insert.

NAME

Simoa® SARS-CoV-2 N Protein Antigen Test

INTENDED USE

The Simoa SARS-CoV-2 N Protein Antigen Test is an automated paramagnetic microbead-based immunoassay intended for the qualitative detection of the nucleocapsid protein (N protein) antigen from SARS-CoV-2 in nasopharyngeal swab specimens collected in Huachenyang iClean Viral Transport Medium (VTM), CDC's formulation of VTM, normal saline, or phosphate buffered saline (PBS) from individuals who are suspected of COVID-19 by their healthcare provider within 14 days of symptom onset. Testing is limited to laboratories certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C. §263a, that meet the requirements to perform moderate or high complexity tests.

The Simoa SARS-CoV-2 N Protein Antigen Test does not differentiate between SARS-CoV and SARS-CoV-2.

Results are for the detection of SARS-CoV-2 nucleocapsid protein antigen. Antigen is generally detectable in nasopharyngeal swab specimens during the acute phase of infection. Positive results indicate the presence of viral antigens, but clinical correlation with patient history and other diagnostic information is necessary to determine infection status. Positive results do not rule out co-infection with other viruses. The antigen detected may not be the definite cause of disease. Laboratories within the United States and its territories are required to report all results to the appropriate public health authorities.

Negative results should be treated as presumptive and confirmation with a molecular assay, if necessary, for patient management, may be performed. Negative results do not rule out SARS-CoV-2 infection and should not be used as the sole basis for treatment or patient management decisions, including infection control decisions. Negative results should be considered in the context of a patient's recent exposures, history and the presence of clinical signs and symptoms consistent with COVID-19.

The Simoa SARS-CoV-2 N Protein Antigen Test is intended for use by trained clinical laboratory personnel specifically instructed and trained in vitro diagnostic procedures. The Simoa SARS-CoV-2 N Protein Antigen Test is intended for use on the Simoa HD-X Analyzer. The Simoa SARS-CoV-2 N Protein Antigen Test is only for use under the Food and Drug Administration's Emergency Use Authorization.

SUMMARY AND EXPLANATION OF TEST

Severe Acute Respiratory Syndrome Coronavirus 2 (SARS-CoV-2) is a recently identified coronavirus strain responsible for the Coronavirus Disease 2019 (COVID-19) and pandemic. SARS-CoV-2 emerged in China in December 2019 and is transmitted mainly through droplets and surface contact routes. The virus infects human cells through interaction between angiotensin converting enzyme 2 (ACE2) on respiratory cells and spike or S-protein on the outer envelope of the virion particle. COVID-19 affects people in different ways. Symptoms can include signs and symptoms of acute respiratory illness, such as fever, cough, shortness of breath, but the infection can also be asymptomatic. Symptomatic, pre-symptomatic and asymptomatic infected individuals can all be sources for viral transmission. The current gold standard for diagnosis of SARS-CoV-2 infection is real-time reverse transcription polymerase chain reaction (rRT-PCR), which detects the presence of SARS-CoV-2 nucleic acid material in upper respiratory specimens, such as nasopharyngeal swab and oropharyngeal swab. In contrast, SARS-CoV-2 protein antigen detection directly detects the presence or absence of proteins from the SARS-CoV-2 virus. One of these proteins, nucleocapsid protein (or N protein) is elevated in respiratory fluids during the initial acute phase of the infection. This N protein elevation can be used to detect a SARS-CoV-2 infection in its earliest stages, prior to the development of a significant anti-SARS-CoV-2 antibody response.¹

PRINCIPLES OF THE PROCEDURE

The Simoa SARS-CoV-2 N Protein Antigen Test is a 2-step microbead-based sandwich ELISA that uses single molecule array (Simoa) technology.² In the first step, anti-N protein antibody coated paramagnetic capture beads, sample, and biotinylated anti-N protein detector antibody are combined. Nucleocapsid protein molecules present in the sample are captured by the anti-N protein capture beads and labeled with biotinylated detector. After washing, a conjugate of streptavidin-β-galactosidase (SBG) is mixed with the capture beads. SBG binds to the biotinylated detector antibodies, resulting in enzyme labeling of captured nucleocapsid protein. Following a second wash, the capture beads are resuspended in a resorufin β-D-galactopyranoside (RGP) substrate solution for signal generation. Digital processing occurs when beads are transferred to the Simoa array disc which is composed of microarrays of femtoliter reaction wells. Individual capture beads are then sealed within microwells in the array through the addition of oil, which forms a liquid seal trapping the labeled immunocomplexes and RGP within the wells. If nucleocapsid protein from the sample has been captured and labeled, the β-galactosidase hydrolyzes the RGP substrate into a fluorescent product that provides the signal for digital counting. The fraction of bead-containing microwells counted with an enzyme is converted into 'average enzymes/bead' (AEB). AEB values are converted into N protein concentration in unknown samples by interpolation from a calibration curve obtained by 4-parameter logistical regression fitting. Total time to first result on a single sample is 80 minutes. Time to perform 96 tests is approximately two hours and 30 minutes.

For additional information on system and assay technology, refer to the *Simoa HD-X Analyzer User Guide (EUA Edition)*.

REAGENTS

Reagent Kit – Materials Provided

Simoa SARS-CoV-2 N Protein Antigen Test kit

Bead Reagent	1 bottle (4.4 mL)	Anti-N protein (mouse monoclonal) antibody coated capture beads in Tris buffer with a protein stabilizer (bovine) and a surfactant. Preservative: ProClin 300.
Detector Reagent	1 bottle (4.0 mL)	Biotinylated anti-N protein (mouse monoclonal) antibody in phosphate buffer with a protein stabilizer (bovine) and a surfactant. Preservative: ProClin 300.
SBG Reagent	1 bottle (12.3 mL)	Conjugate of streptavidin-β-galactosidase (SBG) in phosphate buffer with a protein stabilizer (bovine). Preservative: ProClin 300.
Sample Diluent	2 bottles (14.5 mL ea)	Phosphate buffer with a protein stabilizer (bovine), a heterophilic blocker, and a surfactant. Preservative: ProClin 300.
RGP Reagent	2 bottles (3.4 mL)	Resorufin β-D-galactopyranoside (RGP) in phosphate buffer with a surfactant.
Calibrators A-H (0 plus 7 levels)	8 vials (1 mL each)	Recombinant N protein in phosphate buffer with a protein stabilizer (bovine), a surfactant. Preservative: ProClin 300.
Positive Control 1	1 vial (0.5 mL each)	Recombinant N protein in Sample Diluent.
Positive Control 2	1 vial (0.5 mL each)	Recombinant N protein in Sample Diluent.
Negative Control 3	1 vial (0.5 mL each)	Phosphate buffered saline with sodium azide as a preservative.

Materials Required but Not Provided

- Simoa HD-X Analyzer, Simoa software v3.0 (Item # 103385)
- Simoa HD-X System Wash Buffer 1 (Item # 101619)
- Simoa HD-X System Wash Buffer 2 (Item # 100205)
- Simoa HD-X Sealing Oil (Item # 100206)
- Simoa cuvettes (Item # 103160, 500 ct box)
- Simoa disposable pipettor tips (Item # 101681)
- Simoa Discs (Item # 100001)
- Simoa microplate shaker (Item # 102899)
- X-Pierce™ XP-100 plate seals (K-1080, Sigma Aldrich item # Z722502)
- Simoa 96 well assay plates (Item # 101457)
- Standard nasopharyngeal swab collection kit — Users should obtain samples in accordance with the IFU procedures on sample collection.

WARNINGS AND PRECAUTIONS

For in vitro diagnostic and laboratory professional use. For emergency authorization use only.

- This product has not been FDA cleared or approved; but has been authorized by FDA under an Emergency Use Authorization (EUA) for use by laboratories certified under the CLIA that meet the requirements to perform moderate or high complexity tests.
- This product has been authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens.
- This test is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostic tests for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

Safety Precautions

- **CAUTION:** This product requires the handling of human specimens. It is recommended that all human-sourced materials be considered potentially infectious and be handled in accordance with the OSHA Standard on Bloodborne Pathogens.³ Biosafety Level 2⁴ or other appropriate biosafety practices should be used for materials that contain or are suspected of containing infectious agents.
- Simoa reagents contain methylisothiazolones, which are components of ProClin and are classified per applicable European Community (EC) Directives as: Irritant (Xi). The following are the appropriate Risk (R) and Safety (S) phrases.



R43 May cause sensitization by skin contact.

S24 Avoid contact with skin.

S35 This material and its container must be disposed of in a safe way.

S37 Wear suitable gloves.

S46 If swallowed, seek medical advice immediately and show this container or label.

- For a detailed discussion of safety precautions during instrument operation, refer to the *Simoa HD-X Analyzer User Guide*.

Handling Precautions

- Proper specimen collection, storage, and transport are critical to the performance of this test. Inadequate or inappropriate specimen collection, storage, and transport may yield false negative test results.
- The reagents are one-time use; any remaining material should be discarded after the completion of assay run.
- Calibrators and controls are one-time use; any remaining material should be discarded after completion of assay run.

- Fresh RGP reagent (prepared as described in the Assay Procedure) should be used with every run.
- Do not use reagent kits beyond the expiration date. When stored and handled as directed, reagents and calibrator are stable until the expiration date
- Do not pool reagents within a kit or between reagent kits.
- Do not attempt to reuse tips, cuvettes, or Simoa Discs, as this will cause significant data quality deterioration and may result in incorrect results.
- Dispose of test components and clinical specimens in accordance with Federal, State and Local regulatory requirements. The Safety Data Sheets (SDS) for kit components are available upon request. Contact Quanterix Customer Service (1-877-786-8749).

Shipping and Storage Instructions

Simoa SARS-COV-2 N Protein Antigen Test reagents are shipped on cold packs. If the reagents arrive at room temperature or frozen, reagent integrity may be suspect. Contact Quanterix Customer Service (1-877-786-8749).



- Simoa SARS-COV-2 N Protein Antigen Test reagents must be stored at 2–8°C in an upright position.



- Simoa SARS-COV-2 N Protein Antigen Test Calibrators and Controls are shipped on dry ice. If these components arrive in a partially frozen or unfrozen state, their integrity may be suspect. Contact Quanterix Customer Service.



- Simoa SARS-COV-2 N Protein Antigen Test Calibrators and Controls must be stored at -80°C and should be kept upright.
- When stored and handled as directed, reagents and calibrators are stable until the expiration date.

- Total time for RGP reagent on board the instrument should be limited to no more than eight hours.

Indications of Reagent Deterioration

If a control sample returns a concentration value out of the expected range, this may indicate deterioration of reagents or errors in technique. Associated test results may be invalid and may require retesting. Assay recalibration may be necessary. Refer to the Calibration and Quality Control Procedures section of this document.

SPECIMEN COLLECTION AND PREPARATION FOR ANALYSIS

- Use sterile rayon, foam, polyester or flocked flexible-shaft NP swabs to collect a nasopharyngeal sample.
- When collecting any type of sample, follow universal collection precautions and guidelines according to your organization. For specimen collection of nasopharyngeal swabs, follow the Centers for Disease Control and Prevention (CDC) Swab Collection Guidelines and swab manufacturers' recommendations.
- To ensure proper collection, the swab should be passed a distance that is halfway from the nose to the tip of the ear. This is about half the length of the swab. The swab should travel smoothly with minimal resistance; if resistance is encountered, withdraw the swab a little without taking it out of the nostril. Then elevate the back of the swab and move it forward into the nasopharynx.
- Immediately place the swab in 3 mL of transport medium and cap the specimen tube. Huachenyang iClean Viral Transport Medium, Centers for Disease Control Viral Transport Medium, normal saline, and phosphate buffered saline have been validated as transport media for the Simoa SARS-COV-2 N Protein Antigen Test.
- After collection, store the specimen at 2–8°C until ready for testing.

- Stability of nasopharyngeal specimens in saline has been established for up to four hours at 30°C, three days at 2-8°C, and 14 days at -80°C. Nasopharyngeal specimens in saline are stable to up to seven freeze/thaw cycles.
- Specimens thawed after frozen storage must always be mixed THOROUGHLY by low speed vortexing or inverting 10 times. Visually inspect the specimens to confirm homogeneity.
- If specimens are to be transported prior to testing, transportation must comply with all applicable regulations for the transport of etiologic agents. Samples to be transported outside the collection facility must be transported in insulated cold shipping containers with frozen gel packs sufficient to ensure cold conditions through the duration of shipment. The receiving laboratory must confirm the integrity of the samples upon receipt by confirming the frozen state of the gel packs. Upon receipt, transfer the samples to 2-8°C storage. Samples to be transported within a facility may be maintained at room temperature (\leq 30°C) if transport time does not exceed 4 hours. If longer transport times are anticipated, follow the instructions above for shipment outside the facility.
- If samples are frozen, thaw completely prior to processing. Centrifuge all specimens prior to assay. Centrifugation conditions should be sufficient to efficiently remove particulate matter and to clarify the sample, for example 5 minutes at 10,000 g. Remove the swab prior to centrifugation.
- For optimal results, inspect all samples for bubbles immediately before placing on the instrument. Remove bubbles with an applicator stick prior to analysis. Use a new applicator stick for each sample to prevent cross-contamination.
- Use caution when handling patient specimens to prevent cross-contamination. Use of disposable pipettes or pipette tips is recommended.
- Specimens with obvious microbial contamination should not be used.

PROCEDURE

Assay Procedure

- The Simoa SARS-CoV-2 N Protein Antigen Test assay definition must be downloaded from the Quanterix customer portal website (portal.quanterix.com) and installed on the Simoa HD-X Analyzer prior to performing the assay. Do not open or attempt to edit the file or assay results could be compromised.
- Calibrators, controls, and samples must be allowed to come to room temperature and mixed thoroughly before loading onto the Simoa HD-X Analyzer.
- Simoa reagents excluding RGP and Bead Reagent can be loaded onto the HD-X straight from refrigerated storage or after equilibration to room temperature.
- Solubilize RGP fully by heating at 30–37°C with constant vigorous shaking on the Simoa microplate shaker set to 800 rpm for a minimum of 30 minutes and a maximum of four hours. One RGP bottle is required for up to 48 tests. Two RGP bottles are required for 49–96 tests.
- Immediately before loading on the HD-X Analyzer, the Bead Reagent bottle must be mixed to resuspend the capture beads that may have settled. To resuspend the beads, vortex for a minimum of 30 seconds. Re-vortex if there is a delay of more than five minutes loading the Bead Reagent on the HD-X.
- Note:** The bead diluent is formulated with an antifoam agent, but vortexing can still create foaming. If the foam does not dissipate within a few minutes, remove excess foam with a pipette prior to loading bead reagent onto the Simoa HD-X Analyzer.
- Set up the assay run on the instrument (see the *Simoa HD-X Analyzer User Guide*).
 - Remove caps from reagent bottles and load the Simoa SARS-CoV-2 N Protein Antigen Test reagents (Bead Reagent, Detector Reagent, SBG Reagent, Sample Diluent) into the reagent bay.
 - Load samples, calibrators, controls, and RGP into the sample bay. (Note: Two bottles of RGP are needed per 96-test run. One bottle of RGP may be used for multiple smaller runs up to 48 tests, but total time on board the instrument should be limited to no more than eight hours.

- In the run set up screen, specify “neat” protocol for calibrators and “dilution” protocol for controls and samples. The “dilution” protocol performed a 1:4 dilution of the controls and samples.
- Replenish consumables and system resources as needed prior to initiating the run, as described in the *Simoa HD-X Analyzer User Guide*.
- Initiate the run.
- Upon completion of the run, generate and export the Clinical Results Report and place into the patient record as per the laboratory's standard operating procedures. Only qualitative results should be reported in the Clinical Results Report, with limitations as indicated in the Interpretation of Results Table included below.

Sample Inputs and Volumes

- The SARS-CoV-2 N Protein Antigen Test was validated with Quanterix-supplied 96-well plates as the sample input. Each replicate test consumes 25 µL of sample. Multiple replicate tests of the same sample may be run from one well. The minimum sample volume in the well depends on the number of replicates to be run and the required dead volume of the well. The dead volume of a well is 30 µL, so the volume of sample to pipette for a single replicate is 55 µL. The maximum recommended sample volume for 96-well plates is 300 µL.
- If samples are left on board the instrument or pipetted samples are left on a lab bench for more than an hour, evaporation effects may influence the results depending on the volume of sample. Use of plate seals is recommended to prevent evaporation. X-Pierce Sealing Films (cat # XP-100) are the only plate seals compatible with the Simoa HD-X Analyzer. When placing seals, care must be taken to center the black circular marks over all plate wells. Once seal is placed, do not tip the plate or sample well contents maywick up onto the seal and cause cross-contamination of wells.

Calibration

- To perform a SARS-CoV-2 N Protein Antigen Test calibration, test Calibrators A through H in duplicate. All levels of SARS-CoV-2 N Protein Antigen Test Controls should be tested in duplicate to evaluate the assay calibration. All assay control values should be within their expected concentration ranges for a valid calibration curve. Refer to the Certificate of Analysis for expected control values and ranges. If the results from one or more of the controls are outside their expected ranges, the assay calibration may not be valid, and a re-calibration is recommended.
- The HD-X Analyzer is capable of storing a calibration curve for analysis of samples on subsequent batch runs. Calibration curve storage has not been validated for the SARS-CoV-2 N Protein Antigen Test.
- The SARS-CoV-2 N Protein Antigen Test utilizes a 4 Parameter Logistic Curve fit data reduction method to generate a calibration curve. Specimen results are interpolated from the calibration curve.

Preparing Calibrators

- Calibrators should be brought to room temperature prior to pipetting. Do not heat the vial to accelerate thawing.
- When the solution is fully thawed, THOROUGHLY mix by multiple gentle inversions or vortexing. Frozen protein solutions can partition during freezing, so complete mixing of thawed material is critical for accurate calibrators.
- Each replicate test consumes 100 µL of calibrator. Multiple replicate tests of the same calibrator may be run from the same well. The minimum calibrator volume in the well depends on the number of replicates to be run and the required dead volume for the well. The dead volume of the well is 30 mL, so replicates of two require 230 µL of each calibrator level.

Preparing Controls

- Controls should be brought to room temperature prior to pipetting. Do not heat the vial to accelerate thawing.
- When the solution is fully thawed, THOROUGHLY mix by multiple gentle inversions or vortexing. Frozen protein solutions can partition during freezing, so complete mixing of thawed material is critical for accurate controls.

QUALITY CONTROL PROCEDURES

- Follow the specific quality control procedures in your laboratory.
- The SARS-CoV-2 N Protein Antigen Test Negative Control and Positive Controls should be included in every run to assess run validity.
- Control values and ranges listed in the lot-specific Certificate of Analysis should be considered as guides. It is recommended that each laboratory establish statistically based control values and ranges with a sufficiently powered study. For guidance, it is recommended to consult with Clinical and Laboratory Standards, Institute (CLSI) Guideline C24, 4th ed., or other published guidelines for quality control recommendations.
- If quality control results do not meet acceptance criteria defined by your laboratory, sample results may be suspect. Follow the established quality control procedures for your laboratory. For troubleshooting information, refer to the *Simoa HD-X Analyzer User Guide*.

Interpretation of Results

- The assay clinical cutoff for the presence of N protein is 3.0 (arbitrary units). The following table provides guidance in interpreting the information displayed in the Clinical Results Report following completion of the Simoa SARS-CoV-2 N Protein Antigen Test. (Note: 'AEB' refers to 'average enzymes/bead,' the digital signal of the assay).

#	Result Displayed	Reason for Result	Interpretation
1	Numerical result <3.0	Undetectable SARS-CoV-2 N protein antigen. AEB within calibration range.	Negative for SARS-CoV-2 N protein antigen.
2	Numerical result ≥3.0	SARS-CoV-2 N protein antigen detected. AEB within calibration range.	Positive for SARS-CoV-2 N protein antigen.
3	"Undetectable"	AEB below calibration range.	Negative for SARS-CoV-2 N protein antigen.
4	"Positive, out of range"	AEB above calibration range.	Positive for SARS-CoV-2 N protein antigen.
5	"Positive, too much fluor"	Fluorescence exceeds optical range.	Positive for SARS-CoV-2 N protein antigen.
6	"NaN" (Not a Number)	Invalid.	Retest

Non-numerical result scenarios #3-5 are valid for determining the presence or absence of viral antigens in the sample. "NaN" results indicate one or more system errors (such as an aspiration error) during processing. Inspect the Flags field of the Clinical Results Report for additional information and address the cause of any flags where required before repeating a sample. If non-numerical results recur after repeat testing with recommended sample dilutions, contact Quanterix Customer Support.

Note: Only the qualitative results from the Interpretation table should be reported. Semi-quantitative numerical results have not been clinically or analytically validated and may not correlate with patient disease status, duration of illness or severity of illness. Semi-quantitative results have not been demonstrated to correlate with the success or failure of any therapeutic interventions and should not be used to guide clinical management. **Note:** Results of this assay should always be interpreted in conjunction with the patient's medical history, clinical presentation, and other findings.

Note: For a description of the messages that may appear in the Flags field, refer to the *Simoa HD-X Analyzer User Guide*.

LIMITATIONS OF THE PROCEDURE

- For in vitro diagnostic use on the HD-X under Emergency Use Authorization Only. Other Quanterix instrument platforms are not authorized for diagnostic testing.
- Only qualitative results should be reported. Semi-quantitative numerical results have not been clinically or analytically validated and may not correlate with patient disease status, duration of illness or severity of illness. Semi-quantitative results have not been demonstrated to correlate with the success or failure of any therapeutic interventions and should not be used to guide clinical management.

- This test will indicate the presence of SARS-CoV-2 nucleocapsid protein antigen in the specimen from both viable and non-viable SARS-CoV-2 virus. Test performance depends on the amount of virus (antigen) in the sample and may or may not correlate with viral culture results performed on the same sample.
- Failure to follow the instructions for use may adversely affect test performance and/or invalidate the test result.
- Test results should be considered in the context of all available clinical and diagnostic information, including patient history and other test results.
- Positive test results do not differentiate between SARS-CoV and SARS-CoV-2.
- Positive test results do not rule out co-infection with other pathogens.
- Negative test results are not intended to rule in other non-SARS viral or bacterial infections.
- Negative results should be treated as presumptive and confirmation with a molecular assay, if necessary, for patient management, may be performed.
- A negative test result may occur if the level of antigen in a sample is below the detection limit of the test or if the sample was collected or transported improperly.
- If the differentiation of specific SARS viruses and strains is needed, additional testing, in consultation with state or local public health departments, is required.
- Performance has only been established with nasopharyngeal specimens. Other specimen types have not been evaluated and should not be used with this assay.
- Bacterial contamination of specimens may affect the test results.
- Therapeutic doses of biotin can interfere with assays that utilize biotinylated reagents.

Conditions of Authorization for the Laboratory

The Simoa SARS-CoV-2 N Protein Antigen Test Letter of Authorization, along with the authorized Fact Sheet for Healthcare Providers, the authorized Fact Sheet for Patients, and authorized labeling are available on the FDA website: <https://www.fda.gov/medical-devices/coronavirus-disease-2019-covid-19-emergency-use-authorizations-medical-devices/vitro-diagnostics-euas>.

Authorized laboratories* using the Simoa SARS-CoV-2 N Protein Antigen Test, must adhere to the Conditions of Authorization indicated in the Letter of Authorization as listed below:

- Authorized laboratories using the Simoa SARS-CoV-2 N Protein Antigen Test must include with test result reports, all authorized Fact Sheets. Under exigent circumstances, other appropriate methods for disseminating these Fact Sheets may be used, which may include mass media.
- Authorized laboratories using the Simoa SARS-CoV-2 N Protein Antigen Test must use the Simoa SARS-CoV-2 N Protein Antigen Test as outlined in the Instructions for Use. Deviations from the authorized procedures, including the authorized instruments, authorized clinical specimen types, authorized control materials, authorized other ancillary reagents and authorized materials required to use the Simoa SARS-CoV-2 N Protein Antigen Test are not permitted.
- Authorized laboratories that receive the Simoa SARS-CoV-2 N Protein Antigen Test must notify the relevant public health authorities of their intent to run the Simoa SARS-CoV-2 N Protein Antigen Test prior to initiating testing.
- Authorized laboratories using the SARS-CoV-2 N Protein Antigen Test must have a process in place for reporting test results to healthcare providers and relevant public health authorities, as appropriate.
- Authorized laboratories must collect information on the performance of the Simoa SARS-CoV-2 N Protein Antigen Test and report to DMD/OHT7-OIR/OPEQ/ CDRH (via email: CDRH-EUA-Reporting@fda.hhs.gov) and Quanterix Corporation (at customerservice@quanterix.com) any suspected occurrence of false reactive or false non-reactive results and significant deviations from the established performance characteristics of the Simoa SARS-CoV-2 N Protein Antigen Test of which they become aware.
- All laboratory personnel using the Simoa SARS-CoV-2 N Protein Antigen Test must be appropriately trained in performing and interpreting the results laboratory and personal protective equipment when handling this kit and

use the product in accordance with the authorized labeling. Quanterix Corporation, authorized distributors, and authorized laboratories using the Simoa SARS-CoV-2 N Protein Antigen Test must ensure that any records associated with this EUA are maintained until otherwise notified by FDA. Such records will be made available to FDA for inspection upon request.

* The letter of authorization refers to, "Laboratories certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C. §263a, that meet requirements to perform moderate or high complexity tests" as "authorized laboratories."

PERFORMANCE CHARACTERISTICS

Limit of Detection (Analytical Sensitivity)

The Limit of Detection (LoD) was determined by evaluating different dilutions of gamma-inactivated SARS-CoV-2 virus added to an eluate pool from PCR-confirmed negative NP swabs. The LoD range was determined by assaying a 2-fold dilution series of virus (20 dilutions). The LoD was then determined by repeatedly assaying three dilutions of virus near the LoD estimated from the range finding study. The LoD was confirmed as the lowest concentration where at least 95% of 20 replicates (19/20) read positive relative to the clinical cutoff. The results of this study are summarized in the table below.

	Dilution 1	Dilution 2	Dilution 3
Number positive replicates/total	20/20	20/20	15/20
TCID ₅₀ /mL	0.35	0.31	0.28

The LOD in nasal pharyngeal specimens was confirmed as 0.31 TCID₅₀/mL.

Cross Reactivity and Microbial Interference

Cross reactivity and potential interference were evaluated by testing 30 commensal and pathogenic microorganisms (16 viruses, 12 bacteria, 2 fungi and pooled human nasal fluid) that may be present in the nasal cavity. Each of the organisms, viruses, and fungi were tested in triplicate in the absence or presence of gamma-inactivated SARS-CoV-2 virus in a pool of eluates from PCR-confirmed negative NP swabs. No cross-reactivity or interference was seen with the microorganisms presented in the table below with the exception of SARS-Coronavirus. When tested at a concentration of 1.0 E+05 CFU/mL in negative matrix, SARS-Coronavirus gave an assay signal equivalent to 57.1 units.

Potential Cross-Reactant		Test Concentration
Virus	Enterovirus Type 68 (2007 Isolate)	3.78E+05 TCID ₅₀ /mL
	Adenovirus Type 07 (Species B)	3.53E+04 TCID ₅₀ /mL
	Human Metapneumovirus 16 Type A1	9.50E+05 TCID ₅₀ /mL
	Influenza A H1N1 (New Cal/20/99)	2.88E+06 TCID ₅₀ /mL
	Influenza B (Florida/02/06)	3.53E+04 TCID ₅₀ /mL
	Parainfluenza Virus Type 1	2.28E+08 TCID ₅₀ /mL
	Parainfluenza Virus Type 2	2.88E+06 TCID ₅₀ /mL
	Parainfluenza Virus Type 3	1.65E+06 TCID ₅₀ /mL
	Parainfluenza Virus Type 4A	7.05E+06 TCID ₅₀ /mL
	Respiratory Syncytial Virus Type A	9.50E+05 TCID ₅₀ /mL
	Rhinovirus Type 1A	8.88E+04 TCID ₅₀ /mL
	Coronavirus (Strain: 229E)	1.04E+05 TCID ₅₀ /mL
	Coronavirus (Strain: OC43)	2.63E+05 TCID ₅₀ /mL
	Coronavirus (Strain: NL63)	4.25E+04 TCID ₅₀ /mL
	MERS-Coronavirus (Strain: FL/USA-2 Saudi 2014)	1.0 E+05 TCID ₅₀ /mL
	SARS-Coronavirus	1.0 E+05 PFU/mL
Bacteria	Bordetella pertussis A639	1.13E+09 CFU/mL
	Legionella pneumophila Philadelphia	1.88E+09 CFU/mL
	Mycobacterium tuberculosis H37Ra-1	6.86E+06 CFU/mL
	Mycoplasma pneumoniae M129	3.16E+07 CFU/mL
	Pseudomonas aeruginosa	8.44E+08 CFU/mL
	Streptococcus pyogenes Z018	1.64E+08 CFU/mL
	Streptococcus salivarius Z127	4.19E+07 CFU/mL
	Chlamydia pneumoniae (TWAR Strain)	1.49E-01 mg/mL
	Haemophilus influenzae type b	5.43E+06 CFU/mL
	Streptococcus pneumoniae 19F	2.26E+07 CFU/mL
	Staphylococcus aureus	1.48E+06 CFU/mL
	Staphylococcus epidermidis MRSE	1.21E+09 CFU/mL
Fungus	Pooled Human Nasal Fluid	15% v/v
	Candida albicans Z006	6.27E+07 CFU/mL
	Pneumocystis jirovecii	3.45E+07 CFU/mL

To estimate the likelihood of cross-reactivity with SARS-CoV-2 of organisms that were not available for wet testing, in silico analysis using the Basic Local Alignment Search Tool (BLAST) managed by the NCBI was used to assess the degree of protein sequence homology.

- No protein sequence homology was found between *M. tuberculosis* and *P. jirovecii* thus homology-based cross-reactivity can be ruled out.
- Homology for Coronavirus-HKU1 and MERS-Coronavirus was 39.1% and 50.3% across 76% and 88% of the sequence respectively. Cross-reactivity of the assay with Coronavirus-HKU1 and MERS-Coronavirus cannot be ruled out. For reference, sequence homology with SARS-Coronavirus was 90.5% across 100% of the sequence, leading to the cross reactivity observed in the wet testing described in the table above.

High Dose Hook Effect

High dose hook was evaluated with gamma-inactivated SARS-CoV-2 virus added to an eluate pool from negative NP swabs. No high dose hook effect was observed with increasing antigen concentration across the optical range of the instrument (to 137 TCID₅₀/mL). The highest concentration of SARS-CoV-2 virus tested to yield a positive qualitative result was 2.8x10⁵ TCID₅₀/mL.

Endogenous Interference

Potentially interfering substances, naturally present in respiratory specimens or that may be artificially introduced into the nasal cavity or nasopharynx, were evaluated in the absence or presence of gamma-inactivated SARS-CoV-2 virus in an eluate pool from negative NP swabs. The substances listed in the table below were found not to affect test performance.

Substance	Active Ingredient	Test Concentration
Afrin – nasal spray	Oxymetazoline	15%
Alkalol Allergy Relief	N/A	1:10 dilution
Blood (human)	Blood	5%
Chloraseptic	Benzocaine, Menthol	1.5 mg/mL
CVS Nasal Spray	Phenylephrine	15%
Flonase	Fluticasone	5%
Halls Relief Cherry Flavor	Menthol	0.8 g/mL
Homeopathic (Alkalol Nasal Wash)	N/A	1:10 dilution
Mupirocin	Mupirocin	10 mg/mL
Nasacort Allergy 24 hour	Triamcinolone	5%
NasalCrom Nasal Spray	Cromolyn sodium	15%
Naso GEL (NeilMed)	N/A	5%
Neo-Synephrine	Phenylephrine hydrochloride	5%
Oseltamivir	Oseltamivir	2.2 µg/mL
Purified mucin protein	Mucin protein	2.5 mg/mL
Rhinocort	Budesonide (Glucocorticoid)	5%
Saline nasal spray	Saline	15%
Sore Throat Phenol Spray	Phenol	15%
Tobramycin	Tobramycin	2.5 mg/mL
Zanamivir	Zanamivir	400 ng/mL
Zicam Cold Remedy	Galphimia glauca, Luffa operculata, Sabadilla	5%
Zicam Oral Mist	Zinc	5%

Biotin Interference

Biotin was added to pools of VTM eluants from negative NP swabs at test concentrations ranging from 412.5 to 3500 ng/mL. 3 positive pools were prepared by adding inactivated SARS-CoV-2 virus to levels 2-3X the clinical cutoff. The qualitative results from the negative and positive pools were unaffected by the added biotin. One of the 3 positive pools exhibited +18% bias in numerical assay results at biotin concentrations above 1750 ng/mL.

Clinical Performance

- The clinical performance characteristics of the Simoa SARS-CoV-2 N Protein Antigen Test was evaluated in a single-site retrospective study conducted in the U.S. The study evaluated 126 nasopharyngeal swab specimens in viral transport medium from patients suspected of SARS-CoV-2 infection within 14 days of symptom onset. Specimens were transported, stored and tested following validated specimen stability information in the instructions for use. The table below exhibits the positive and negative agreement between the Simoa antigen test and matched results from a highly sensitive EUA RT-PCR.

Simoa SARS-CoV-2 N Protein Antigen Test	Comparator method (RT-PCR)		
	Positive	Negative	Total
Positive	86	0	87
Negative	2	38	39
Total	88	38	126
Positive Agreement: 86/88 97.7% (95% CI: 92.03% - 99.72%)			
Negative Agreement: 38/38 100% (95% CI: 90.75% - 100.0%)			

- Patient demographics (gender, age, time elapsed since onset of symptoms) are available for the 126 samples used in the analysis. The table below shows the positive results by age.

Age	Simoa SARS-CoV-2 N Protein Antigen (n=126)		
	Total #	Positive	Prevalence
≤ 5 years	0	0	0 %
6 to 21 years	4	3	75.0 %
22 to 59 years	69	55	79.7 %
≥ 60 years	53	30	54.7 %

- The table below shows positive results from days since symptom onset.

Days since symptom onset	Cumulative RT-PCR positives	Cumulative Simoa antigen positives	PPA
0	8	8	100%
1	23	22	93.3%
2	34	33	100%
3	44	43	100%
4	53	52	100%
5	57	56	100%
6	61	60	100%
7	69	68	100%
8 to 14	88	86	94.7%

BIBLIOGRAPHY

- Porte L, Legarraga P, Vollrath V, et al. Evaluation of novel antigen-based rapid detection test for the diagnosis of SARS-CoV-2 in respiratory samples. *Int J of Infectious Diseases.* 2020; DOI:<https://doi.org/10.1016/j.ijid.2020.05.098>
- Rissin DM, Kan CW, Campbell TG, et al. Single-molecule enzyme-linked immunosorbent assay detects serum proteins at subfemtomolar concentrations. *Nat Biotech* 2010; 28:595–99.
- US Department of Health and Human Services. Biosafety in micro-biological and biomedical laboratories, 4th ed. Washington, DC: US Government Printing Office, May 1999.
- World Health Organization. Laboratory biosafety manual. Geneva: World Health Organization, 2004.
- Clinical and Laboratory Standards Institute. Protection of laboratory workers from occupationally acquired infections: Approved guideline, 3rd ed. CLSI Document M29-A3. Wayne, PA: Clinical and Laboratory Standards Institute, 2005.
- Clinical and Laboratory Standards Institute. Statistical Quality Control for Quantitative Measurement Procedures: Principles and Definitions: Approved guideline, 4th ed. CLSI Document C24. Wayne, PA: Clinical and Laboratory Standards Institute, 2017.
- Clinical and Laboratory Standards Institute. Evaluation of Precision of Quantitative Measurement Procedures; Approved Guideline—Third Edition. Wayne, PA: Clinical and Laboratory Standards Institute; 2014. CLSI Document EP05-A3.
- Clinical and Laboratory Standards Institute. Evaluation of Detection Capability for Clinical Laboratory Measurement Procedures, 2nd Edition; Approved Guideline. Wayne, PA: Clinical and Laboratory Standards Institute; 2017. CLSI Document EP17-A2.
- Clinical and Laboratory Standards Institute. Evaluation of the Linearity of Quantitative Measurement Procedures: A Statistical Approach; Approved Guideline. Wayne, PA: Clinical and Laboratory Standards Institute; 2003. CLSI Document EP06-A.

For additional information on Simoa technology, instrument operation maintenance, and data analysis, refer to the *Simoa HD-X Analyzer User Guide*. Simoa and Quanterix are trademarks of Quanterix Corporation. © 2014 Quanterix Corporation. All rights reserved.

ProClin is a trademark of Rohm & Haas Company.

Nalgene is a registered trademark of Thermo Fisher Scientific or its subsidiaries.

Quanterix Corporation
900 Middlesex Turnpike
Building 1
Billerica, MA 01821
www.quanterix.com
Email: Techsupport@quanterix.com
1-877-786-8749

January 2021

For use under an Emergency Use Authorization (EUA) Only.

Prescription Use Only. For In Vitro Diagnostic Use Only.

The Simoa SARS-CoV-2 N Protein Antigen Test is an automated paramagnetic microbead-based immunoassay intended for the qualitative detection of the nucleocapsid protein (N protein) antigen from SARS-CoV-2 in nasopharyngeal swab specimens collected in Huachenyang iClean Viral Transport Medium (VTM), CDC's formulation of VTM, normal saline, or phosphate buffered saline (PBS) from individuals who are suspected of COVID-19 by their healthcare provider within 14 days of symptom onset.

Study the HD-X Analyzer User Guide and Simoa SARS-CoV-2 N Protein Antigen Test Instructions for Use thoroughly before using these Quick Reference Instructions or performing a test. This is not a complete product insert. Refer to the Simoa SARS-CoV-2 N Protein Antigen Test Instructions for Use for Warnings, Precautions and Limitations.

The Simoa SARS-CoV-2 N Protein Antigen Test is intended for use by trained clinical laboratory personnel specifically instructed and trained in vitro diagnostic procedures. The Simoa SARS-CoV-2 N Protein Antigen Test is intended for use on the Simoa HD-X Analyzer. The Simoa SARS-CoV-2 N Protein Antigen Test is only for use under the Food and Drug Administration's Emergency Use Authorization.

EMERGENCY USE AUTHORIZATION – WARNING AND PRECAUTIONS

In the USA, this product has not been FDA cleared or approved; but has been authorized by FDA under an EUA for use by authorized laboratories; use by laboratories certified under the CLIA, 42 U.S.C. §263a, that meet requirements to perform moderate or high complexity tests. This product has been authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens. In the USA, - this test is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of the virus that causes COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

Please see the Simoa SARS-CoV-2 N Protein Antigen Test Instructions for Use for the full instructions for use, including Limitations, Warnings and Precautions

(https://www.quanterix.com/wp-content/uploads/2021/01/Simoa_SARS-CoV2_N_Protein_Antigen_Test_Instructions_for_Use.pdf).

1. Prepare Instrument

Note: The system must be shut down once every 24 hours.

- Turn on the HD-X Analyzer: Turn on the PC > Start the HD-X Analyzer (switch on the right side of instrument) > Open the Simoa® software. When instrument **initialization** is complete the system will say **Ready**.

- Pre-Run Maintenance:** Maintenance Tab > Check the **Start of Day Task** > click **Run Task** (about 20 minutes).

Note: If Start of Day Task was performed, but system has been idle for over 4 hours, run **Idle System Prime**.

- Adjust the Lot-Specific Calibrator Concentrations:** For each new lot of kits to be used, in the **Custom Assay Tab** > select the SARS-CoV-2 N Protein Antigen Assay > deselect **Read Only** > select **Plexes** > select the **SARS-CoV-2 N Protein Antigen Plex** > Modify the concentration value for each calibrator level according to the lot-specific CoA.

2. Set Up Assay

- Solubilize RGP:** Place RGP vial(s) on a temperature-controlled shaker set at 30°C for ≥30 minutes at a shaking speed of 800 rpm. One vial is sufficient for up to 48 tests. Do not shake for > 4 hours prior to use. Do not re-use RGP vials.
- Prepare Samples:** Equilibrate Calibrators, Controls, and Specimens to room temperature, then mix thoroughly by low speed vortexing or inverting 10 times. Visually inspect the specimens to confirm homogeneity. Centrifuge samples at 10,000g for 5 minutes to clear any debris.
- Determine volume required per well:** Each well of the assay plate requires 30 uL dead volume + the volume consumed during testing. Calibrators are run neat, and 100 uL is consumed per test. Controls and Specimens are run with an on-board 4X dilution, and 25 uL is consumed per test.
- Prepare Assay Plate:** Pipette required volumes of calibrators, controls and specimens into a Simoa 96-well microplate. Cover with an X-Pierce XP-100 sealing film.

Table 1: Sample volumes required. A minimum of two replicates should be run for calibrators and controls.

Replicates	Calibrators No Dilution Required	Controls and Specimens 4X Dilution on Instrument
Tests/well	Minimum volume per well (uL)	Minimum volume per well (uL)
1	not recommended	55
2	230	80
3	prep in separate wells	105

3. Load Instrument and Run Assay

Note: Only the RGP reagent needs to be warmed before use. Other assay reagents can go immediately from 2 – 8 °C storage on to the instrument.

1. **Vortex Bead Reagent:** Vortex Bead Reagent for about 30 seconds immediately before loading. If beads are not loaded within 5 minutes, vortex again.
2. **Load Reagents:**
 - a. **Load Bead, Sample Diluent, Detector and SBG reagents into Rack:** Remove all caps. Make sure beads are in one of the three shaking positions on the rack. Other reagents should be in non-shaking positions. Position bottles so barcodes are readable on the right side of the rack.
Select **Load Reagent Tab** > Select reagent lane > Enable the on-board barcode scanner > Insert reagent rack. Confirm that reagent barcodes have scanned correctly.
 - b. **Load RGP:** Select an RGP lane > Make sure on-board barcode scanner is enabled > Place RGP bottle(s) in RGP rack (marked with an O) and remove cap > Insert RGP rack. Confirm RGP barcode has scanned correctly.
 - c. Touch **Done Loading Reagents.**
3. **Create Plate Layout:**
 - a. **Setup Run Tab** > Assign **Batch name** > Assign **Plate Barcode** > Click **Enter** on your keyboard.
 - b. **Assign Calibrators:** Select **Assign Calibrators Tab** > Highlight the well containing Calibrator A > Select **CVDNAg Assay Definition** > Select Calibrator A from the **Select Calibrator** pop up > Click **Ascending/Descending** to populate the remaining calibrators in the plate column > Highlight all calibrator wells and set the **Replicates per well** to 2.
 - c. **Assign Samples:** Select **Assign Sample Tab** > Highlight all wells that contain controls/samples > Select **CVDNAg** in the **Assay** section > Set the desired **Replicates per well**.
 - d. If desired, select individual wells and enter sample IDs with the onscreen keyboard, USB keyboard, or handheld **barcode scanner**.
 - e. When setup is complete, click on **List View** to confirm selections. Ensure dilutions are set properly for calibrators (neat), and controls/specimens (dilution). Ensure calibrator concentrations match those on the lot-specific CoA.
 - f. Slide the prepared sample plate onto the plate rack with well A1 in the indicated position and the plate skirt under the catch screw. Clip plate into place. Insert plate rack, then touch **Done with Setup.**
4. **Check Liquid and Solid Consumables and Empty Waste**
 - a. **Load Liquid Consumables:** Check the System Resources tab. If indicated by yellow or red color coding, fill the DI Water and Wash Buffer 1 secondary containers and the Wash Buffer 2 container.
 - b. **Load Cuvettes, Tips, and Discs, if Required:**

Load Cuvettes: Cuvettes are added by placing a full stack of 50 in the chute. Do not allow stacks to twist while loading.

Note: After loading a stack, wait for the system to indicate **Ready** prior to loading the next stack.

Load Tips: **System Resource Tab** > Select **Solid Resources** > Click **Unlock Drawers** > Place tip racks carefully into place and ensure they sit level > In the software, tap twice on each position where a new tip rack was loaded. The tip positions in the rack diagram turn light blue. Only load full racks of tips. Insert Drawer.

Load Discs: **System Resource Tab** > Select **Solid Resources** > Click **Unlock Drawers** > Use the handheld scanner to scan the barcode on the wrapper > Remove blue base plate (from the old stack) from the disc pole in the drawer > Load the new stack on an empty disc pole > remove the wrapper and the top disc with the Quanterix logo. Insert Drawer.

c. Empty Solid and Liquid Waste: If necessary, empty solid and liquid waste containers located in the system bay.

Note: Only remove one liquid waste container at a time. Ensure the sensor has been re-inserted in the first waste container prior to removing the sensor from the second waste container.

5. **Start Run:** **System Resource Tab** > Select **All Resources** > Click **Start Run**. If button is not active, check for flags in **Resource Details**.
6. **Current Run Tab:** Use this screen to monitor progress of run. Run is finished when this tab reads 00:00 and status line at bottom left corner says **READY**.

4. Post-Run: Data Export, Maintenance, and Shut Down

1. **Remove Sample Plate, Reagents and RGP and dispose appropriately.** All kit reagents are single use only.
2. **Export the Clinical Results Report:** **History & Reports Tab** > Select **Reports** > Select **Clinical Results Report** > select appropriate batch name from pull down list > **Preview Report** > **Export**. Reports may be exported as .pdf files and saved to the instrument S:\ drive, a network location, or a USB. Reports should be exported promptly at the completion of each run.
3. **End of Day Maintenance:** After completing the final run of the day, navigate to **Maintenance Tab** > Check the **End of Day** task > Click **Run Task** (about 20 minutes).
4. **Shut down:** Shut down software > Turn off instrument > Shut down PC.

5. Interpretation of Results

Prior to accepting assay results ensure that the control values are within the acceptance criteria guided by the lot-specific CoA and established by your laboratory. Results for controls and samples will be automatically corrected for the dilution applied during the run.

The assay clinical cutoff for the presence of N protein is 3.0 (arbitrary units). The following table provides guidance in interpreting the Clinical Results Report.

Note: ‘AEB’ refers to ‘average enzymes/bead,’ the digital signal of the assay.

#	Result Displayed	Reason for Result	Interpretation
1	Numerical result <3.0	Undetectable SARS-CoV-2 N protein antigen. AEB within calibration range.	Negative for SARS-CoV-2 N protein antigen.
2	Numerical result ≥3.0	SARS-CoV-2 N protein antigen detected. AEB within calibration range.	Positive for SARS-CoV-2 N protein antigen.
3	“Undetectable”	AEB below calibration range.	Negative for SARS-CoV-2 N protein antigen.
4	“Positive, out of range”	AEB above calibration range.	Positive for SARS-CoV-2 N protein antigen.
5	“Positive, too much fluor”	Fluorescence exceeds optical range.	Positive for SARS-CoV-2 N protein antigen.
6	“NaN” (Not a Number)	Invalid.	Retest

Note: Results of this assay should always be interpreted in conjunction with the patient’s medical history, clinical presentation, and other findings.

Only the qualitative results from the Interpretation table should be reported. Semi-quantitative numerical results have not been clinically or analytically validated and may not correlate with patient disease status, duration of illness or severity of illness. Semi-quantitative results have not been demonstrated to correlate with the success or failure of any therapeutic interventions and should not be used to guide clinical management.

Negative results from patients with should be treated as presumptive and confirmation with a molecular assay, if necessary, for patient management, may be performed.

Simoa HD-X Analyzer® User Guide (EUA Edition)

Software Version 3.1

USER-0105 05

04 Jan 2021

simOaTM



Customer Support

Customer support is available 8 AM to 5 PM, EDT.

[email: techsupport@quanterix.com](mailto:techsupport@quanterix.com)

Contact Us

Quanterix® Corporation

900 Middlesex Turnpike

Billerica, MA 01821

[email: techsupport@quanterix.com](mailto:techsupport@quanterix.com)

Simoa® Semi-Quantitative SARS-CoV-2 IgG Antibody Test

For use under an Emergency Use Authorization (EUA) for In Vitro Diagnostic Use Only.

Prescription Use Only.

The intended use of the Simoa® Semi-Quantitative SARS-CoV-2 IgG Antibody Test and Simoa HD-X Analyzer is for use under the Food and Drug Administration's Emergency Use Authorization (EUA) only. The tests have not been FDA cleared or approved; the tests have been authorized by the FDA under an EUA for use by laboratories certified under CLIA, that meet requirements to perform moderate or high complexity tests. Diagnostic use is authorized while the EUA remains in place.

This product has not been FDA cleared or approved, but has been authorized for emergency use for use by laboratories certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C 263a, that meet requirements to perform moderate or high complexity tests.

This product has been authorized for detecting the presence of IgG antibodies against SARS-CoV-2, not for any other viruses or pathogens.

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug, and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

Simoa® SARS-CoV-2 N Protein Antigen Test

For use under an Emergency Use Authorization (EUA) for In Vitro Diagnostic Use Only.

Prescription Use Only.



The intended use of the Simoa® SARS-CoV-2 N Protein Antigen Test and Simoa HD-X Analyzer is for use under the Food and Drug Administration's Emergency Use Authorization (EUA) only. The tests have not been FDA cleared or approved; the tests have been authorized by the FDA under an EUA for use by laboratories certified under CLIA, that meet requirements to perform moderate or high complexity tests. Diagnostic use is authorized while the EUA remains in place.

This product has not been FDA cleared or approved, but has been authorized for emergency use for use by laboratories certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C 263a, that meet requirements to perform moderate or high complexity tests.

This product has been authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens.

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug, and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

This product is protected by US and Foreign patents and patent filings.

Quanterix Corp. provides this document to its customers with a product purchase to use in the product operation. This document is copyright protected and any reproduction of the whole or any part of this document is strictly prohibited, except with the written authorization of Quanterix Corp.

The contents of this document are subject to change without notice. All technical information in this document is for reference purposes only. System configurations and specifications in this document supersede all previous information received by the purchaser.

Quanterix Corp. makes no representations that this document is complete, accurate, or error-free and assumes no responsibility and will not be liable for any errors, omissions, damage, or loss that might result from any use of this document, even if the information in the document is followed properly.

This document is not part of any sales contract between Quanterix Corp. and a purchaser. This document shall in no way govern or modify any Terms and Conditions of Sale, which Terms and Conditions of Sale shall govern all conflicting information between the two documents.

© 2020 Quanterix Corporation. All rights reserved.



Contents

1	About This Document	1
	Document Purpose	1
	Intended Audience.....	1
	Customer Support.....	1
	Purchasing Supplies	2
	Abbreviations Used in This Document.....	2
	Notices Used in This Document.....	2
2	Safety Information	3
	Read the Safety Notices.....	3
	General Safety Notices	3
	Electrical Safety Notices.....	4
	Chemical/Biological Safety Notices	5
	Laser and Electrical Safety Labels	6
	Safety Certifications.....	7
	Labeling.....	7
	Restriction of Hazardous Substances (RoHS2).....	7
3	Overview of the Simoa® HD-XAnalyzer.....	9
	Introduction	9
	Key Instrument Parts	10
	External View	10
	Reagent Bay	11
	Sample Bay.....	12
	System Resources Drawers.....	13
	System Bay.....	14
	System Fluid Containers	15
	System Waste Containers.....	16
	Overview of Key Instrument Operations	17
4	Quick Start Guide.....	19
5	Operating the Simoa HD-X Analyzer®	22
	Connecting the Simoa HD-X Analyzer to ElectricalPower	23
	Turning the Instrument Computer On and Off.....	24
	Turning the Instrument On and Off	25



Logging In to and Out of the Simoa Software.....	26
Preparing the Instrument for the First Run of the Day.....	28
Shutting Down the Instrument After the Last Run of the Day	29
Operating the Barcode Scanners	29
Using the Handheld Scanner	29
Using the On-Board Scanner.....	30
Opening and Closing the Instrument Bay Drawers.....	31
Inserting and Removing a Reagent Rack	32
Inserting and Removing a Plate Rack.....	33
Inserting and Removing a Tube Rack.....	35
Inserting and Removing an RGP Rack	37
Loading Instrument Consumables	38
Loading Cuvettes.....	38
Loading Disposable Pipettor Tips.....	40
Loading Simoa Discs.....	43
Refilling the System Fluid Containers	48
Changing the Sealing Oil Bag	51
Emptying the Solid Waste Container	53
Emptying the Liquid Waste Containers	55
Monitoring the Instrument's Internal Temperature	58
6 Importing Assay Definitions	59
Understanding Assay Definitions.....	59
Checking Imported Assay Definitions	59
Importing an Assay Definition	61
Configure Assay Settings.....	62
Read-Only	63
Visible Only to Me	63
Archive	64
7 Performing an Assay Run	65
Steps in an Assay Run	67
The Role of Barcodes in an Assay Run	67
Manually Entering a Reagent Barcode	67
Preparing Reagents, Calibrators, Controls, and Samples	68
Use Approved Sample Plates and Sample Tubes.....	69

Preparing Reagents, Calibrators, Controls and Samples	69
Preparing the Instrument for the Run	69
Loading Reagents and RGP into the Instrument.....	71
Using the Barcode Scanners	71
Understanding the Load Reagents Screen.....	72
Loading Reagents	73
Loading RGP	77
Setting Up a Run	80
Understanding the Assign Calibrators Screen	81
Understanding the Assign Samples Screen	83
Understanding the Plate Position Diagram	85
Selecting Processing Order	85
Selecting Plate Wells and Removing Assignments	85
Understanding the Required Sample Volume Field	86
Understanding the Autoname Tool	86
Scanning Sample Barcodes	87
Assigning Calibrators.....	88
Assigning Samples	96
Loading Calibrators/Samples into the Instrument	102
Importing and Exporting a Layout	103
Starting a Run.....	103
Understanding the System Resources Screen	105
Understanding the Solid Resources Screen	106
Monitoring a Run	107
Understanding the In-Process Run Screen for Plates	107
Understanding the In-Process Run Screen for Tubes	108
Checking Instrument Status.....	109
Viewing the Event Log	109
8 Analyzing Run Data.....	111
Data Analysis Overview	111
Working with Run History.....	111
Exporting Results in CSV Format.....	114
Exporting Results in PDF or Microsoft Excel Format.....	115
Analyzing Results in the Data Reduction Tab	115



9	Using Reports, Run History, and the Event Log	119
	The History and Reports Tab	119
	Understanding the Reports Screen.....	120
	Understanding the Run History Screen	121
	Understanding the Event Log Screen	122
	Customizing the Table View in the History & Reports Tab.....	123
	Creating a Report.....	125
	Exporting or Printing a Report	126
	Customizing the Logo for a Report	126
10	Maintaining the Simoa® HD-XAnalyzer.....	127
	Performing Maintenance Tasks	128
	Understanding the Maintenance Screen.....	129
	Determining Whether a Maintenance Task Was Performed	129
	Setting and Viewing Task Details.....	130
	Changing Columns	131
	Performing the Start of Day Task.....	132
	Performing the End of Day Task	133
	Performing the Idle Fluid Prime Task.....	134
	Performing the Database Clean Task.....	135
	Cleaning Exterior Instrument Surfaces	136
	Cleaning the Touchscreen.....	136
	Cleaning the System Resources Drawers.....	136
	Cleaning the Reagent and Sample Bays.....	137
	Cleaning the Liquid Waste Containers.....	137
	Cleaning the System Fluid Containers and Tube Rinsing.....	138
	Maintaining Idle Instruments	143
11	Configuring Software Settings.....	145
	Software Settings.....	145
	Changing System Settings.....	146
	Selecting the Sample Loading Mode	146
	Setting an Automatic Lock Out Time	147
	Enabling Password Policies	147
	Controlling Tab Switching when Bay Doors are Opened	148
	Setting a Password Expiration Time	148

Setting Default Decimal Places	149
Specifying the Customer Defined System Name	149
Locking an Account due to Incorrect Login Attempts	150
Specifying User Settings.....	150
Specifying Account Management Settings	152
Managing Users	153
Managing Roles and Features.....	155
Managing Report Settings	156
Appendix A – Configure Columns Definitions.....	159
Appendix B – Instrument Specifications	167
Physical Specifications	167
Installation & Service Dimension Requirements	167
Power Requirements	168
Instrument	168
Computer and UPS.....	168
Environmental Operating Conditions	168
Noise Emission	168
Lasers	169
Laser Barcode Scanner.....	169
Laser Light Barrier	169
Packaging	169
Environmental Conditions	170
Instrument Installation Dimensions	171

1 About This Document

Document Purpose	1
Intended Audience	1
Customer Support	1
Purchasing Supplies.....	2
Abbreviations Used in This Document.....	2
Notices Used in This Document	2

Document Purpose

This document provides information about how to perform Simoa® assays with the Simoa HD-X Analyzer® (termed “the instrument” subsequently). After reading this document, you will know how to operate the instrument, perform assays, maintain the instrument, and obtain help for instrument or assay problems.

Intended Audience

This guide is for anyone who operates or maintains the instrument. Users should be familiar with standard laboratory practices and with computers running the Windows 10 operating system.

Customer Support

Customer support is available 8 AM to 5 PM EST, Monday through Friday.

Email: techsupport@quanterix.com

When you contact Quanterix®, you may be asked to use the Customer Support Tool to generate a support package that will help troubleshoot the problem. See *Quanterix Customer Support Tool User Guide* for instructions.

Log in to the Quanterix portal to view training videos at
<http://portal.quanterix.com>.

Purchasing Supplies

For information about purchasing reagent kits, consumables, and instruments, contact Sales at Quanterix:

Quanterix Corporation
900 Middlesex Turnpike
Billerica, MA 01821
Phone: 617-301-9400
Email: sales@quanterix.com

Abbreviations Used in This Document

CCD	charge-coupled device (camera imaging system)
ELISA	enzyme-linked immunosorbent assay
GUI	graphical user interface
SBG	Streptavidin-beta-galactosidase
RGP	resorufin-D-galactopyranoside

Notices Used in This Document



Note! This type of notice highlights important information about procedures or provides tips.



CAUTION! This type of notice provides information about minor potential hazards to the instrument or assay performance. Failure to comply with the notice may result in inaccurate results, instrument maintenance problems, or instrument failure.



WARNING! This type of notice provides information about situations that pose a danger to the instrument or to operators. Failure to observe the warnings in this notice may result in instrument damage or a safety hazard for operators.

2 Safety Information

Read the Safety Notices	3
General Safety Notices.....	3
Electrical Safety Notices	4
Chemical/Biological Safety Notices.....	5
Laser and Electrical Safety Labels.....	6
Safety Certifications	7

Read the Safety Notices

Before you operate the Simoa HD-X Analyzer, familiarize yourself with the safety and regulatory notices in this chapter. If improperly operated, the instrument can pose electrical or mechanical safety hazards to operators. You may also encounter chemical or biological hazards, depending on the assays that you perform with the instrument.

Do not perform any instrument procedures that are not described in this guide without the guidance of Quanterix® Customer Support.

The Simoa HD-X Analyzer must be installed by trained service professionals only. Do not interfere with system electrical grounding or move the instrument after installation without the assistance of Quanterix Customer Support.

General Safety Notices



CAUTION! As the Simoa HD-X Analyzer is designed for measuring biological samples, universal laboratory precautions should be used at all times.



WARNING! Do not install any software on the Simoa HD-X Analyzer computer unless instructed to do so by Quanterix Customer Support. Installing third-party software can result in system failure, job loss, and voiding of your warranty. Contact Technical Support for further details and assistance.



CAUTION! Do not move the instrument after the initial installation.



CAUTION! Do not remove the cover of the instrument.

Label located on pipette rail near cuvette loader.



CAUTION – Mechanical hazard! Do not insert your hands into the reagent bay, the sample bay, or the system bay if the instrument is turned on. Doing so may expose you to risk of personal injury. The instrument contains robotic pipettors that access these areas at times, and they can injure your hands.

Label is located on sample bay door, reagent bay door, drawers, waste container drawer (left and right side) and cuvette loader.



CAUTION! Hot surface!

Label is located on LSI unit.



Overheating!

Improper placement of the instrument may cause fire or serious instrument damage in case of overheating.

- Do not block or cover ventilation slots.
- The air shall be able to circulate.

Electrical Safety Notices



CAUTION! The appliance inlet is intended to be used as a disconnecting device in case of an emergency. Set up the instrument so that the appliance inlet is always easily accessible.



Improper connection of mains supply.

Improper connection of the instrument and the peripheral devices to the mains supply can cause serious personal injury with potentially deadly consequences and material damage (e.g. fire).

- Only use grounded connection and extension cables with sufficient capacity (voltage and current) to connect the instrument and any peripheral devices to the mains power supply.
- Never remove ground connections.
- Grounding of the instrument and its peripheral devices to the same protective earth potential shall be ensured.
- The use of a multi-outlet power strip is not allowed!
- Only use power cables that fulfill the minimum requirements for this instrument.

Chemical/Biological Safety Notices



CAUTION! Biohazard materials may be used on this instrument. Proper personal protective equipment (PPE) should be worn at all times to meet the level of biohazard materials used.




Risk of infection! Monthly maintenance and waste removal procedures put the user in close contact with biohazard materials. Proper personal protective equipment (PPE) should be worn at all times to meet the level of biohazard materials used. Label is located on drawers.

- Strictly follow the local and national provisions, legislation and laboratory regulations.
- Use appropriate gloves!
- Use an appropriate lab coat!
- Use an appropriate eye protection (e.g. protective glasses)!
- Avoid contact between skin/mucous membrane and samples/test reagents or parts of the instrument.
- Clean, disinfect and decontaminate the instrument immediately if potentially infectious material has been spilled.
- Do not use broken or chipped tubes or bottles.
- Observe the instructions in the package inserts for correct use of reagents.
- Observe the legal regulations for the handling of infectious material.

-
- Never use bio-hazardous liquids for testing the instrument!
 - The instrument shall be cleaned, disinfected and decontaminated before servicing!
-

Laser and Electrical Safety Labels

The following regulatory labels are affixed to the Simoa HD-X Analyzer.

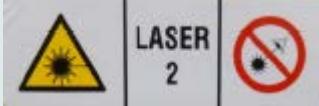Label Image	Description of Label Information
	<p>Laser safety label. The label is located on the side of the reagent bay (right) and the sample bay (left), adjacent to the laser inside of the bay.</p>
	<p>Laser information label. Provides the user with information on the class, wavelength, output power, and regulatory compliance status of the laser. The label is for the barcode readers next to the reagent bay and sample bay.</p>
	<p>Biohazard label. The label is located in five places:</p> <ul style="list-style-type: none"> • The reagent bay door, in front of the reagent bay lanes. • The sample bay door, in front of the sample bay lanes. • The front of liquid waste bottles 1 and 2. • The solid waste intermediate trap door to the waste bin. • The system bay waste drawer, inside rear surface.
 <p>Vor Öffnen des Gehäuses Netzstecker ziehen. Before opening disconnect mains. Avant d'ouvrir l'appareil retirez la fiche mâle.</p>	<p>Electrical hazard label. The label is located above and adjacent to the electrical power switch, on the side of the instrument.</p>

Label Image	Description of Label Information
	<p>Heavy weight. The label is located on the inside of the reagent and sample bay doors. Do not place heavy objects on the doors. When the system bay drawers are open, do not step on them or apply additional weight.</p>


Safety Certifications

This product complies with the EMC standards listed below for class A devices and is suitable for use in commercial environments. In a domestic environment, this product may cause radio interferences, in which case the user may be required to take adequate measures.

- EN ISO 14971:2012 (medical devices – Application of risk management to med devices)
- EN 61010-1 and EN 61326-1 (International) Safety of operation in a laboratory environment

Labeling

The following label is affixed to the instrument:

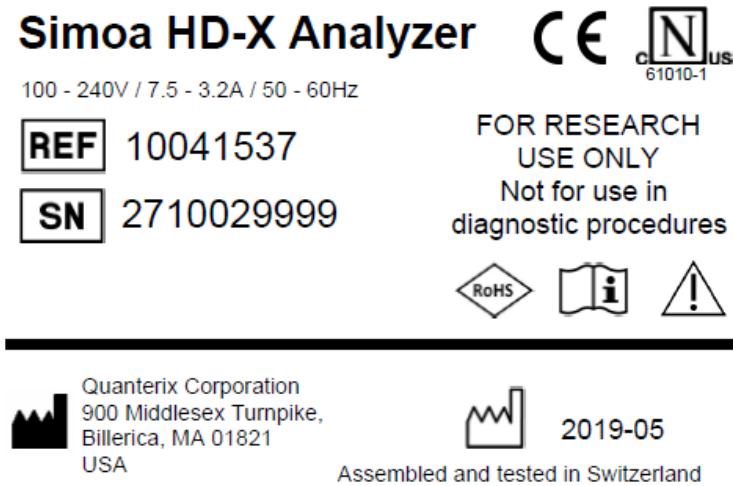

Restriction of Hazardous Substances (RoHS2)

The instrument system complies with the RoHS2 Directive:



- Directive 2011/65/EU – Restriction of Hazardous Substances in Electrical and Electronic Equipment

3 Overview of the Simoa® HD-X Analyzer

Introduction9
Key Instrument Parts.....	10
External View	10
Reagent Bay.....	11
Sample Bay.....	12
System Resources Drawers	13
System Bay	14
System Fluid Containers.....	15
System Waste Containers	16
Overview of Key Instrument Operations	17

Introduction

The Simoa HD-X Analyzer is a fully automated instrument that performs ELISA immunoassays. At the heart of the instrument is a precision-engineered, disposable disc with twenty-four ~238,000 femtoliter-sized reaction chambers, termed single molecule arrays (Simoa®), that can isolate and detect single enzyme-labeled protein molecules. Because the array volumes are approximately 2 billion times smaller than those required for conventional ELISA, a rapid buildup of fluorescent product is generated if a labeled protein is present. With diffusion limited, this high local concentration of product can be readily observed.

Only a single molecule is needed to reach the detection limit.

Quanterix provides a wide range of assay kits for the Simoa HD-X Analyzer, containing ready-to-load reagents.


Performing a Simoa assay of any type is simple:

- Load assay reagents into racks in the reagent bay.
- Load plates or tubes of samples and a fluorescent substrate (RGP) into racks in the sample bay.
- Initiate the assay with the instrument's touch-screen software.

- View and analyze assay results in the instrument software or export them to analysis packages or a laboratory information management system (LIMS).

Key Instrument Parts

External View

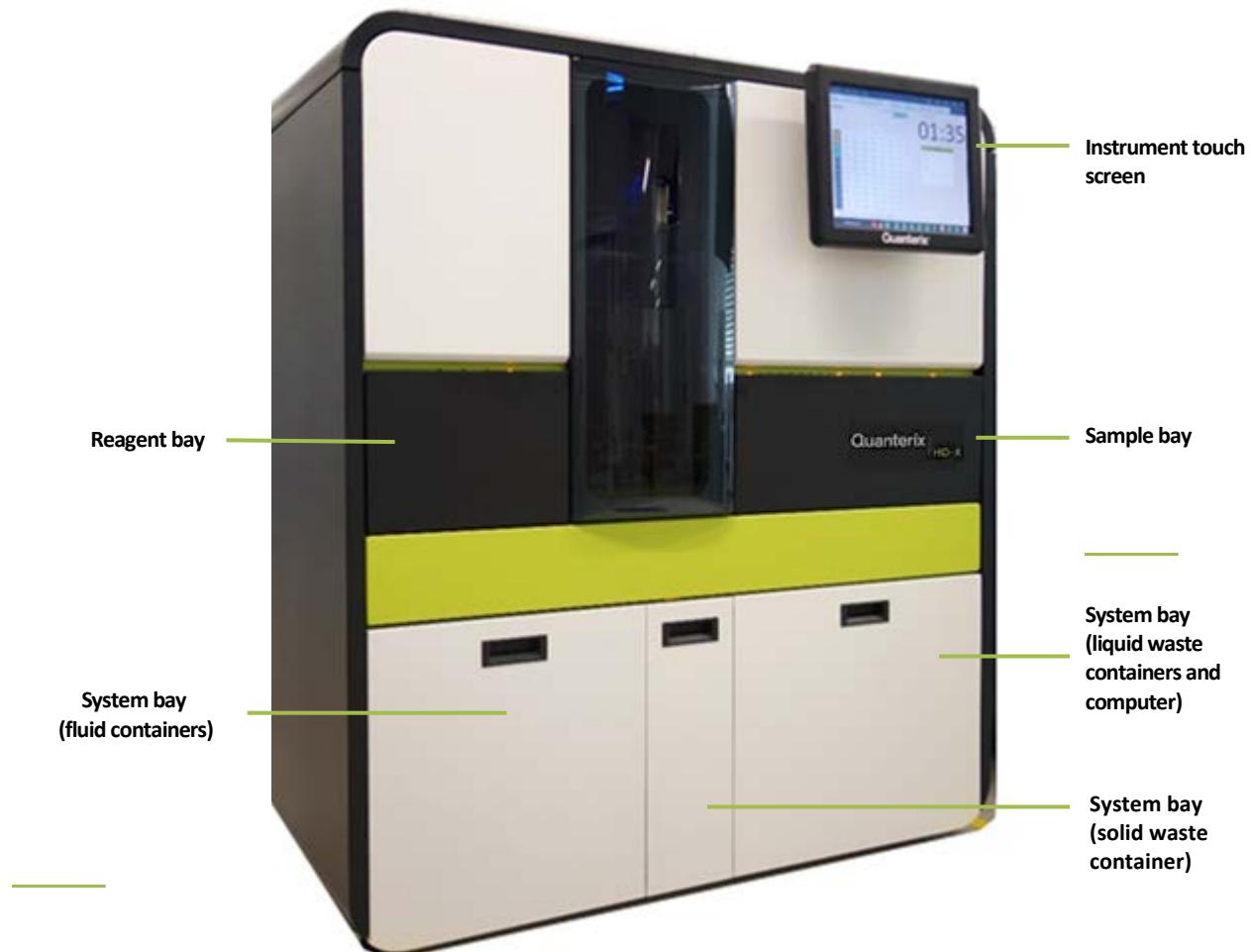

Reagent Bay

Reagent bay showing four reagent rack lanes



Lane 1 2 3 4

Loading a reagent rack
Into reagent bay lane 4




Sample Bay

Sample bay, showing eight lanes for plate racks or tube racks and two lanes for RGP racks




Loaded plate rack being loaded into sample bay lanes 1-4 and RGP rack loaded into RGP lane 1

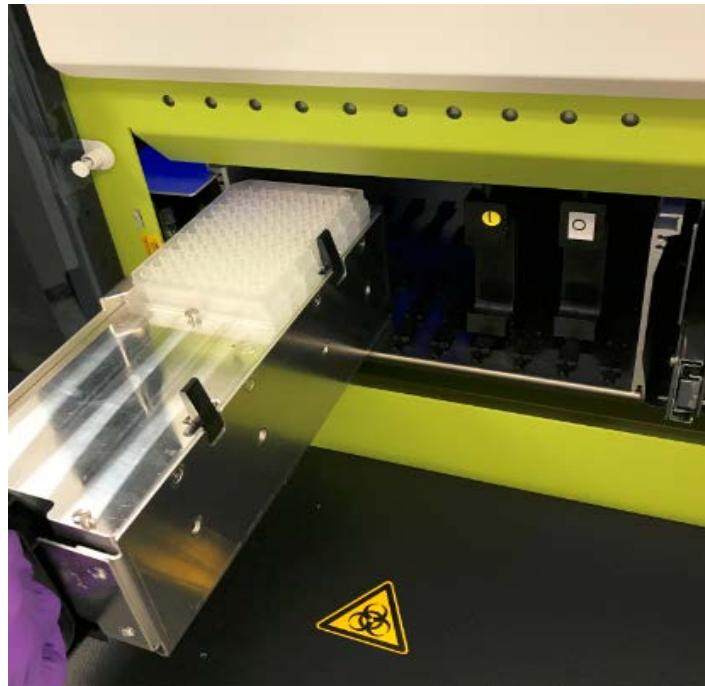

System Resources Drawers

Stack of Simoa discs on mounting pole

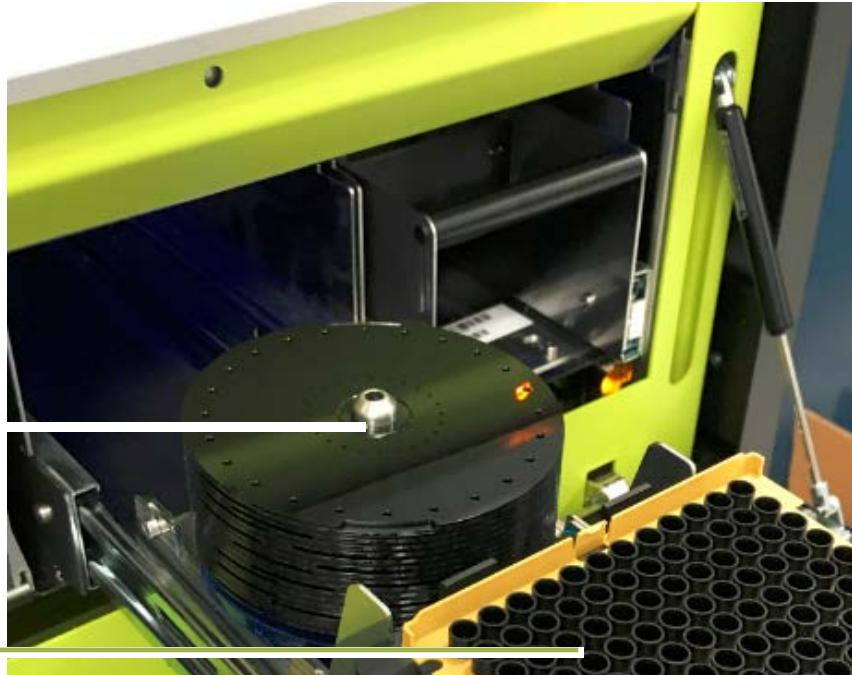

Disposable pipettor tip rack

System Bay



System fluid containers, page 15

Solid waste container, page 16

Liquid waste containers, page 16



Labeling for System fluid (left) and Liquid waste containers (right)


System Fluid Containers

The system uses three system fluids:

- Deionized water (system fluid and needle wash)
- Wash buffer 1 (main assay wash buffer)
- Wash buffer 2 (bead wash pre-loading)

The deionized water and wash buffer 1 each have a secondary container that feeds into a primary container. The instrument draws fluid from the primary containers, which in turn replenish themselves passively from their secondary containers. The container caps connect to tubing that delivers system fluids where needed in the instrument. The secondary containers connect to the primary containers via a plug.

Each secondary system fluid container feeds into a primary container via a short connector. The instrument monitors fluid levels in the containers. The software notifies operators when fluids need replenishing via color-coded icons.

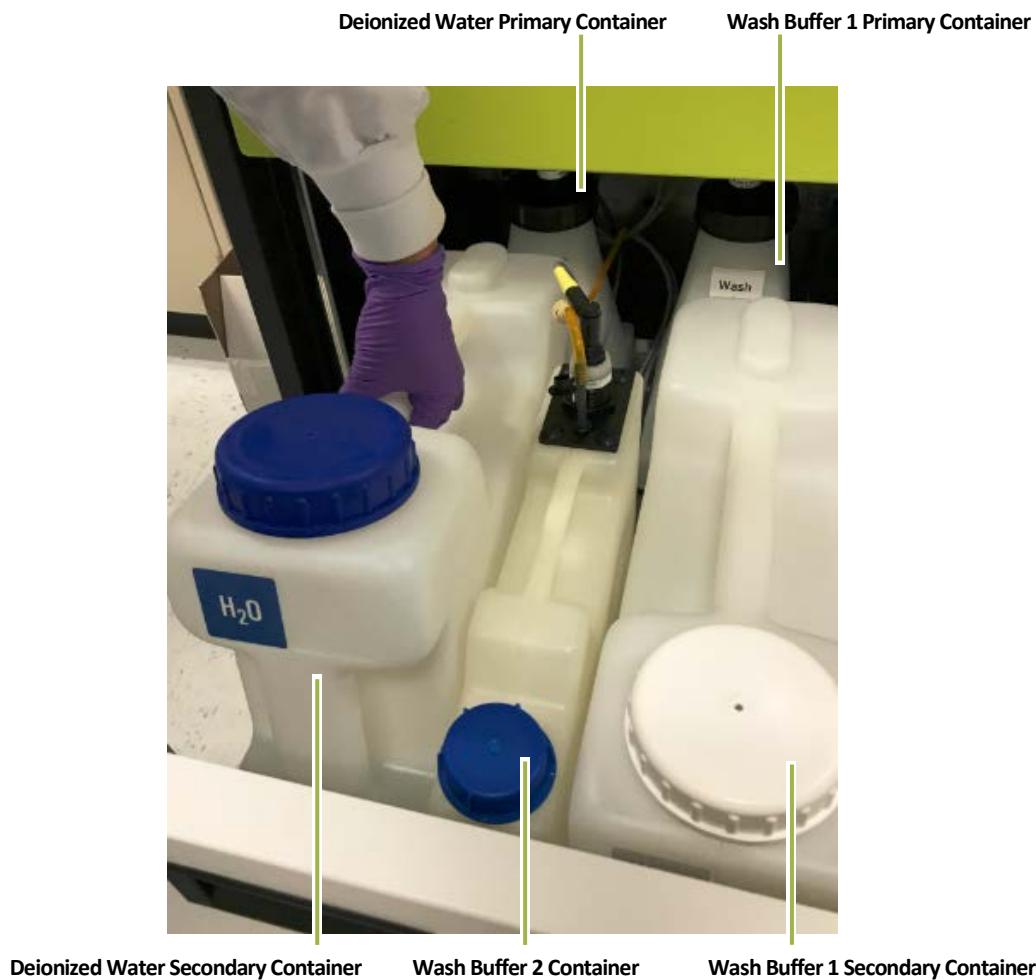

System Waste Containers

The instrument discards used cuvettes, tips, and Simoa Discs into the solid waste container and the effluent from sample washing and incubation into the liquid waste container.

The caps on the liquid waste containers hold liquid-level sensors that allow the instrument to determine when the containers are full and need to be emptied.




Solid waste container



Liquid waste containers

Overview of Key Instrument Operations

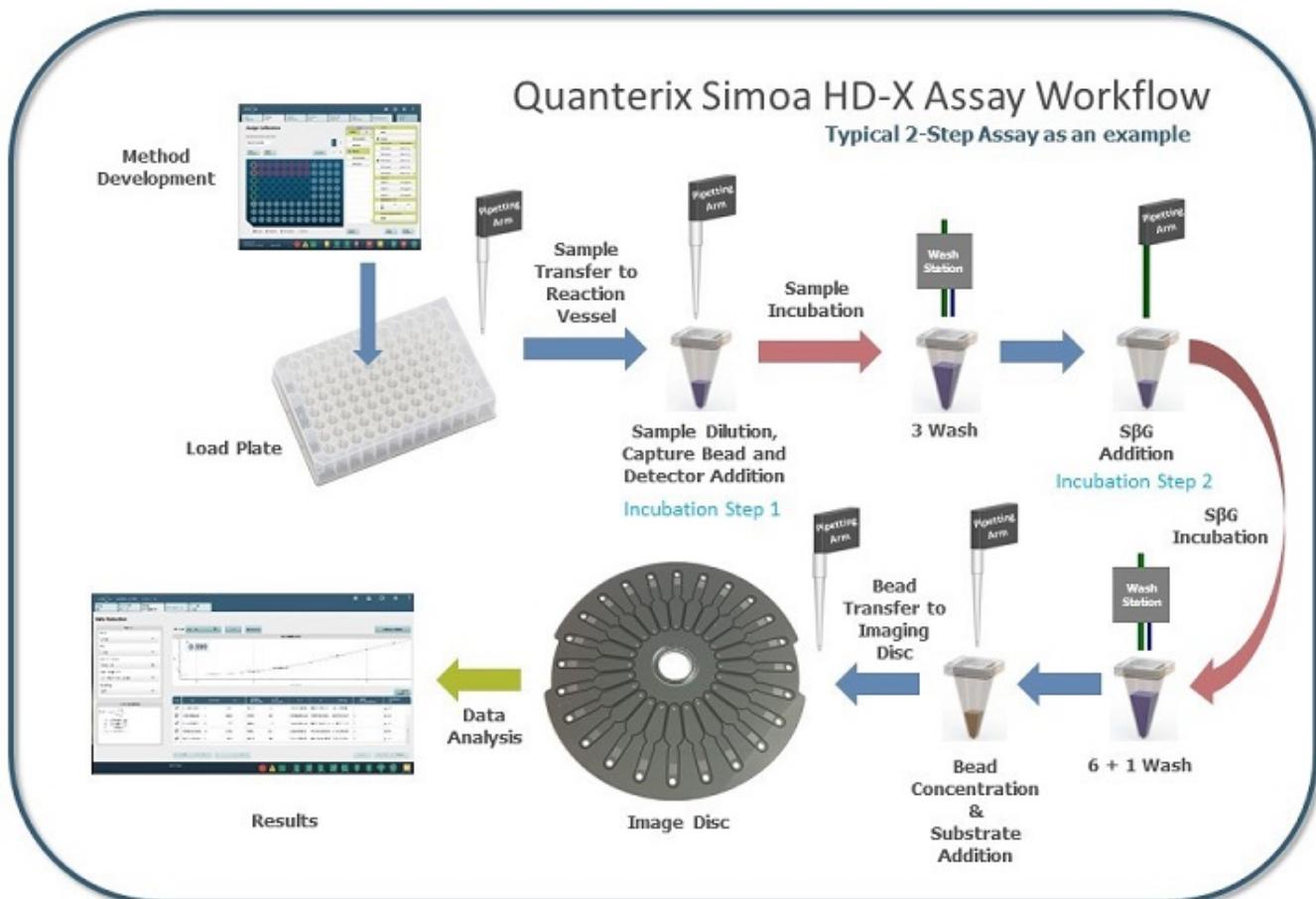

- **Method Development** – Set up your assay using the HD-X touch screen.
- **Load Plate** – Load reagents, samples, and consumables.
- **Sample Transfer to Reaction Vessel** – Robotic pipettors deliver samples into the cuvettes.
- **Sample Dilution, Capture Bead and Detector Addition** – Samples are diluted; capture beads and detector reagents are added to the cuvette.
- **Sample Incubation** – The system implements a one-, two-, or three-step incubation series by moving the cuvette within the washer/incubator ring per the details of the user-specified assay protocol.
- **3 Washes/6+1 Washes** – The washer/incubator ring moves the cuvettes through a series of pelleting, mixing, and washing, according to the protocol. All washes, except the final wash, are performed using Wash Buffer 1. Wash Buffer 2 is used for the final wash preceding substrate addition.

- **SBG Addition / SBG Incubation** – SBG is added to the cuvette and the washer/incubator ring implements the protocol.
- **Bead Concentration & Substrate Addition** – Beads are pelleted, wash buffer is removed and beads are re-suspended via mixing with enzyme substrate, RGP.
- **Bead Transfer to Imaging Disc** – The disposable tip pipettor transfers the mixture of beads and RGP into the sample inlet of a microarray on the Simoa Disc. The instrument loads and seals individual beads into microwells within the array. Loose beads that have not settled into microwells are flushed away.
- **Image Disc** – The instrument moves the disc to the imaging station, where a camera images the sealed wells, capturing the growth of fluorescent signal generated by enzyme labeled beads.
- **Data Analysis/Results** – The Simoa software analyzes the image, determines the average enzymes per bead (AEB), and generates a curve.

4 Quick Start Guide

This quick start guide walks you through performing an assay run on the Simoa® HD-X Analyzer. Be sure to read the entire user guide before you perform your first run.

After that, you can use this guide as a ready reference and to locate explanations of instrument procedures.

Log in to the Quanterix customer portal to view training videos for additional detail at <http://portal.quanterix.com/>.



Note! The Instrument and the computer must be power-cycled once every 24 hours to permit internal systems to reinitialize.

1. Turn on the Simoa HD-X Analyzer computer.	Chapter 5, page 24
2. If the instrument is not running, turn it on. If the instrument has been running continuously overnight, turn it off, wait 15 seconds, and turn it on again.	Chapter 5, page 25
3. Log in to the Simoa software.	Chapter 5, page 26
4. If performing the first run of the day: <ul style="list-style-type: none">Perform the End of Day task if it was not performed the previous night.Perform the Start of Day task.	Chapter 10, page 133 Chapter 10, page 132
5. Prepare reagents, samples, and calibrators.	Chapter 7, page 68
6. If needed: <ul style="list-style-type: none">Empty the solid and liquid waste containers.Refill system fluid containers (deionized water, system wash buffer 1 and 2).Load cuvettes, disposable pipettor tips, and Simoa discs.	Chapter 5, page 53 Chapter 5, page 48 Chapter 5, page 38

7.	Mix the bead reagent according to the kit instructions. Re-mix the bead reagent if it sits for >5 minutes before being loaded.	Chapter 7, page 73
8.	Set up the software for sample plates or sample tubes.	Chapter 11, page 146
9.	Load bead, detector, and SBG reagents into the instrument. <ul style="list-style-type: none"> Place bead, detector, and SBG in a reagent rack. 	Chapter 7, page 73
	<ul style="list-style-type: none"> If using the handheld barcode scanner: Identify the rack position of each reagent bottle and the lane into which you plan to load them. Scan the bottle barcodes and insert the reagent rack. 	
	<ul style="list-style-type: none"> If using the on-board scanner: Enable the on-board scanner in the Load Reagents screen and select the lane into which you will load the reagent rack. Insert the rack. 	
10.	Load RGP reagent into the instrument. <ul style="list-style-type: none"> Place RGP reagent in the RGP rack. 	Chapter 7, page 77
	<ul style="list-style-type: none"> If using the handheld barcode scanner: Identify the rack position of the RGP reagent and the lane into which you plan to load it. Scan the RGP bottle barcode and insert the RGP rack. 	
	<ul style="list-style-type: none"> If using the on-board scanner: Enable the on-board scanner and select the lane into which you will load the RGP rack. Insert the rack. 	
11.	Load calibrators and samples: <ul style="list-style-type: none"> Assign calibrators in the software. Assign samples in the software. 	Chapter 7, page 88
	<ul style="list-style-type: none"> Insert calibrator and sample plates or tubes. 	Chapter 7, page 96
		Chapter 7, page 102

12. Start the run.	Chapter 7, page 103
13. Monitor the assay in the Current Run tab.	Chapter 7, page 107

5 Operating the Simoa HD-X Analyzer®

Connecting the Simoa HD-X Analyzer to Electrical Power	23
Turning the Instrument Computer On and Off	24
Turning the Instrument On and Off	25
Logging In to and Out of the Simoa Software	26
Preparing the Instrument for the First Run of the Day	28
Shutting Down the Instrument After the Last Run of the Day	29
Operating the Barcode Scanners	29
Using the Handheld Scanner	29
Using the On-Board Scanner	30
Opening and Closing the Instrument Bay Drawers	31
Inserting and Removing a Reagent Rack	32
Inserting and Removing a Plate Rack	33
Inserting and Removing a Tube Rack	35
Inserting and Removing an RGP Rack	37
Loading Instrument Consumables	38
Loading Cuvettes	38
Loading Disposable Pipettor Tips	40
Loading Simoa Discs	43
Refilling the System Fluid Containers	48
Changing the Sealing Oil Bag	51
Emptying the Solid Waste Container	53
Emptying the Liquid Waste Containers	55
Monitoring the Instrument's Internal Temperature	58

Connecting the Simoa HD-X Analyzer to Electrical Power

The Simoa HD-X Analyzer must be connected to a dedicated electrical outlet that conforms to the instrument specifications for electrical power (see *Simoa HD-X Analyzer Site Requirements*). Your instrument will be equipped with the correct type of electrical plug for your country.



CAUTION! If you must use an extension cord to connect the instrument to its dedicated electrical outlet, use only a cord with a three-prong connector on both ends. Using a two-prong extension cord defeats the grounding circuit of the outlet, creating a hazard of shock.



CAUTION! Quanterix recommends that the dedicated outlet that you use for the instrument be equipped with a ground fault circuit interrupter.



Note! Do not disconnect the instrument from its electrical power source when the instrument is running an assay. Doing this will cause the instrument to discard the Simoa Disc and the cuvettes for the assay when you reconnect it to power.

Turning the Instrument Computer On and Off

The instrument computer is inside of the system bay. Complete all assay operations and log out of the Simoa software before you turn off the computer.

❖ To turn the computer on

Press the on/off power switch at the top of the computer, wait three minutes, power up the instrument (see “Turning the Instrument On and Off” on page 25), and open the HD-X software (see “Logging In to and Out of the Simoa Software” on page 26).

❖ To turn the computer off

Close the software and shut down the computer through Windows Shut down.



Turning the Instrument On and Off

❖ To turn the instrument on


Move the power switch to the on position: Push down the end of the switch with the dash (–). The instrument initializes itself over the next three minutes.

❖ To turn the instrument off

Move the power switch to the off position: Push down the end of the switch with the circle (O).



Note! Do not turn off the instrument during an assay run. Doing this will cause the instrument to discard the Simoa Disc and the cuvettes for the assay when you reconnect it to power.

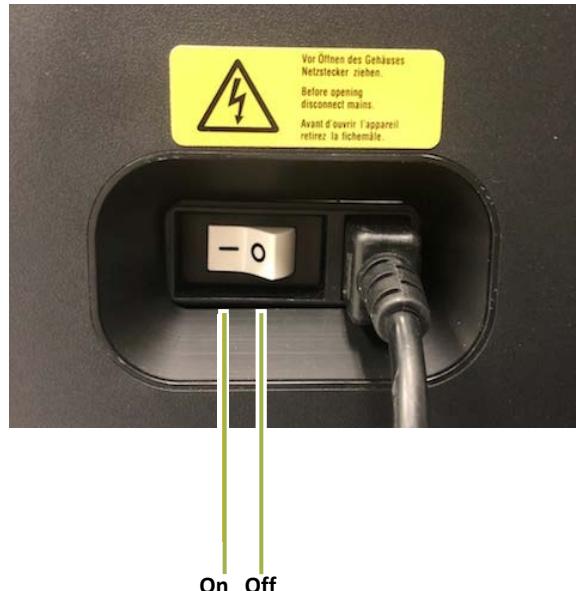

Note! You must turn the instrument on and off once a day to power cycle instrument systems. During power cycling, the instrument performs a series of hardware and software initialization processes. Wait 15 seconds after you power down before powering up again.

If the instrument has not been power cycled in over 24 hours, the software will prompt you to turn it off with the following message:




Select “Yes” to close the software and shut down the instrument computer. You must then manually turn off the instrument by moving the power switch to the off position.

Logging In and Out of the Simoa Software

To log in to the software, you need a user name and a password. Only a user with administrative privileges can add users to the system. For information on managing users, see “Managing Users” on page 153.

❖ To log in to the software

- 1 Double-touch the Simoa software icon on the instrument desktop.

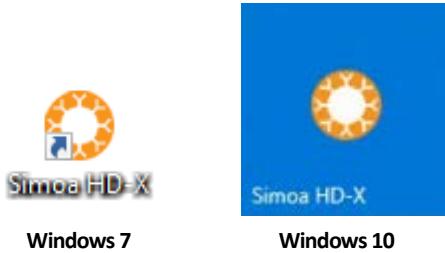

- 2 Enter your user name and password > touch **Enter**.
- 3 Touch **Accept** in the Terms of Use dialog. Depending on how your system is configured, the first time you open the software, you may see the startup screen.



If you check the “Do not show this screen next time” checkbox, the next time you open the software, the tabbed user interface appears.



- To minimize the software GUI, double-touch the letter O in simoa.




❖ To Log Out of the software

- Touch the Logout icon at the top right of the software screen.



- 2 Click Logout to log out the current user (without closing the software) or Shutdown to log out and exit the software.
Click the lock to lock the screen.
- 3 Perform maintenance tasks, exit, or cancel the log out.

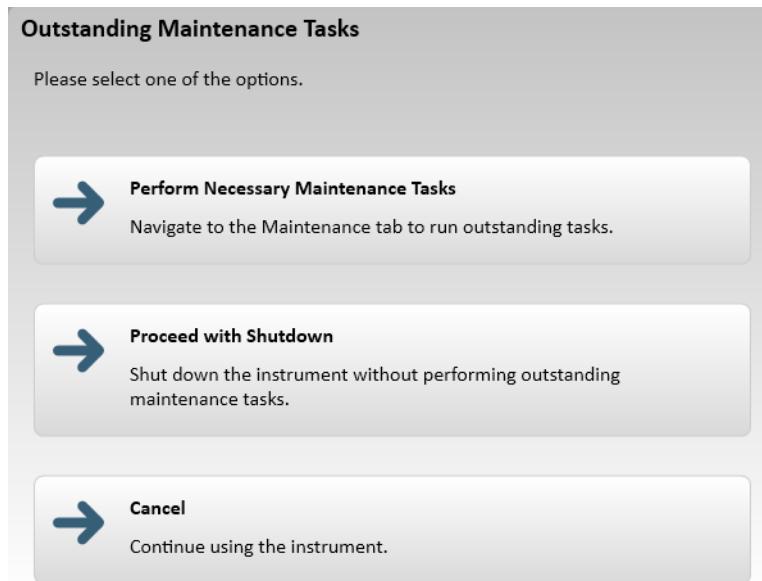

Note! Do not shut down the software during an assay run or while data analysis is in progress, as this will result in run failure or the cancellation of data analysis.

Preparing the Instrument for the First Run of the Day

Before the first run of the day:

- Power cycle the instrument if it was not turned off the preceding night by turning the instrument off, waiting for 15 seconds, and turning it on again (see page 25).
- If the End of Day task was not performed the preceding night, perform it now (see page 133).
- Run the Start of Day task (see page 132).
- Refill the system fluid containers if necessary (see page 48).
- Empty the waste containers if necessary (see page 53).

Shutting Down the Instrument After the Last Run of the Day

After the last run of the day:

- Remove any leftover reagents and samples from the instrument. Recap the reagents and store at 2–8° C or as instructed in the assay kit instructions.
- Empty the solid and liquid system waste containers (see page 53).
- Perform the End of Day task (see page 133).

Operating the Barcode Scanners

You have the option to use the handheld or on-board barcode scanners to scan reagent label barcodes.

Using the Handheld Scanner

The handheld scanner is used to scan reagent, RGP, and Simoa disc barcodes.

The barcode scanner is connected to the instrument computer by default. If you need to connect it, plug it into one of the USB ports on the front of the instrument computer, which is in the system bay. The scanner beeps once when it becomes active.

When not in use, store the scanner in the holster that is mounted on the right side of the instrument.

To operate the scanner, point it at the barcode and press the scan button. The scanner transmits the barcode data to the computer, which displays it in the software.



Using the On-Board Scanner

You can use the on-board scanner as an alternative to the handheld device for scanning reagent, RGP, and sample barcodes. It cannot be used to scan Simoa discs.

The on-board scanner is located on the inside of the reagent and sample bays and reads the barcoded reagents, RGP and samples as you insert the racks into the instrument. To use the on-board scanner, you must enable it by turning on the “Enable On-Board Scanner” toggle in the Load Reagents or Setup Run tab. For more information about using the on-board scanner, see Loading Reagents and RGP into the Instrument on page 71.





Note! The on-board barcode scanner will only work if the red scanner light is visible on the inside of the reagent or sample bay. Make sure that the red light is visible before you load barcoded items into the instrument.



Red scanner light
in sample bay

Opening and Closing the Instrument Bay Drawers

Pull the door/drawers to open the reagent, sample, and system bays. Keep all instrument drawers closed except when you need to open them for access into the instrument or the system bay.

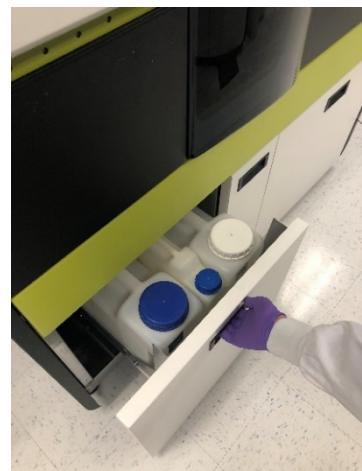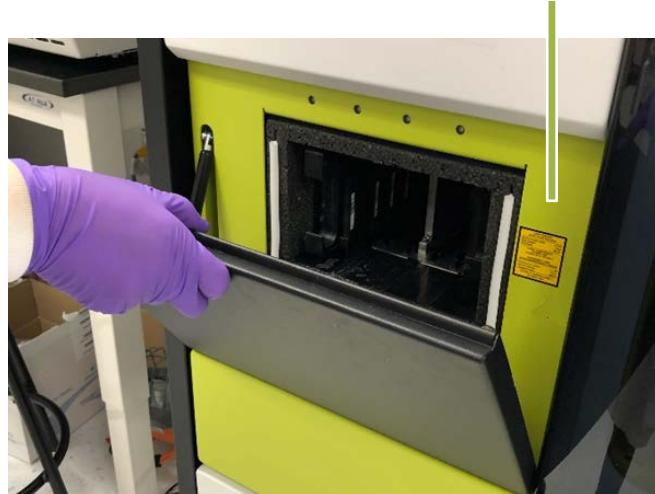

Inserting and Removing a Reagent Rack

Each of the four lanes in the reagent bay can hold one reagent rack. Each reagent rack can hold six 15-mL bottles of assay reagents, for a maximum of 24 reagent bottles.


The first three positions in the reagent racks oscillate during runs and are intended to hold reagents that must be held in suspension (for example, capture bead reagent).



❖ To insert and remove a reagent rack

To insert a rack into a reagent bay lane, line it up between the lane guides and push the rack into the instrument until the catch engages.

To remove a rack, grasp the handle, push the rack in slightly to release the catch, and pull the rack straight out.

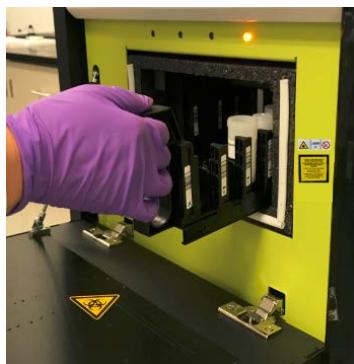

Inserting and Removing a Plate Rack

Lanes 1 to 8 in the sample bay can hold two plate racks. Each plate rack can hold two 96-well plates, for a maximum of four loaded plates. Plate rack 1 goes in lanes 1–4, and plate rack 2 in lanes 5–8.

If you are loading only one plate, place it in position 1 in plate rack 1.



CAUTION! Use only Quanterix-supplied sample plates. These plates have been verified to be nonreactive with Simoa assay reagents and to be correctly sized to work with the instrument pipettors. Quanterix cannot guarantee system performance if you use any other type of plate. For information on ordering plates, see “Purchasing Supplies” on page 2.


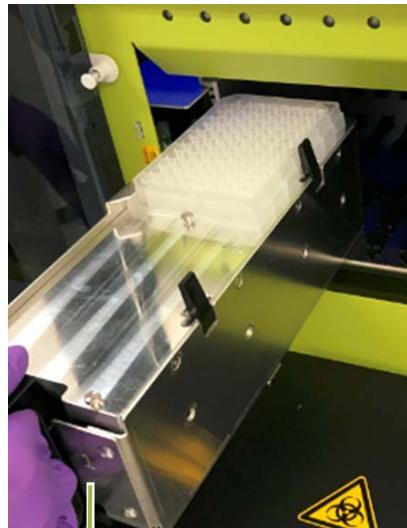


Plate rack number engraved here

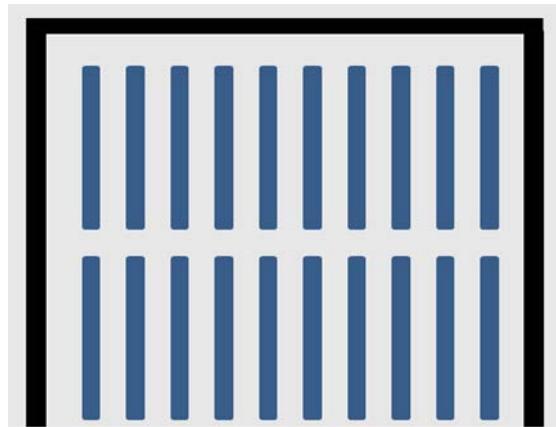


1 2 3 4 5 6 7 8 RGP

Plate rack 1 Plate rack 2

❖ To load a plate onto a rack

- 1 Slide the plate onto the rack.
- 2 Ensure the lip of the plate sits below the head of both positioning screws.
- 3 Lift the tabs to secure the plate.

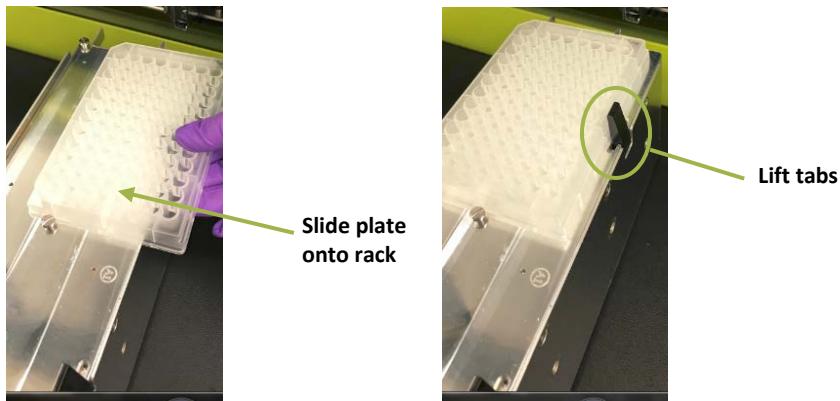


❖ To insert and remove a plate rack

The plate racks are numbered and must be inserted into the instrument in a specific position. Plate rack 1 is inserted in lanes 1 to 4, and plate rack 2 in lanes 5 to 8.

To insert a plate rack, line it up over four lane guide rails in the sample bay and push the rack into the instrument until the catch engages.

To remove a plate rack, grasp the handle, push the rack in slightly to release the catch, and pull the rack straight out.


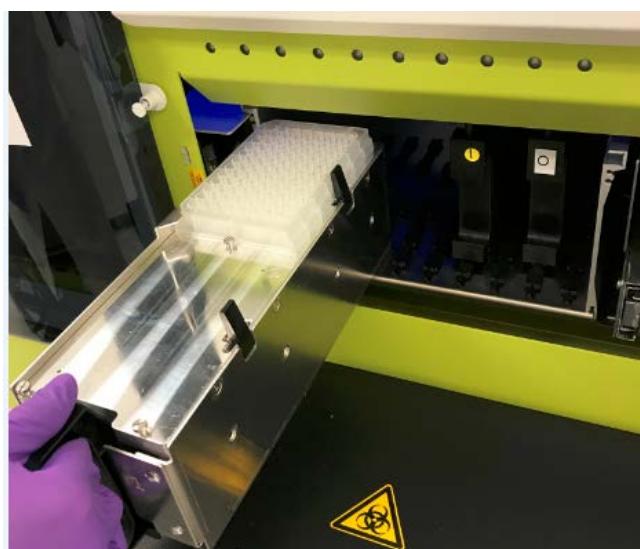


Inserting and Removing a Tube Rack

Lanes 1 to 8 in the sample bay can hold tube racks.

One tube rack can hold 12 tubes, for a maximum of 96 tubes in the sample bay. The tube racks accommodate tubes up to 100 mm tall and 10 to 16 mm in diameter.

1 tube rack per lane

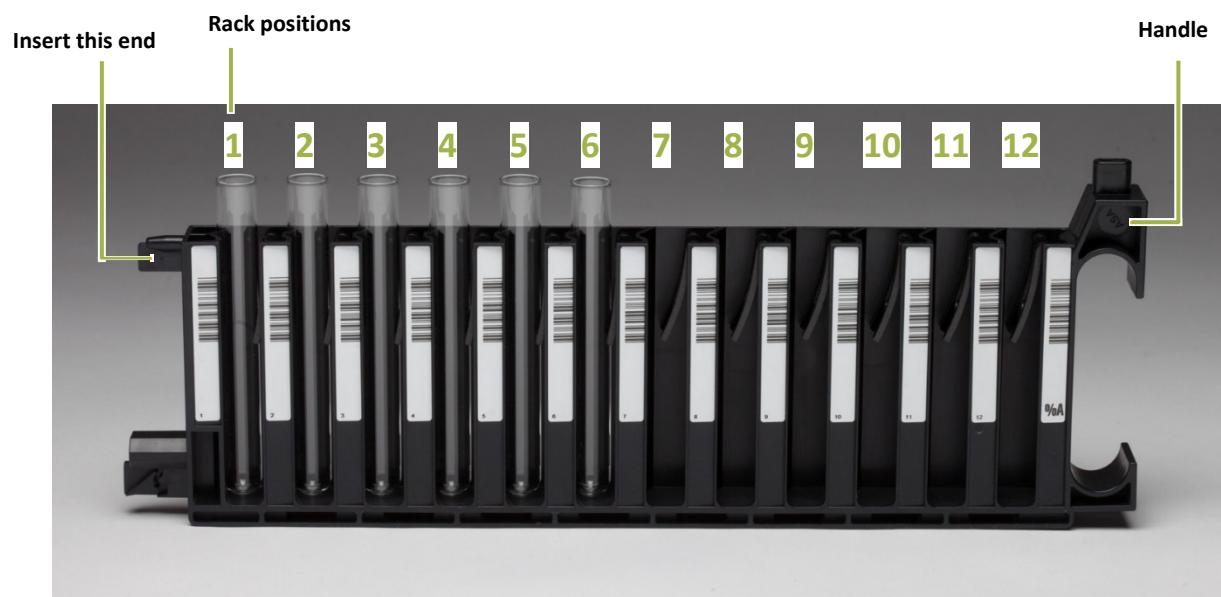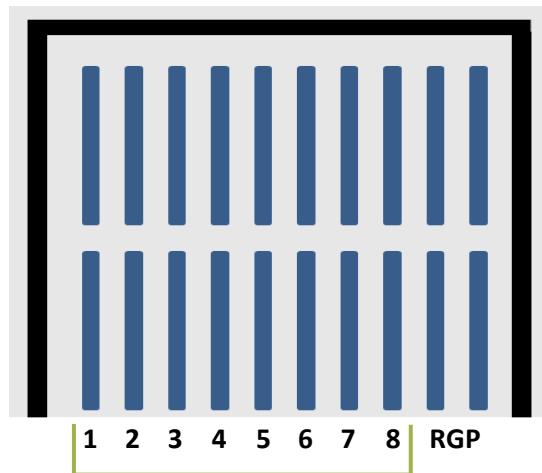

❖ **To insert and remove a tube rack**

To insert a tube rack, center the rack over a lane guide rail in the sample bay and push the rack into the instrument until the catch engages.


To remove a tube rack, grasp the handle, push the rack in slightly to release the catch, and pull the rack straightout.




Inserting and Removing an RGP Rack

The RGP racks are similar to tube racks. They have been customized to hold the short RGP reagent vials.

The RGP racks must be inserted into lanes RGP 1 and RGP 2 only (also termed the “substrate lanes”).

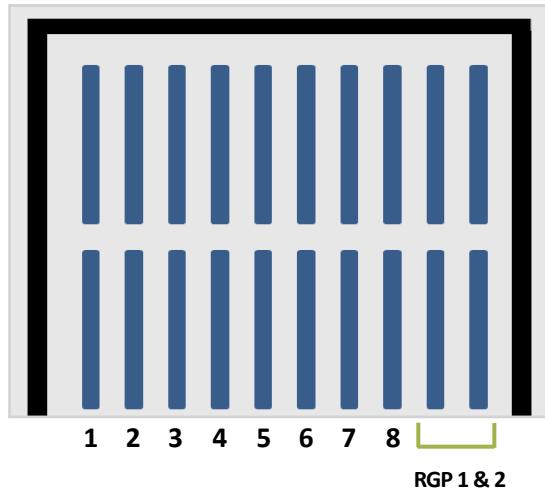

RGP racks are inserted and removed in the same manner as the tube racks (see “Inserting and Removing a Tube Rack” on page 35).

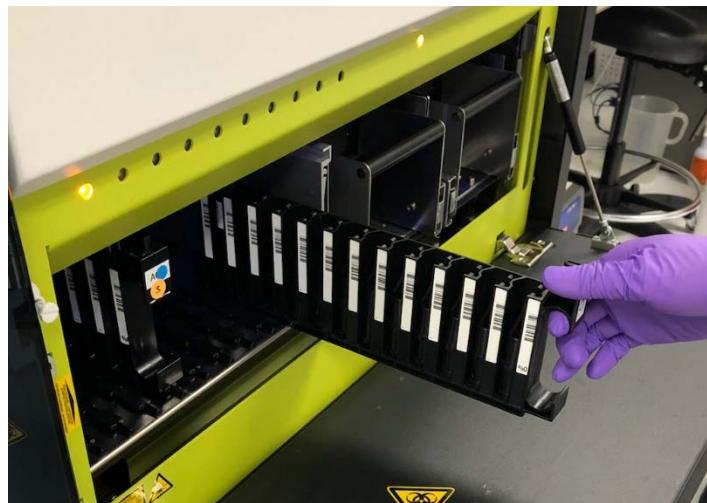

Loading Instrument Consumables



CAUTION! Use only Simoa cuvettes and Simoa disposable pipettor tips. Using any other type of cuvette or disposable pipettor tip may cause unreliable assay results or damage to the instrument. Quanterix cannot guarantee system performance if you use any other type of cuvette or disposable tip. For information on ordering consumables, see page 2.

Loading Cuvettes

The instrument uses 1-mL plastic cuvettes as reaction vessels. Cuvettes are snapped together into stacks of 50 to make them easy to load.

The cuvette holder has 10 chutes into which you insert stacked cuvettes. The instrument monitors the number of cuvettes in each chute. When a chute contains only five cuvettes, that chute will be considered empty. This chute will be refilled the next time by rotating to the front of the instrument when the cuvettes are loaded.



Note! Always load full stacks of 50 cuvettes. The instrument counts partial stacks as a single cuvette.



Single Cuvette



Cuvette Stack


❖ To determine whether you need to load cuvettes

Touch the System Resources tab. If there are insufficient cuvettes, the cuvette diagram is red.




❖ To load cuvettes

- 1 Touch the System Resources tab > touch the cuvette diagram. The Solid Resources screen appears. For information on reading this screen, see page 106.
- 2 Ensure the instrument status in the lower left corner shows the system is in a "Ready" state.
- 3 Unwrap a stack of cuvettes and insert it into the cuvette holder, as illustrated below, with the point down. The cuvette number in the Solid Resources screen will update.

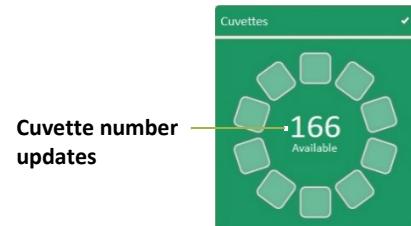



CAUTION! When loading the chute, ensure the stack does not become twisted (see below). Twisting can occur if the stack is dropped too quickly into the chute and may cause a jam in the cuvette loader.

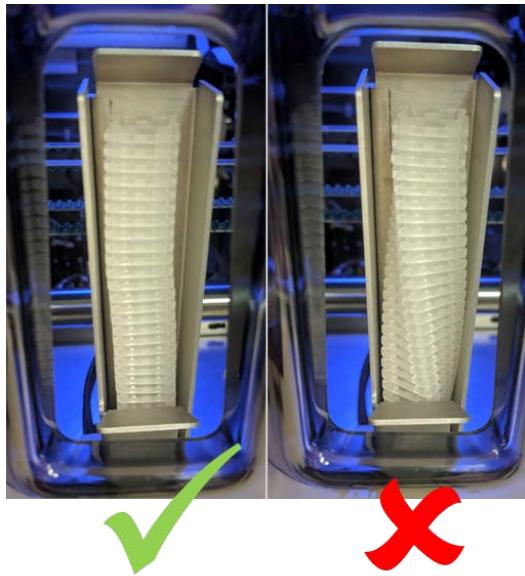

- 4 The cuvette loader rotates and the cuvette stack drops into the loader. “Maintenance” appears in the lower left corner of the screen. **Wait until the rotation is finished and the system status changes from “Maintenance” to “Ready” before loading a new stack.**

Loading Disposable Pipettor Tips

Each of the two System Resources drawers holds up to three disposable pipette tip racks. A full rack holds 96 tips, so the maximum number of tips that may be loaded at one time is 576, or 6 full racks.

Use only Quanterix-supplied tips to ensure proper instrument functioning and transfer of liquids. The instrument uses two tips per sample – one when the sample is picked up for processing and another when the sample is picked up for imaging.

When you load tips, you indicate in the software the position of each loaded tip rack. Doing so specifies the number of loaded tips and allows the software to count down from this number and to notify you if you need to load more tips.

❖ To determine whether you need to load disposable pipettortips

- 1 Touch the System Resources tab > touch one of the Tip diagrams to view the Solid Resources screen (for information on reading the Solid Resources screen, see page 106.)
- 2 If the System Resources drawer is locked, touch the **Unlock** button located under the System Resources Drawer diagram (see page 106).



Unlock

- 3 Pull out both System Resources drawers and verify that enough tips are loaded to complete the next planned assay run.
- 4 Double the total number of tests to be run in the assay to determine the number of tips required.


❖ To load disposable pipettor tips



CAUTION! The software only allows loading of full tip racks to prevent contamination that could result from handling tips.


-
- 1 Open the System Resources drawers > remove empty tip racks.

- 2 Remove the cover from new tip racks and place them in the drawer. Fill all rack positions.

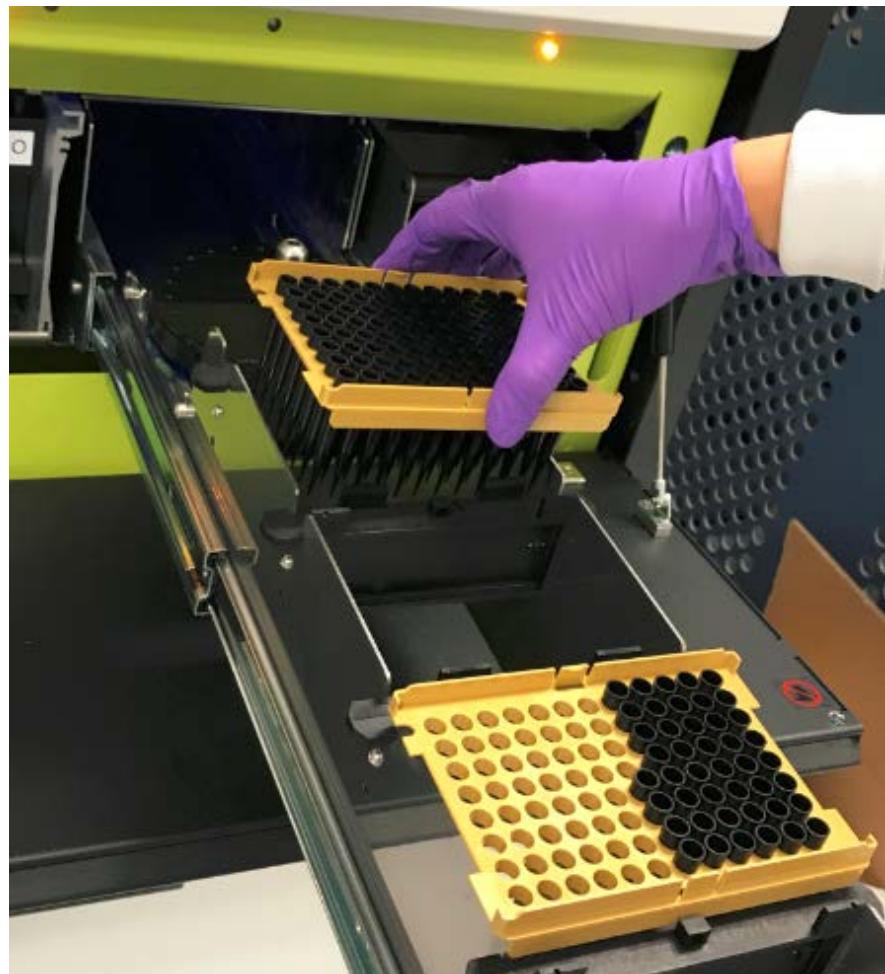

- 3 In the System Resources tab, touch one of the tip diagrams to open the Solid Resources screen.

- 4 In the Solid Resources screen, tap twice in the positions where you loaded new tip racks. The tip positions in the rack diagram turn light blue.

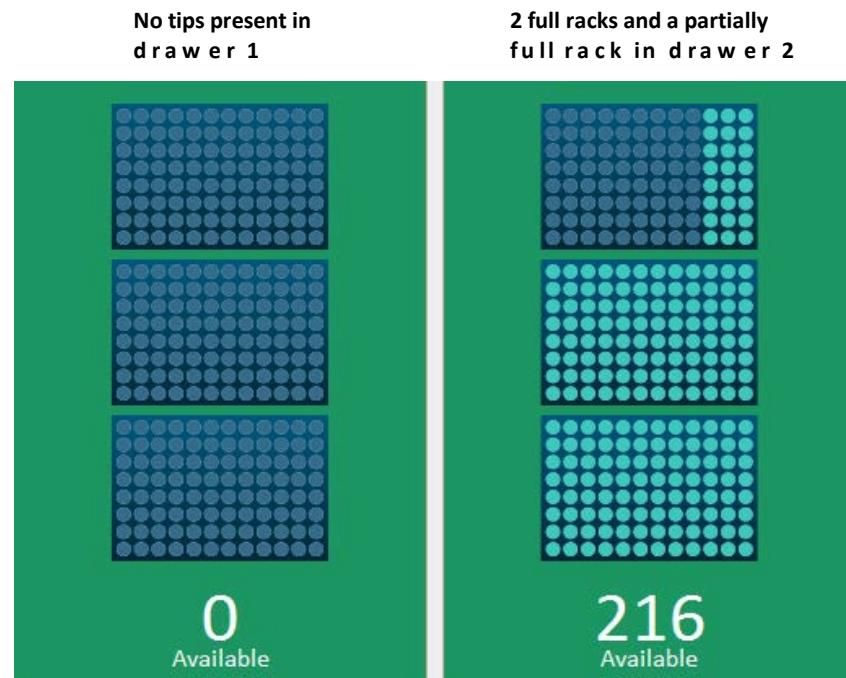

Loading Simoa Discs

Simoa Discs are loaded onto mounting poles at the rear of the System Resources drawers. One disc stack contains 16 discs. The maximum disc capacity is 32.



When you load a stack of discs, you scan the package barcode, which tells the software that you have loaded 16 discs. The instrument counts down from this number every time you consume a disc and notifies you if you need to load more discs.



CAUTION! Only load whole packages of Simoa Discs. Do not attempt to top off a partially used stack with new discs. To prevent disc contamination, do not remove the plastic wrapper from a stack of discs before you load the stack into the instrument.

❖ **To calculate the number of discs needed for a run**


Each disc has 24 arrays, and one array is required per job. Therefore, divide the number of jobs in the run by 24 and then round up to the nearest whole disc. For example: 96 jobs/24 = 4 discs; 78 jobs/24 = 3.25, round up to 4 discs.

To calculate the number of jobs in the run, multiply the number of samples by the number of replicates per sample.

- 1 Touch the System Resources tab > touch one of the disc diagrams to view the Solid Resources screen (for information on reading the Solid Resources screen, see page 106).
- 2 View the Solid Resources screen to determine if there are enough discs for your run.


❖ **To load Simoa Discs onto an empty disc pole**

- 1 Obtain one package of discs.
- 2 When the Solid Resources screen indicates there are zero discs available, touch **Unlock** to open the appropriate System Resource drawer.

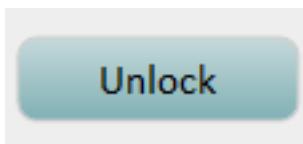

- 3 Pull out the System Resources drawer.

- 4 Remove the blue base plate from the disc pole. This plate sits at the bottom of every package of discs.

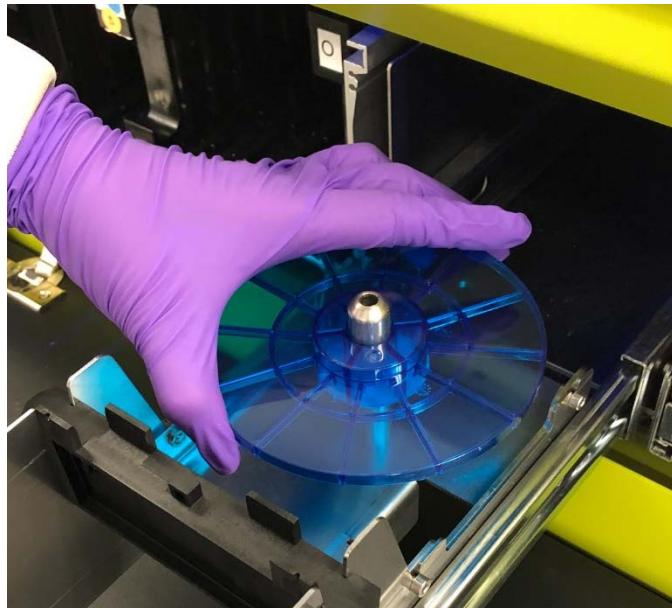


- 5 Load one package of discs:




Note! Only load a complete package of discs onto an empty disc pole. Do not add discs from the other stack.

- a Tap the disc barcode box.
- b If an existing barcode appears in the box, delete it.

- c Scan the package barcode with the barcode scanner. The barcode and the number of discs in the package (16) appear in the disc diagram.

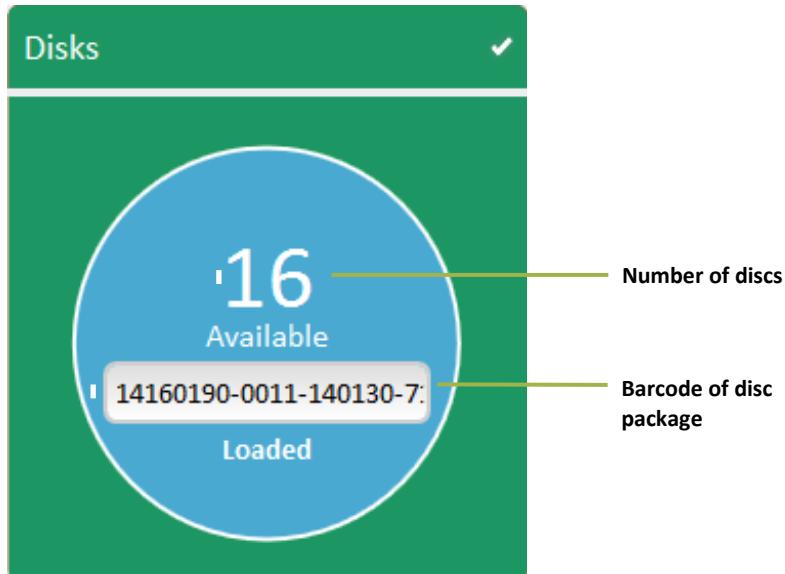

- d Press the package onto the mounting pole. The package label should be on top.

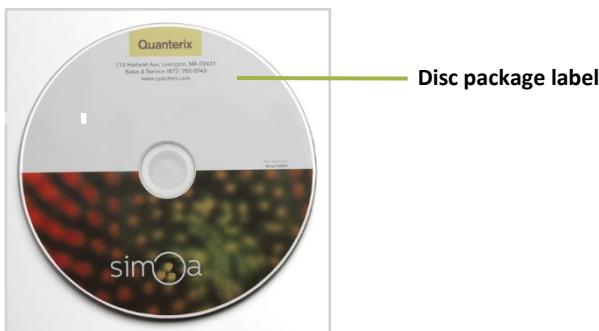

- e Line up the perforation with the top of the empty disc pole, slide the disc stack onto the pole, and orient the red pull tab toward the front.



- f Pull the red tear strip and remove the plastic disc covering completely.



- g Remove the top disc imprinted with the Simoa logo from the spindle.



Note! Take care not to shift discs against each other as they are very sensitive to scratching, which can affect image quality and data.

- 6 Push the drawer into the instrument until it catches.
- 7 Verify that the Resources screen reports the correct number of discs in both drawers.
- 8 Repeat “To load Simoa discs onto an empty disc pole” on page 44 if you need to load a second package of discs.

Refilling the System Fluid Containers

Refill the secondary system fluid containers for wash buffer 1 and deionized water and the system fluid container for wash buffer 2 before the first run of the day and whenever the container alert icons at the bottom of the software screen turn red.



Note that for wash buffer 1 and the deionized water, you only need to fill the secondary containers, not the primary containers. The primary containers passively refill themselves from the secondary containers. Wash buffer 2 only has one container, and it must be refilled.



CAUTION! Be sure to refill the system fluid containers with the correct solution. Filling with the wrong solution may require multiple flushes with water and/or the appropriate fluid in order to clean the lines.

- ❖ To refill a secondary system fluid container or wash buffer 2 container

- 1 Remove the container as described below.

- To remove the secondary containers for the deionized water and the system wash buffer 1, lift the containers out of their seats. As you lift the front end of the container, the plug connecting it to the primary container disconnects.



- To remove the system wash buffer 2 container, disconnect the sensor and fluid port and lift the container.



- 2 Unscrew the cap and refill the container with the correct fluid. Screw the container cap back on when you are finished.
- 3 Replace the container, as described below. **Do not force a container into its seat.** The connectors should mate easily. If you experience resistance, remove the container and try again.
 - To replace the deionized water and the system wash buffer 1 secondary containers, tip the connector end of the container down slightly and insert it into the seat.
 - To replace the system wash buffer 2 container, place it in the seat between the secondary deionized water and secondary wash buffer 1 containers. Reconnect the sensor and fluid port.

Changing the Sealing Oil Bag

The sealing oil bag hangs from the ceiling of the lower system bay and connects to the instrument by tubing. Change the bag when the sealing oil icon at the bottom of the software screen is red.




Note! Prior to running the Start of Day maintenance task, the sealing oil icon will appear red. An oil prime is performed as part of the Start of Day maintenance and the sealing oil icon status will be updated accordingly.




❖ To change the sealing oil bag

- 1 Close the tubing clamp just below the bag.

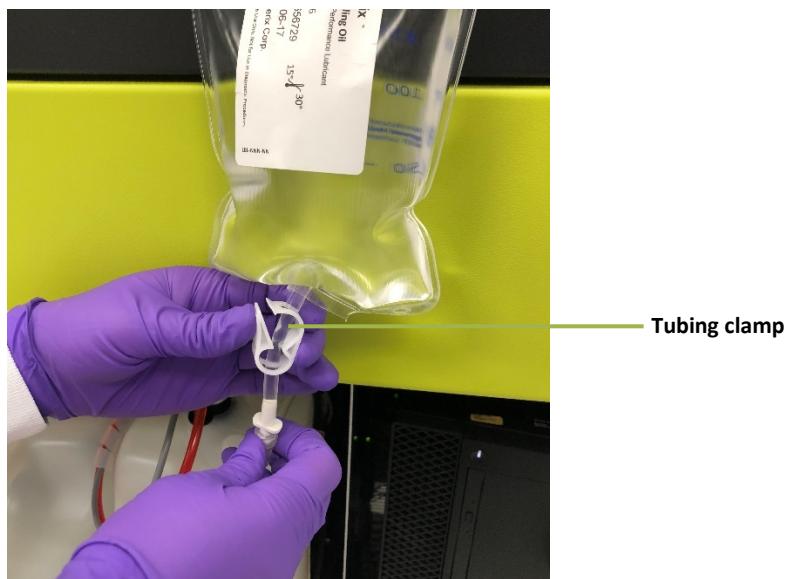

- 2 Unscrew the oil tubing from the bag.
- 3 Unhook the bag and remove from the system bay. Dispose of the bag according to the waste disposal policies of your organization.
- 4 Connect the hook to the new bag by piercing through the perforated circle at the top of the bag. The magnetic hook may be temporarily removed from the system bay ceiling to accomplish this.

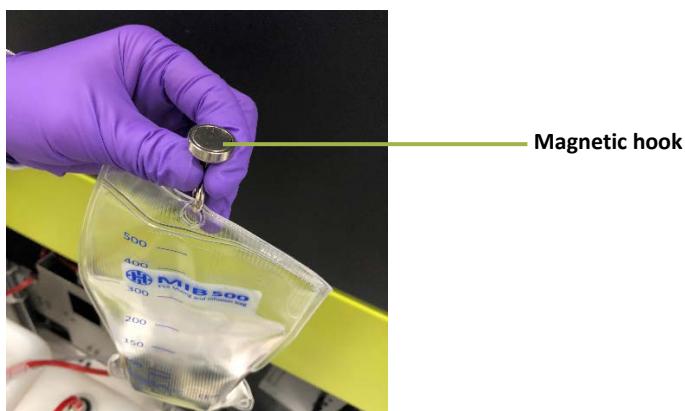

- 5 Confirm that the tubing clamp is closed on the new bag, then unscrew the yellow cap.



- 6 Connect the oil tubing to the new bag. Before connecting, twist the connectors in the opposite direction to allow “self-threading.”
- 7 Open the tubing clamp that is below the bag.
- 8 Touch the Maintenance tab > touch the Replenish Oil task checkbox > touch **Run Task**. The task takes approximately 15 minutes to complete. During the task, the instrument software is inactive. When the software screen reactivates, the task is complete.

Emptying the Solid Waste Container

Empty the solid waste container whenever the solid waste alert icon at the bottom of the software screen turns yellow or red.



The solid waste container holds discarded cuvettes, disposable pipettor tips, and Simoa Discs.

❖ To empty the solid waste container

- 1 Pull the solid waste drawer out of the system bay.




CAUTION! Do not place your hand in the opening in the shelf above the solid waste container when you pull the container drawer open. As the drawer slides out, it carries the container with it. If your hand is in the opening, your fingers or hand could be pinched between the edge of the opening and the rear top edge of the container.

- 2 Remove the waste bag > insert a new waste bag. Dispose of the waste bag according to the waste disposal policies of your organization.



CAUTION! The bag contains biohazard waste materials and should be handled with care.

- 3 Touch the System Resources tab > touch one of the disc diagrams to view the Solid Resources screen (for information on reading the Solid Resources screen, see page 106).
- 4 Push the solid waste drawer into the system bay.
- 5 Select “Yes” when you see the following dialog box to confirm that you emptied the solid waste container.

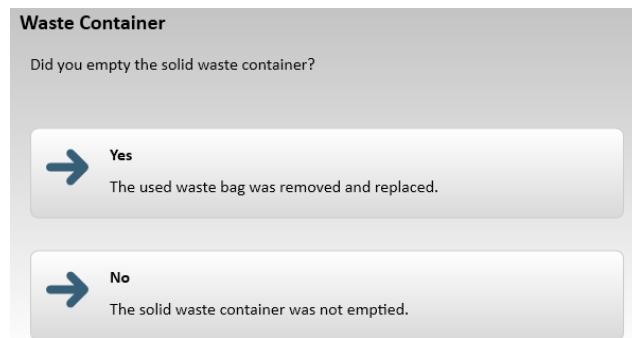

Emptying the Liquid Waste Containers

The instrument discards liquids into only one liquid waste container at a time. When a container is in use, a light in its cap turns on.


Empty both liquid waste containers before the first run of the day.

During a run, the instrument senses the fluid level in the containers with a floating sensor inside the containers.

To ensure decontamination of liquid waste, follow the instructions below to clean the waste container with the cleaner specified in your in-house decontamination procedures.

❖ **To empty the liquid waste container**

- 1 Disconnect the fluid tubing and the container sensor by pulling them straight up.
- 2 Once disconnected, place the container sensor cables into the holsters behind the waste bottles.

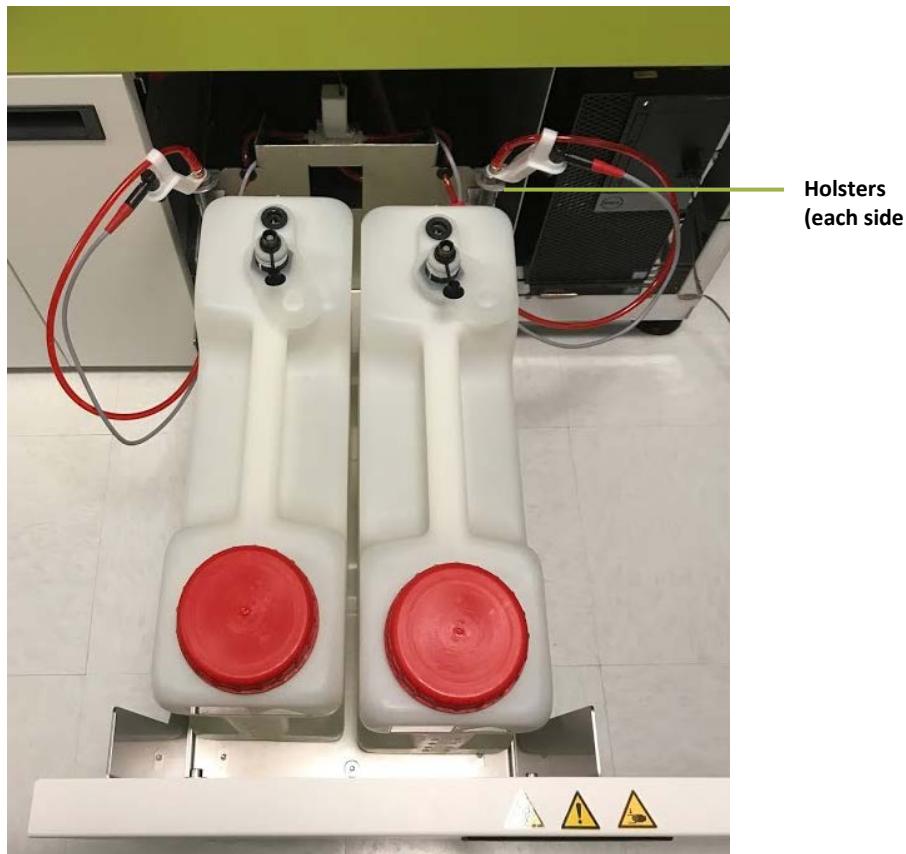

- 3 Remove the container and dispose of the liquid waste according to the waste disposal policies of your organization.



CAUTION! The liquid waste container holds biohazard waste and should be handled with care.

- 4 Place cleaner into the container as specified by your in-house cleaning procedures and place it back in the system bay.
- 5 Reconnect the tubing and the container sensor. When done correctly, you will hear a click when you reinsert each sensor.



Monitoring the Instrument's Internal Temperature

The instrument's internal temperature is displayed in the upper left corner of the System Resources screen. See Understanding the System Resources Screen on page 105 for more information.

The recommended internal temperature range is between 18.0°C (64.4°F) and 25.0°C (77.0°F). The system will notify you if the temperature falls outside of this range. If this happens, read the notification and follow the instructions on the screen to restore the instrument to an acceptable temperature.



CAUTION! Do not start a run while the internal instrument temperature is above or below the recommended range. Assay data may be compromised and any individual jobs processed while the temperature is out of range will be flagged.

6 Importing Assay Definitions

Understanding Assay Definitions	59
Checking Imported Assay Definitions	59
Importing an Assay Definition.....	61
Configure Assay Settings.....	62
Read-Only.....	63
Visible Only to Me	63
Archive	64

Understanding Assay Definitions

The system software provides a range of pre-defined 1-, 2-, and 3-step incubation protocols that specify reagents, calibrator levels and concentrations, reagent volumes, and incubation times. (For information on the protocol steps, see the instructions in your assay kit.)

An **assay definition** is an XML document that specifies the pre-defined protocol that a given Simoa assay uses and all of the settings for the protocol. It also includes information on curve fitting and concentration result generation.

To perform a Simoa assay, you must import an assay definition in the system software.

Checking Imported Assay Definitions


Before you set up a run, make sure that the assay definition you need has been imported on your Simoa HD-X Analyzer.

❖ To check imported assay definitions


- Touch the **Custom Assay tab** > touch the **Assay** field under **General**.
The **Select Assay** dialog opens, showing all assays that have been defined for the system.



Note! For this specific application, enter SARS-CoV-2 IgG for the assay shown below. This applies to subsequent images in this chapter that reflect IL-10.

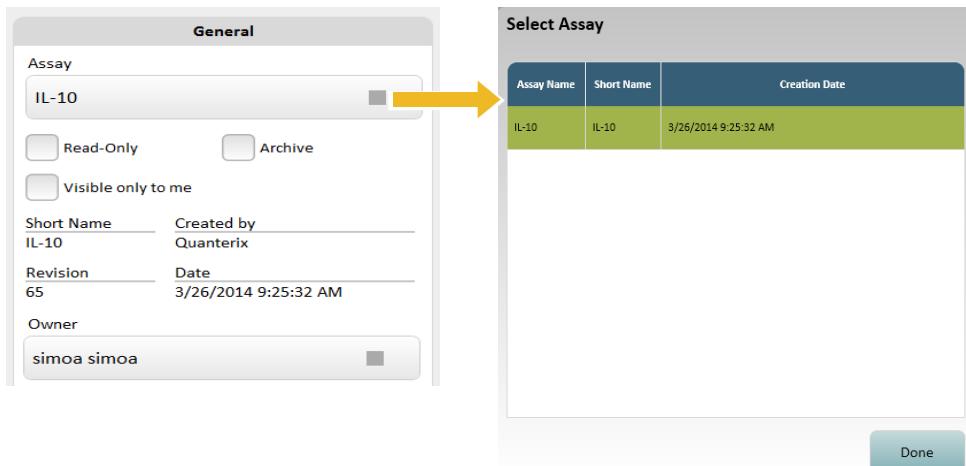

- Alternatively, you can check imported assay definitions in the **Setup Run tab**. Touch the **Setup Run tab** > turn off the **Assays Loaded only** toggle to display all assay definitions that have already been imported into the system. Note that the **Loaded only** toggle is set to off by default.

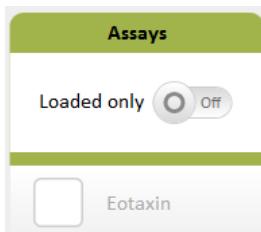


Importing an Assay Definition

You can download assay definitions from the Quanterix Customer Support website (<http://portal.quanterix.com>).

If you are importing an assay definition that shares reagents with existing assays, you can replace the existing reagent definition with the imported assay definition, or you can leave the existing definition for those shared reagents unchanged.

❖ To import an assay definition

- 1 Touch the Custom Assay tab.
- 2 In the Assay Overview screen, touch **Import**.
- 3 Navigate to the assay definition you want to import and touch **Open**.
- 4 Choose the importing options:
 - a Select **Update the existing reagent definition** to replace existing shared reagent definitions with the values specified in the imported assay definition.
 - b Select **Keep the existing reagent definition** to leave existing shared reagent definitions unchanged. The imported assay uses the values of the previously specified reagents instead of the ones in the imported assay definition.
 - c Select **Cancel** to cancel the import.

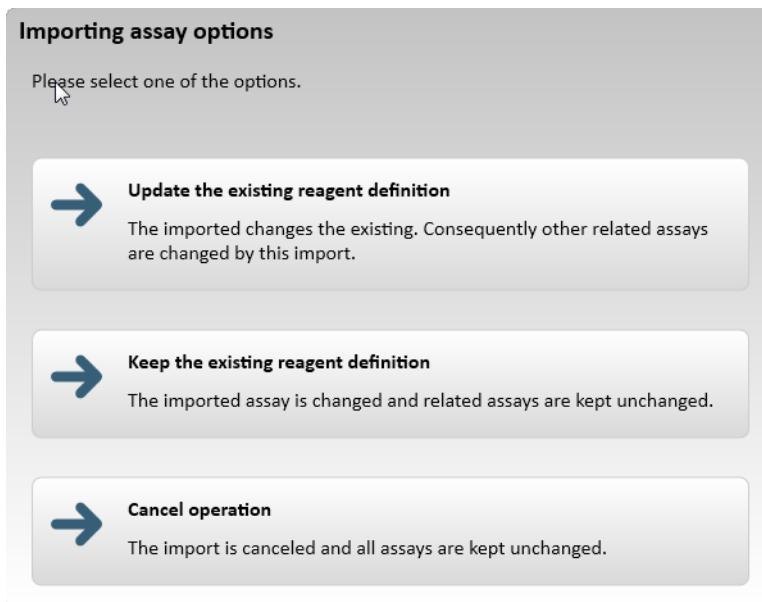

5 Optionally, edit the Assay Name or the Assay Short Name.




6 Touch **Save**.

Configure Assay Settings

You can specify the following assay settings that control how assays are displayed and accessed:

- Read-Only
- Visible Only to Me
- Archive

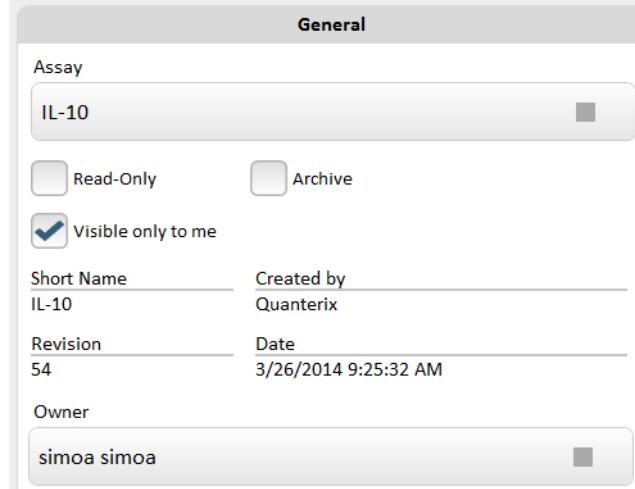

The 'General' settings screen for assay 'IL-10'. It includes fields for 'Assay' (IL-10), 'Read-Only' (unchecked), 'Archive' (unchecked), and 'Visible only to me' (checked). Other fields show 'Short Name' (IL-10), 'Created by' (Quanterix), 'Revision' (54), 'Date' (3/26/2014 9:25:32 AM), and 'Owner' (simoa simoa).

Read-Only

Depending on how roles and features are configured for your system, you may or may not see a Read-Only button on the Custom Assay screen.

Typically, administrators hide this button from users with non-administrator roles.


By default, imported assays are read-only and cannot be edited. Administrators can disable the Read-Only button to allow editing of assay definition values.

For more information, see “Managing Roles and Features” on page 155.

Visible Only to Me

Your assay definitions are associated with your user profile to prevent other users from accidentally modifying them and to limit the number of assays displayed on the Setup Run tab.

To protect your assays, check the Visible Only to Me checkbox.

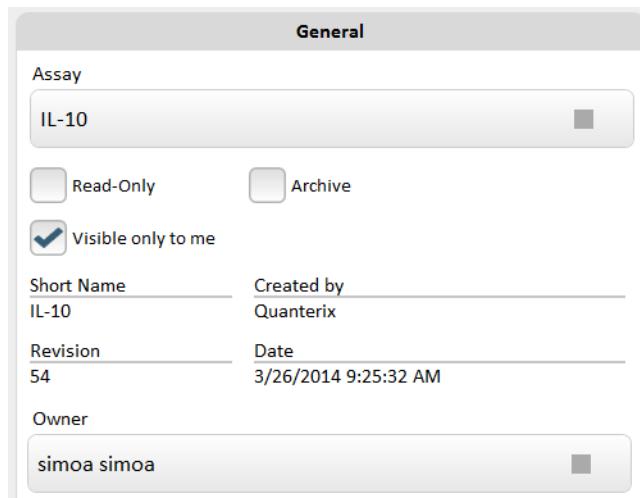

The screenshot shows the 'General' tab of an assay configuration dialog. The 'Assay' section contains the name 'IL-10'. Below it are three checkboxes: 'Read-Only' (unchecked), 'Archive' (unchecked), and 'Visible only to me' (checked). Underneath these are two rows of information: 'Short Name' (IL-10) and 'Created by' (Quanterix); 'Revision' (54) and 'Date' (3/26/2014 9:25:32 AM); and 'Owner' (simoa simoa).

When you log in, you will see all assays that you defined and all assays created by others, except for those configured as Visible Only to Me when created by other users.

Archive

Archive assays that are not currently in use to limit the number of assays displayed in the other tabs. Check the **Archive** checkbox to hide an assay or uncheck it to display it again. When you archive an assay, it is stored in the database but is not visible in the **Setup Run** tab under **Assays**.

General

Assay	
IL-10 <input type="button" value="..."/>	
<input type="checkbox"/> Read-Only	<input checked="" type="checkbox"/> Archive
<input type="checkbox"/> Visible only to me	
Short Name	Created by
IL-10	Quanterix
Revision	Date
56	3/26/2014 9:25:32 AM
Owner	
simoa simoa <input type="button" value="..."/>	

7 Performing an Assay Run

Steps in an Assay Run.....	67
The Role of Barcodes in an Assay Run	67
Manually Entering a Reagent Barcode.....	67
Preparing Reagents, Calibrators, Controls, and Samples.....	68
Use Approved Sample Plates and Sample Tubes	69
Preparing Reagents, Calibrators, Controls and Samples.....	69
Preparing the Instrument for the Run	69
Loading Reagents and RGP into the Instrument.....	71
Using the Barcode Scanners.....	71
Understanding the Load Reagents Screen	72
Loading Reagents.....	73
Loading RGP.....	77
Setting Up a Run.....	80
Understanding the Assign Calibrators Screen.....	81
Understanding the Assign Samples Screen.....	83
Understanding the Plate Position Diagram.....	85
Selecting Processing Order	85
Selecting Plate Wells and Removing Assignments.....	85
Understanding the Required Sample Volume Field.....	86
Understanding the Autoname Tool	86
Scanning Sample Barcodes.....	87
Assigning Calibrators.....	88
Assigning Samples	96
Loading Calibrators/Samples into the Instrument.....	102
Importing and Exporting a Layout.....	103
Starting a Run	103
Understanding the System Resources Screen	105
Understanding the Solid Resources Screen	106
Monitoring a Run	107



Understanding the In-Process Run Screen for Plates	107
Understanding the In-Process Run Screen for Tubes	108
Checking Instrument Status	109
Viewing the Event Log.....	109

Steps in an Assay Run

Although the examples in this chapter show plates of samples and calibrators, the instructions apply to tubes as well.

Performing a Simoa assay run involves the following general steps:

- 1 Prepare reagents, samples, and calibrators (if you will be running them).
- 2 Select Rack or Plate mode for samples. Prepare the instrument for the run.
- 3 Place bead, detector, and SBG reagents into a reagent rack, identify the reagents by position and barcode in the software, and load the rack into the instrument.
- 4 Place one or more vials of RGP into an RGP rack, identify the vials by position and barcode in the software, and load the rack into the instrument.
- 5 Place samples, calibrators, and controls in plate or tube racks, select an assay to run on them, name the run, name the rack or plate, identify the samples, calibrators, and controls in the software, and load the racks into the instrument.
- 6 Start the run and monitor it.

When the run finishes, the results appear in the Run History table in the History & Reports tab. For information on analyzing results, see “Analyzing Run Data” on page 111.

The Role of Barcodes in an Assay Run

The labels on all assay reagents and Simoa Discs have barcodes that provide information such as reagent or item name, fill volume (as specified in the reagent definition), lot number, and expiration date. When you load barcoded reagents or RGP into the instrument, you can scan them with the handheld scanner or use the on-board scanner. When loading barcoded Simoa Discs, you must manually scan their barcodes using the handheld device.

The software uses the scanned barcode data to:


- Identify loaded reagents in screens.
- Monitor expiration dates.
- Calculate reagent and consumable consumption during an assay run and notify you when more reagents or consumables should be loaded.

Be sure to follow instructions in this guide regarding barcode scanning. If barcodes are not scanned, the instrument will not be able to perform the assay.

Manually Entering a Reagent Barcode


In the event that a reagent label is damaged and cannot be scanned, you will need

to type the barcode number into the software. The human-readable version of the barcode is printed beneath the barcode.

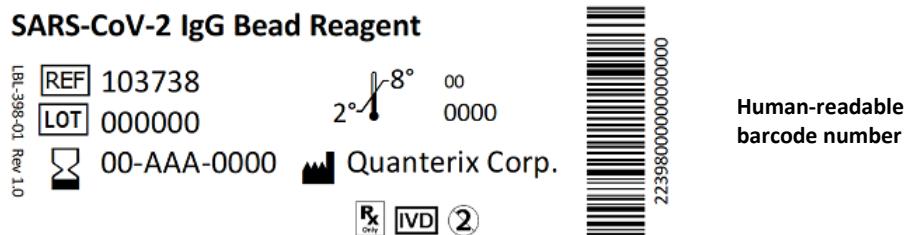

❖ To manually enter a reagent barcode

In the Load Reagents tab, touch the Reagent Barcode field > type the number > touch **Enter**. For information on loading reagents, see “Loading Reagents and RGP into the Instrument” on page 71.

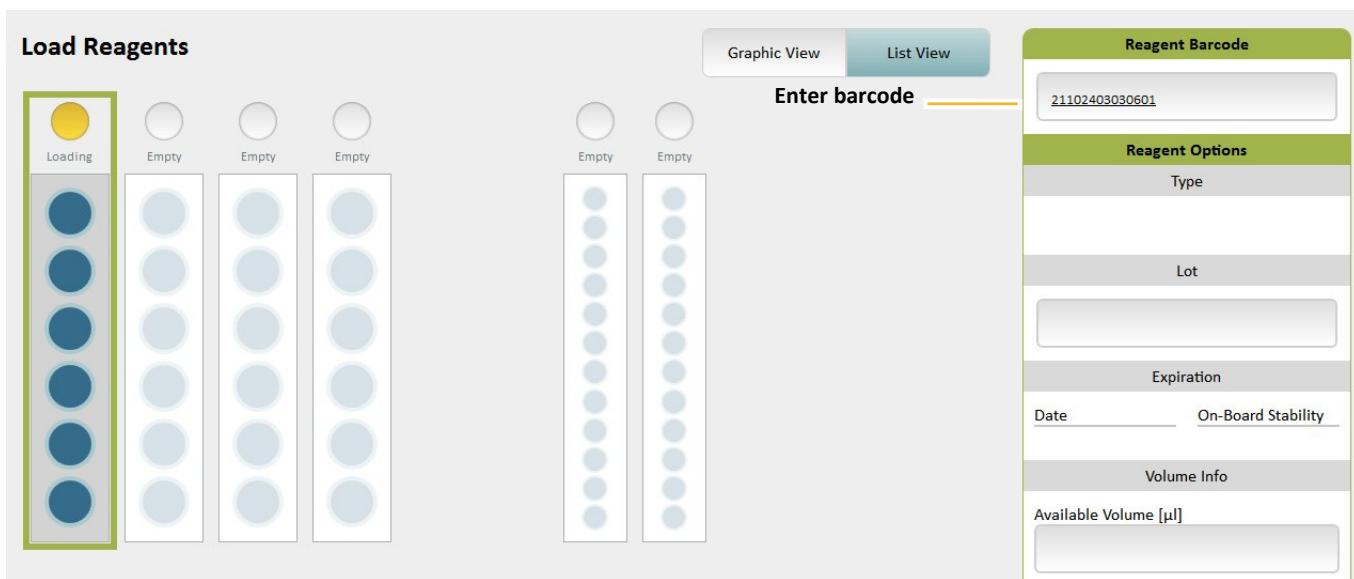

Preparing Reagents, Calibrators, Controls, and Samples

Before preparing any assay materials or performing an assay on the instrument, read the assay kit instructions. The kit instructions contain important safety information and details about the assay procedure that are vital to performing the assay successfully.

Use Approved Sample Plates and Sample Tubes

Use only Quanterix®-supplied sample plates. These plates have been verified to be nonreactive with Simoa assay reagents and to be correctly sized to work with the instrument pipettors. For information on ordering plates, see page 2.

Please contact Quanterix Technical Support for information regarding which sample tubes you may use with your HD-X instrument.



CAUTION! Quanterix cannot guarantee system performance if you use any type of plate other than the ones supplied by Quanterix.

Preparing Reagents, Calibrators, Controls and Samples

Follow the relevant test kit Instructions for Use to prepare reagents, calibrators, controls and samples.

Preparing the Instrument for the Run

For information on how to load reagent racks, plate and tube racks, and RGP racks, see pages 32 to 37.

❖ To prepare the instrument

- 1 If this is the first run of the day, perform the instrument startup procedure (see page 28).
- 2 Load Simoa cuvettes, Simoa disposable pipettor tips, and Simoa Discs into the instrument as needed. For information on how to load these consumables, see page 38.
- 3 Set up the software for sample plates or sample tubes:
 - a Touch the Settings icon at the top right of the software to open the System Settings screen.



Settings Icon

- b Touch **Plate Mode** or **Rack Mode** > touch **Done**.
- 4 Obtain enough reagent racks, RGP racks, and either tube racks or plate racks to perform the assay.

- 5 Make sure that the necessary assay definition is imported on the instrument.
For information on checking this, see page 59.

Loading Reagents and RGP into the Instrument

Loading reagents and RGP involves:

- Placing reagent bottles into a reagent rack.
- Enabling the on-board barcode scanner or scanning reagent label barcodes with the handheld scanner.
- Inserting the rack.

Load as follows:

- Insert the rack into any lane of the reagent bay.
- Load detector and SBG reagents into any position in the rack.
- Load beads into the first three (shaking) positions.
- Load RGP into lanes labeled RGP 1 or RGP 2.

Using the Barcode Scanners

You can use the on-board or handheld scanners to scan reagent and RGP barcodes. When the on-board barcode scanner is enabled, it automatically reads the reagent/RGP barcodes and associates each barcoded bottle with its position in the rack. When the on-board scanner is turned off, you must identify the rack position of the reagent bottles in the Load Reagents screen and then scan their barcodes using the handheld scanner.

You can manually enter one or more reagent barcodes (either by using the handheld scanner or by typing it) and still use the built-in scanner for the remainder of reagents on the rack.


When the on-board scanner is turned on, it reads all reagent barcodes in the rack, including those that have already been scanned with the handheld device or typed in. If a barcode detected by the built-in scanner does not match the previously entered barcode, then the barcode and all associated reagent information will be overwritten by the on-board scanner.

The software associates the barcodes with their positions in the rack so that the instrument pipettors can be directed to the correct bottles in the loaded rack. The barcode data also allows the software to calculate running reagent volume and count down to reagent expiration dates.

Occasionally, a barcode may not be scannable because the label has been damaged in some way. When this happens, manually enter the barcode into the appropriate software field (see page 67).

Understanding the Load Reagents Screen

By default, the screen shows a reagent rack diagram, as shown below. When you touch one of the RGP lanes, the screen changes to show an RGP rack diagram. To view an example of the RGP rack diagram, see page 77.

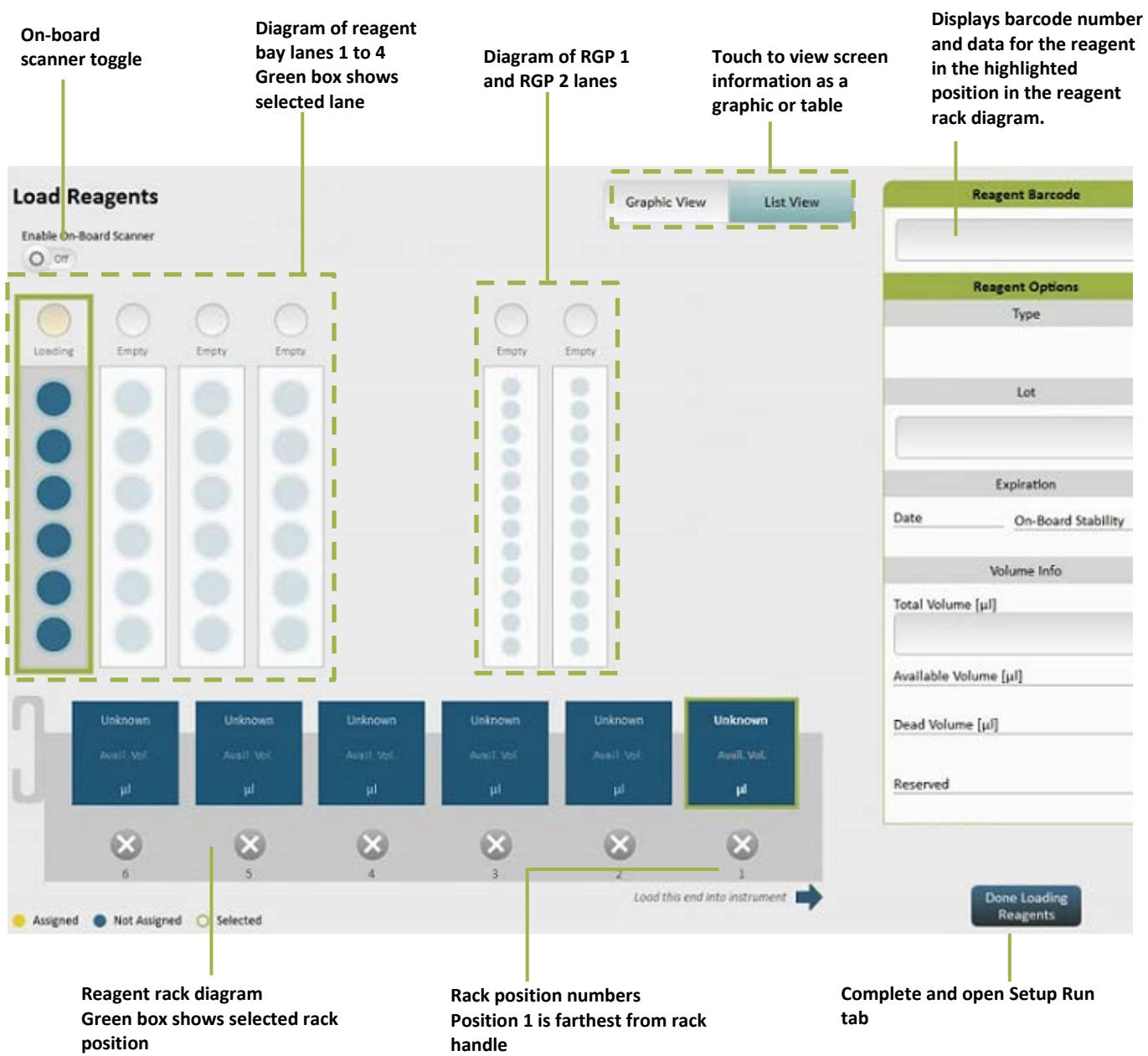

Loading Reagents

Loading reagents involves placing the bead, detector, SBG, and sample diluent reagents in the reagent rack and inserting the rack into the reagent bay.

❖ To load reagents using the handheld barcode scanner

- 1 Mix the bead reagent according to instructions in the assay kit. If the bead reagent is left sitting stationary for more than five minutes before you load it, repeat this step.
- 2 Remove and set aside the caps from the reagent bottles.
- 3 Insert the bead reagent into reagent rack position 1, 2, or 3.



Note! You must insert the bead reagent in rack position 1, 2, or 3. These positions shake during the assay run to keep the beads suspended.


- 4 Insert the other reagents into any other reagent rack position, including unoccupied shaking positions (1, 2, or 3).

- 5 Turn the bottles so that the reagent label barcodes are centered in the notch of the bottle holder. The barcodes must be unobstructed so that they can be scanned accurately.



- 6 Touch the Load Reagents tab > touch a reagent bay lane in the screen. This step identifies the lane into which you will load the reagent rack.

Reagent Bay position 1
has been selected

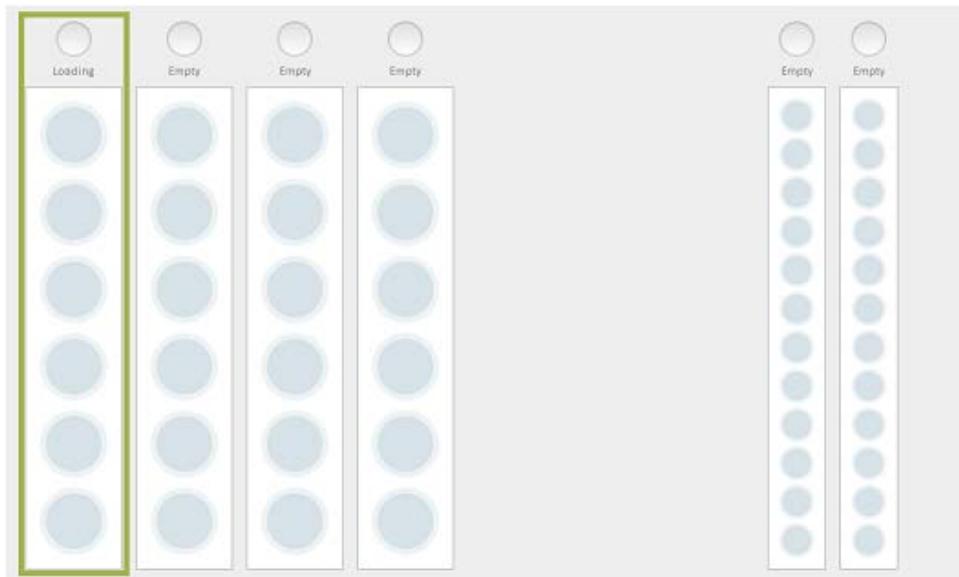

7 Scan the bottle barcodes:

- a In the reagent rack diagram, touch the position that corresponds to the position of the bead reagent bottle in the reagent rack.



- b Touch the reagent barcode field and scan the bead reagent barcode with the handheld scanner. The barcode and data about the bead reagent appear in the green box in the screen. If a barcode does not scan, enter it manually (see “Manually Entering a Reagent Barcode” on page 67).
- c Repeat these two steps for the rest of the reagents in the rack.



Note! If you have not inserted the bead reagent into position 1, 2, or 3, the software displays a warning, and you cannot proceed until you move the bottle into one of these positions.

- 8** Insert the rack in the reagent bay lane that you selected in step 6 until it clicks into place. The rack diagram updates to show the short name of each loaded reagent in its identified rack position. The lane diagram updates to show the positions of reagents in the selected lane. The Reagent Barcode box shows barcode data for the position that is highlighted in the reagent rack diagram.

❖ **To load reagents using the on-board barcode scanner**

- 1 Perform **Steps 1 through 5** as described in the previous section, “To load reagents using the handheld barcode scanner”.
- 2 Touch the Load Reagents tab. On the Load Reagents screen, turn on the “Enable On-Board Scanner” toggle.




- 3 Touch a reagent bay lane in the Load Reagents screen. After you select a lane, the on-board scanner focuses its laser on the lane and a message appears to notify you that the barcode scanner is focusing. Wait until focusing completes and the message disappears before proceeding to the next step.

Barcode scanner is focusing...

- 4 Insert the loaded rack in the reagent bay lane that you selected in step 3 until it clicks into place. As you load it, the on-board scanner reads the barcodes of each reagent bottle in the rack. The rack diagram updates to show the short name of each loaded reagent in its identified rack position. The lane diagram updates to show the positions of reagents in the selected lane. The Reagent Barcode box shows barcode data for the position that is highlighted in the reagent rack diagram.



Note! The on-board scanner will only work if the red scanner light is visible while you are loading the reagent rack. It appears on the left side of the reagent bay when the scanner is turned on.

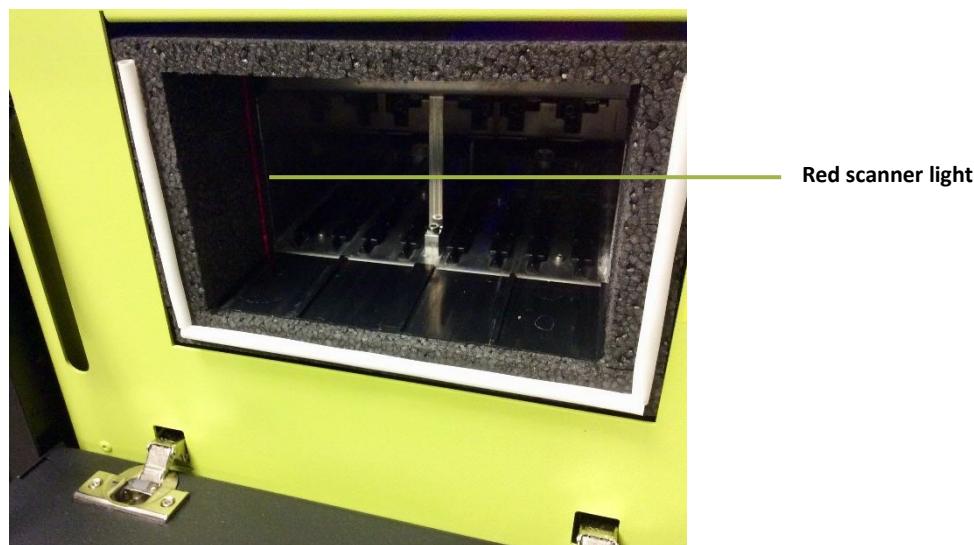

Loading RGP

RGP must be loaded into an RGP rack and placed in sample bay RGP 1 or RGP 2 only (the substrate lanes).


The RGP racks are customized to hold the short RGP reagent vials at optimum position for the instrument pipettors.



CAUTION! Do not load RGP into a sample tube rack. Always use an RGP rack. If you load RGP into a sample rack, the pipettor will strike the bottom of the loaded vials, possibly damaging the instrument.

❖ To load RGP using the handheld barcode scanner

- 1 Prepare RGP according to the assay kit instructions.
- 2 Remove the caps of one or more RGP vials > insert the RGP vials into the RGP rack.
- 3 Touch the Load Reagents tab > touch lane RGP 1 or RGP 2 in the screen. This step identifies the lane into which you load the RGP rack. The Load Reagents screen shows an RGP rack diagram.

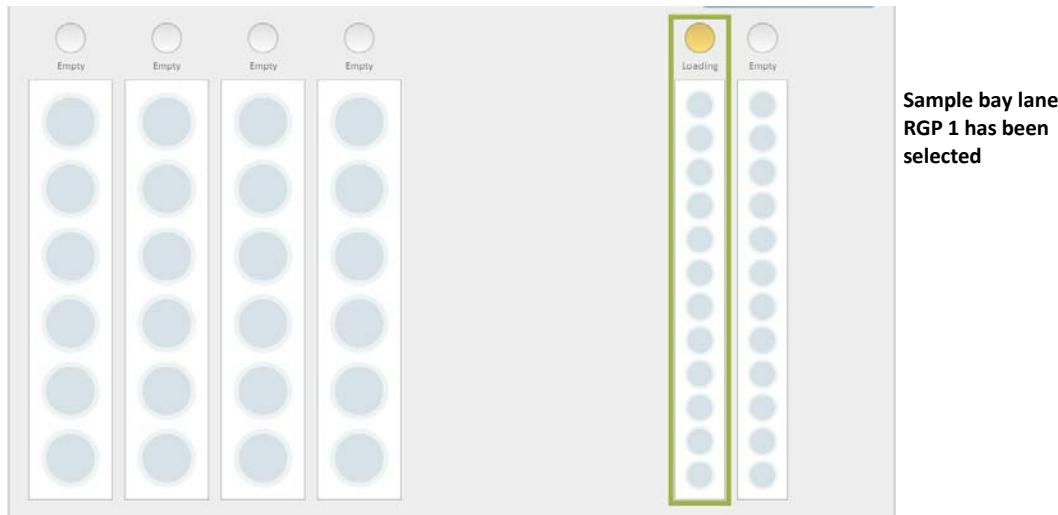

4 Scan the vial barcodes:

- a In the RGP rack diagram, touch the position that corresponds to the position of the RGP vial in the reagent rack.



- b Scan an RGP vial barcode with the handheld scanner and insert the vial in position 1 of the rack. The barcode and data about the scanned vial appear in the Load Reagents screen.
 - c Repeat these two steps for the rest of the RGP vials.
- 5 Insert the rack into the lane that you selected in step 3 until it clicks into place. The RGP rack diagram updates to show each RGP vial in its identified rack position. The lane diagram updates to show the RGP vial positions in the selected lane.

❖ **To load RGP using the on-board barcode scanner**

- 1 Prepare RGP according to the assay kit instructions.
- 2 Remove the caps from one or more RGP vials > insert the RGP vials into the RGP rack.
- 3 Turn the RGP vials in the rack so that the barcodes are centered in the notch of the bottle holder. The barcodes must be unobstructed so that they can be scanned accurately.
- 4 Touch the Load Reagents tab. Turn on the “Enable On-Board Scanner” toggle in the upper left part of the screen.




- 5 Touch lane RGP 1 or RGP 2 in the Load Reagents screen. Once you have selected an RGP lane, the on-board barcode scanner focuses its laser on the lane and a message appears to notify you that the scanner is focusing. Wait until the message disappears before proceeding to the next step.

Barcode scanner is focusing...

- 6 Insert the rack into the lane that you selected in step 4 until it clicks into place. As you insert the rack, the on-board barcode scanner reads the barcode for each RGP bottle. The RGP rack diagram updates to show each RGP vial in its identified rack position. The lane diagram updates to show the RGP vial positions in the selected lane.

Note! The on-board scanner will only work if the red scanner light is visible while you load the RGP rack. The red light appears on the right side of the sample bay when the scanner is turned on.

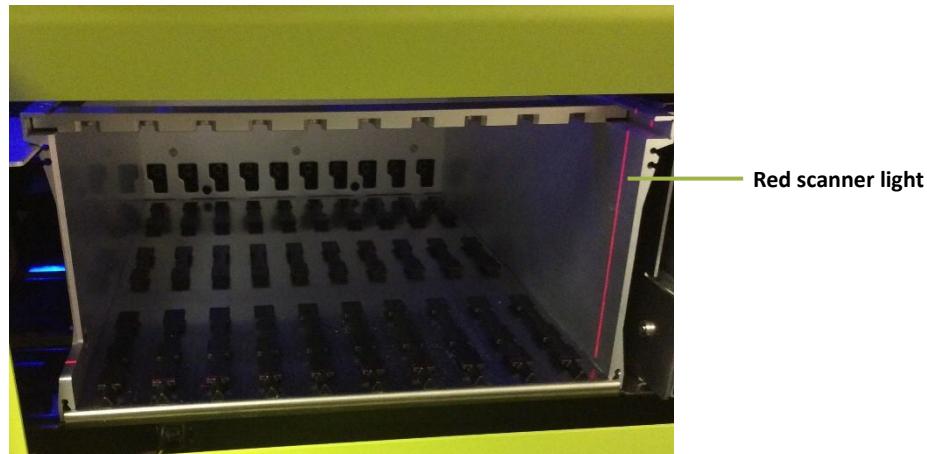

Setting Up a Run

Setting up the run involves identifying calibrators and samples in the Setup Run tab. In this tab, you can toggle between the Assign Calibrators and Assign Samples screen.

In this chapter, the examples describe a single plate that contains all of the calibrators. You can also set up your assay calibrators on multiple plates.

Once you have set up a run, you can export the assignments to a file and import the mappings the next time you run the assay.

If an experiment requires multiple assay runs, you can assign a **Batch Name** to all the plates or tube racks that you run for the experiment. This procedure allows you to group the experiment's results in the Run History screen.

Assigning calibrators and samples involves:

- Selecting an assay.
- Identifying in the software the physical location of calibrators and samples in the plates or tubes that you will be loading into the instrument.
- Associating each calibrator with one of the calibrator concentrations in the assay definition.
- Selecting the number of replicates to run. See “**Error! Reference source not found.**” on page **Error! Bookmark not defined.** for a complete explanation.
- Naming each sample, selecting the number of replicates to run, and, depending on the assay, selecting a dilution factor for the sample.
- Selecting the sample processing order (Rows or Columns).



Note! You will find it easier to assign calibrators and samples if you put them into plate racks or tube racks and then place the racks in the sample bay before you begin. Do not fully insert the racks; simply place them into the bay. The procedures for assigning calibrators and samples include directions to do this.


Understanding the Assign Calibrators Screen

The Assign Calibrators screen shows a plate diagram or a tube diagram depending on how you have set up the instrument (see “Preparing the Instrument for the Run” on page 69).

For information on using the well selection tools, see “Selecting Plate Wells and Removing Assignments” on page 85.

Use Multi Select to assign replicates for multiple calibrators.

The Assign Calibrators screen in plate view appears below.

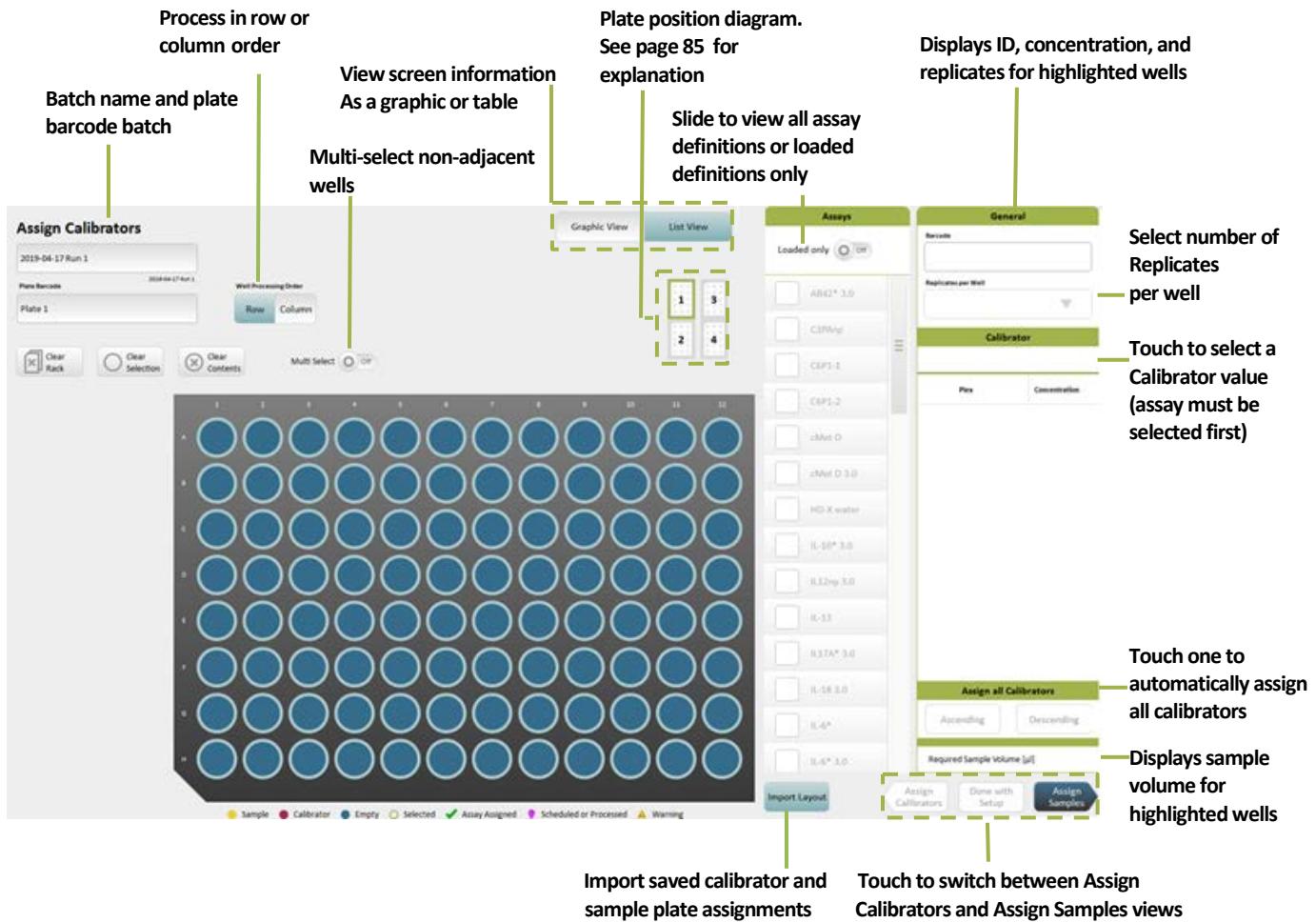

The Assign Calibrators screen in tube rack view appears below. It is identical to the screen for plates except that it shows a tube rack diagram, a location diagram for the eight sample bay lanes, and the Enable On-Board Scanner toggle.

Tube position diagram (sample bay rows 1-8)
Green box shows selected lane

On-board barcode scanner toggle

Touch to specify replicates

Assign Calibrators

Batch [MM/dd/yyyy HH:mm]
8910 03/04/2019 13:40

Rack Barcode

None											
1	2	3	4	5	6	7	8	9	10	11	12

Load this end into instrument

Legend: ● Sample ● Calibrator ● Empty


Icons: ● Selected ✓ Assay Assigned ● Scheduled or Processed ⚠ Warning

Tube rack diagram for selected lane

Tracks volume of calibrator used

Understanding the Assign Samples Screen

The Assign Samples screen in **plate view** appears below.

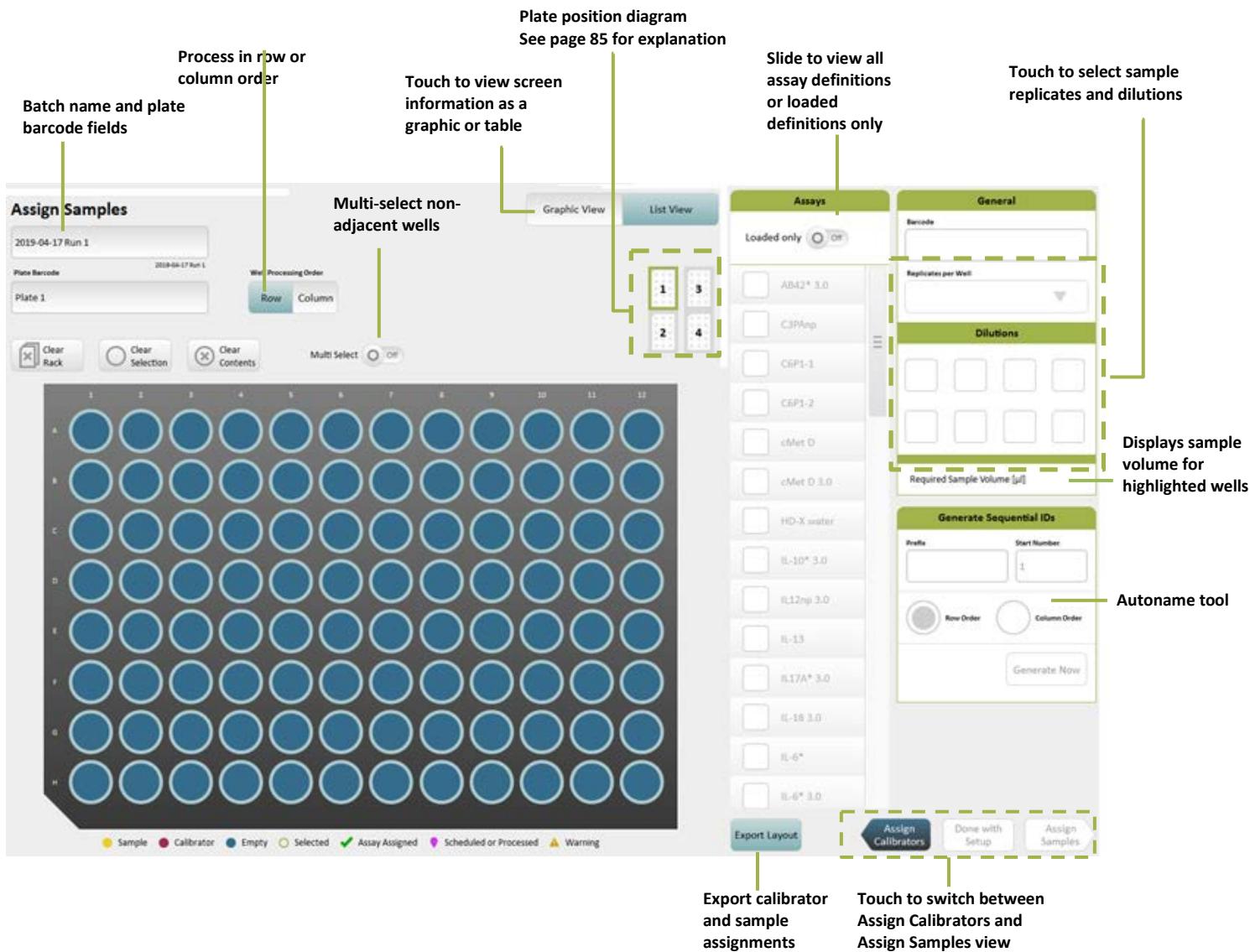

The Assign Samples screen in tube rack view appears below. It is identical to the Assign Samples screen for plates except that it shows a tube rack diagram, a location diagram for the eight sample bay lanes, and the Enable On-Board Scanner toggle.

**Tube position diagram (sample bays rows 1-8)
Green box shows selected lane**

On-board barcode scanner toggle

Graphic View List View

Multi Select Tubes On Off

Enable On-Board Scanner On Off

Touch to view screen information as a graphic or table

Assays General

Barcode

Replicates per Well

Dilutions

Required Sample Volume [μl]

Generate Sequential IDs

Prefix: Start Number: 1

Generate Now

Assign Samples

Batch (MM/dd/yyyy HH:mm)
Batch: 02/04/2018 13:40

Rack Barcode: Clear Rack

On-board barcode scanner toggle: Graphic View (highlighted), List View, Multi Select Tubes, Select All Tubes, Deselect All Tubes, Enable On-Board Scanner (highlighted), On, Off.

Tube position diagram for selected lane: Shows a 4x12 grid of sample bays. The first column (row 1) is highlighted with a green border and contains a yellow circle icon labeled 'Loading'. Other columns are labeled 'Empty'.

Autoname tool: Shows fields for Prefix (empty), Start Number (1), and Generate Now button.


Bottom buttons: Export Layout, Assign Calibrators, Done with Setup, Assign Samples.

Legend: Sample (yellow dot), Calibrator (red dot), Empty (blue dot), Selected (green border), Assay Assigned (green checkmark), Scheduled or Processed (purple dot), Warning (yellow triangle).

Understanding the Plate Position Diagram


The plate position diagram identifies the lanes into which a plate is loaded and the position of a plate in those lanes. Plate rack 1 may only be loaded into lanes 1–4, and plate rack 2 into lanes 5–8. There are two plate positions for each plate rack.

Selections in the diagram apply to both the Assign Calibrators and the Assign Samples screens. For example, if you select plate 2 in the Assign Calibrators screen and then switch to the Assign Samples screen, the plate remains selected.

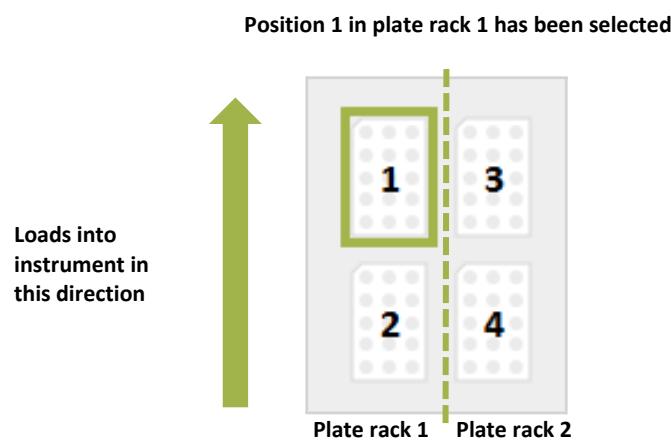

Selecting Processing Order

You can process your calibrators and samples by Row (left to right) or Column (top to bottom). This is a batch-wide setting that applies to all jobs in a run, even if on separate plates. If you have multiple plates, they are processed in ascending numerical order from 1 to 4.

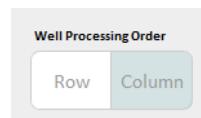

Selecting Plate Wells and Removing Assignments

You must select plate wells to assign calibrators or samples to them. After you have made assignments, you can deselect or select wells without changing their assignments.

- To select a block of wells, touch the plate and draw a marquee around the wells.
- To select individual, nonadjacent wells, turn on **Multi Select** > touch the wells.

- To remove specific selections, touch the plate and draw a marquee around the wells > touch **Clear Selected**.
- To remove all selections and assignments on a plate, touch **Clear Rack**.
- To remove assignments from specific wells, select the wells > touch **Clear Contents**.

Understanding the Required Sample Volume Field

The Required Sample Volume field reports the total amount of sample that will be consumed by all assignments that you have made. This field can help you monitor sample consumption so that you do not assign more sample volume than you have placed in a plate or a tube.


Understanding the Autoname Tool

Use the Autoname tool to quickly name a range of plated samples. This tool creates names of the following format:

[alphanumeric prefix up to nine characters long][first number in series]

For example: PSA02Mar1, PSA02Mar2, PSA02Mar3

Autonames are applied to a set of selected samples in row order or column order.

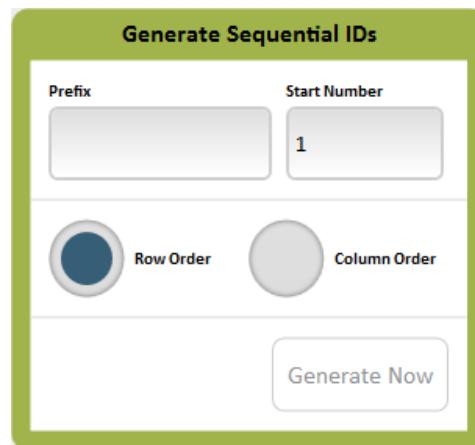

❖ To generate sequential IDs

- 1 Enter a text or numeric prefix in the Prefix field.
- 2 Touch **Row Order** or **Column Order**.
- 3 Optionally, change the Start Number.
- 4 Touch **Generate Now**.


Scanning Sample Barcodes

When loading samples in tube racks, you can use the handheld or on-board scanners to enter sample barcode information.

Note that if you are using the on-board scanner, you must assign all sample properties (including assay definition, replicates, and dilutions) prior to loading the barcoded sample tube rack. These properties cannot be modified once you insert the rack.

❖ To scan samples with the handheld scanner

- In the **Setup Run tab > Assign Samples** screen, select the lane into which you plan to load the sample tube rack.

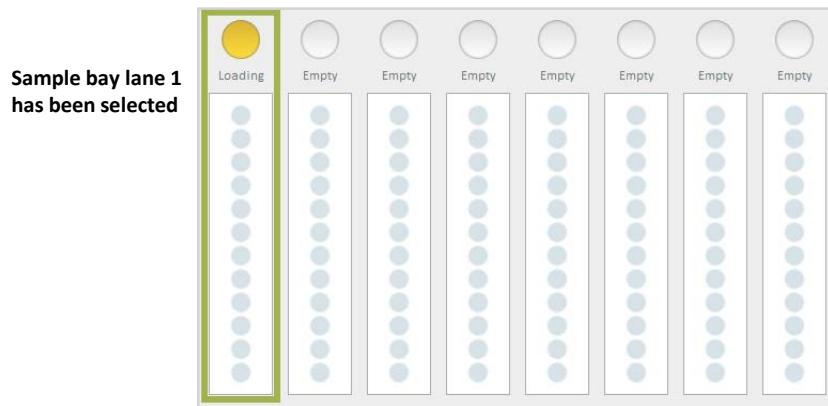

- Touch the position of the sample tube in the tube rack diagram.



- Scan the sample barcode with the handheld scanner.

❖ To scan samples with the on-board scanner

- In the **Setup Run tab > Assign Samples** screen, turn on the **Enable On-Board Scanner** toggle.
- Select the lane into which you plan to insert the tube rack.

3. Wait for the on-board scanner to focus its laser on the lane. You will see a notification informing you that the laser is focusing. Do not navigate to a new tab while the scanner is focusing, as this will halt the process.



Barcode scanner is focusing...

4. Insert the rack into the sample bay lane that you selected. The on-board scanner will read the sample barcodes as you insert the rack.

Assigning Calibrators

To assign calibrators, you identify in the software the physical location of all calibrators in a plate or tube rack, and then you associate each calibrator with a calibrator concentration in the assay definition.




Note! You can import a saved map of calibrator and sample assignments. For information about this option, see “Importing and Exporting a Layout” on page 103. The Import button appears on the Assign Calibrators screen. The Export button appears on the Assign Samples screen.

Assigning Calibrators in Plates

Make sure that you have plated the calibrators correctly. You can assign multiple calibrators to a plate up to the capacity of the plate well. See “Preparing Reagents, Calibrators, Controls, and SamplesError! Reference source not found.” on page 68.

- 1 Open the Assign Calibrators screen:

- If you are in the Assign Samples screen, touch **Assign Calibrators**.
- If you are not in the Assign Samples screen, touch the Setup Run tab > touch **Assign Calibrators**.

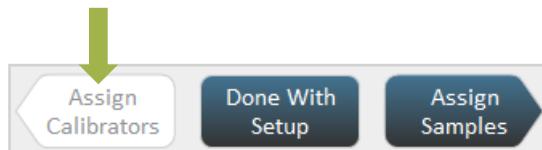

- 2 Recommended: Place prepared plates of calibrators and samples in plate racks and put the racks into the sample bay, but do not fully insert them.



Note! If you choose to insert the plate racks after you finish assigning calibrators and samples, **you must make sure** that the plate containing the calibrators is loaded in the position that you identify for it in the next step of this procedure. Placing the racks into the sample bay now makes it easier to correctly identify plate position.

3 Name the plate that contains the calibrators:

- a Touch the plate position that matches the location of the plate that contains the calibrators. For information on reading the plate position diagram, see page 85.



- b Enter a name for the batch and a plate name > press **Enter**. You can use any combination of letters and numbers.




4 Perform the following steps to associate each set of calibrator replicates with the calibrator values in the assay definition, in ascending or descending order.

- a Identify the first of the plated calibrators:


- ◆ **If you have plated only one set of calibrators:** Touch the well that contains the first calibrator in the set.
- ◆ **If you have plated two or more replicates:** Touch the first calibrator for one of the replicates. (You will repeat the entire assignment procedure to assign the remaining replicates.)

b Identify the assay:


- ♦ **If you have already loaded reagents:** Set the Assays Loaded only button to **On** so that the screen only displays assay definitions that contain the loaded reagents. Select your assay.

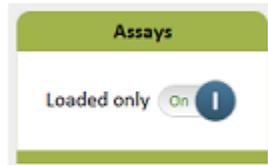

c Touch the highest or lowest calibrator > touch **OK**. The well that you selected turns red and a checkmark appears on it, signifying that an assay has been assigned to the well.

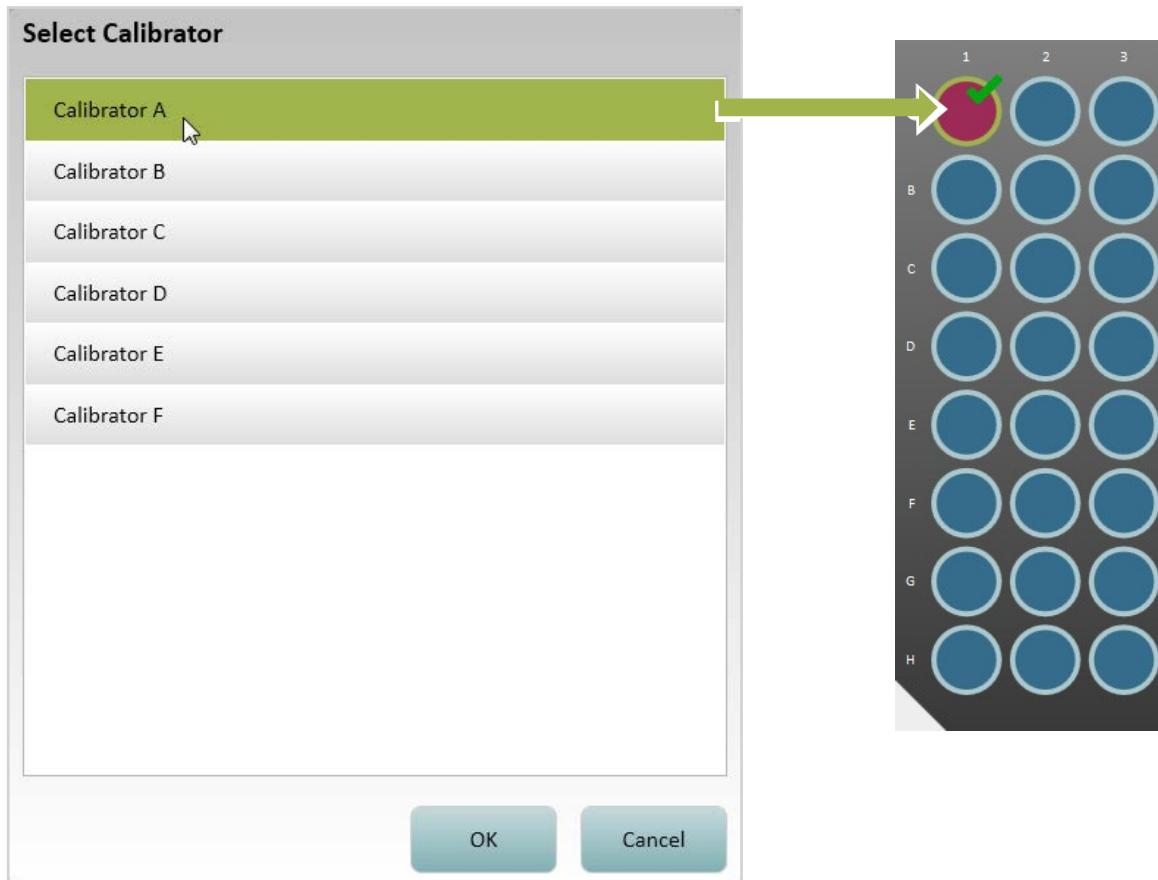


- d Touch **Ascending** or **Descending** under Assign All Calibrators. The software assigns all calibrators that are defined for the assay in the specified order, beginning with the selected well. The plate diagram shows the assignments. To view the concentration that is assigned to a well, touch the well and view the Calibrator box.
- e Repeat to assign all sets of calibrator replicates.

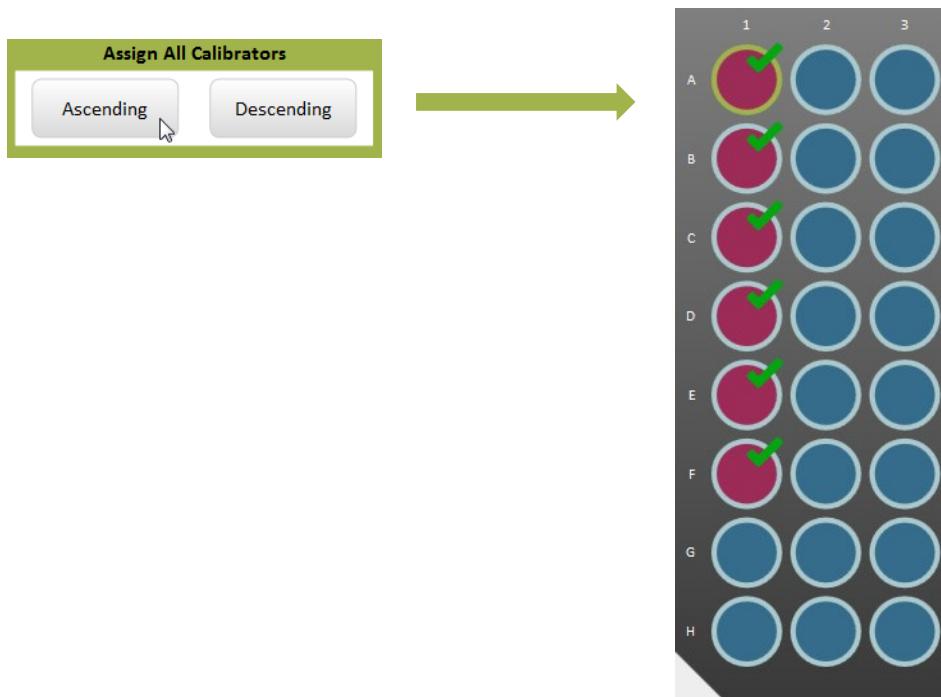

- 5 Select the processing order for the plates. This is a batch-wide setting that applies to all jobs in the run, even if on separate plates.
 - Select **Row** to process the calibrators from left to right
 - Select **Column** to process the calibrators from top to bottom.


If you have multiple plates, they are processed in ascending numerical order from 1 to 4.

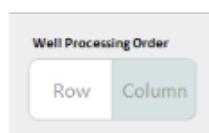

- 6 If the plates contain samples, assign them now. See “Assigning Samples in Plates” on page 88. If you are finished, see “Loading Calibrators/Samples into the Instrument” on page 102.

Assigning Calibrators in Tubes

- 1 Open the Assign Calibratorsscreen:
 - If you are in the Assign Samples screen, touch **Assign Calibrators**.
 - If you are not in the Assign Samples screen, touch the Setup Run tab > touch **Assign Calibrators**.

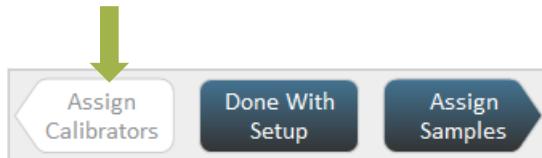


- 2 Recommended: Place tubes of prepared calibrators and samples in tube racks and put the racks into the sample bay, but do not fully insert them.




Note! If you choose to insert the tube racks after you finish assigning calibrators and samples, **you must make sure** that tubes containing the calibrators are loaded in the positions that you identify for them in the next step of this procedure. Placing the racks into the sample bay now makes it easier to correctly identify tube positions.

- 3 In the tube position diagram, touch the lane that contains the calibrators.

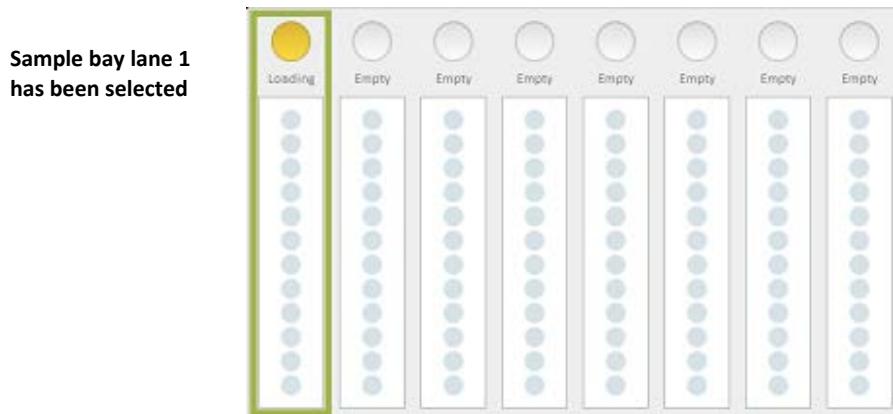

- 4 Enter a batch name and a tube rack name > touch **Enter**. You can use any combination of letters and numbers.

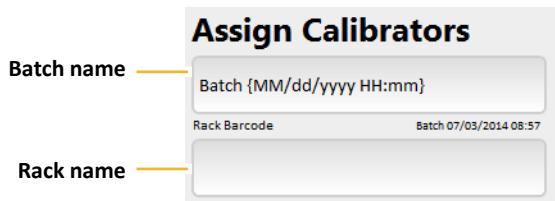


- 5 Perform the following steps to associate tubes of calibrators with the calibrator values in the assay definition, in ascending or descending order.

- a In the tube rack diagram, touch the position that holds the first calibrator in the rack.




- b Identify the assay:

- If you have already loaded reagents: Set the Assays Loaded only button to **On** so that the screen only displays assay definitions that contain the loaded reagents. Select your assay.

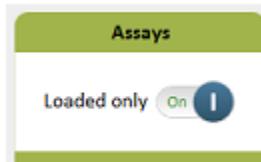


- c In the Select Calibrator dialog, touch the position that contains the highest or lowest calibrator > touch **OK**.

In the rack position diagram, the tube position that you selected turns red and a checkmark appears on it, signifying that an assay has been assigned to the tube.

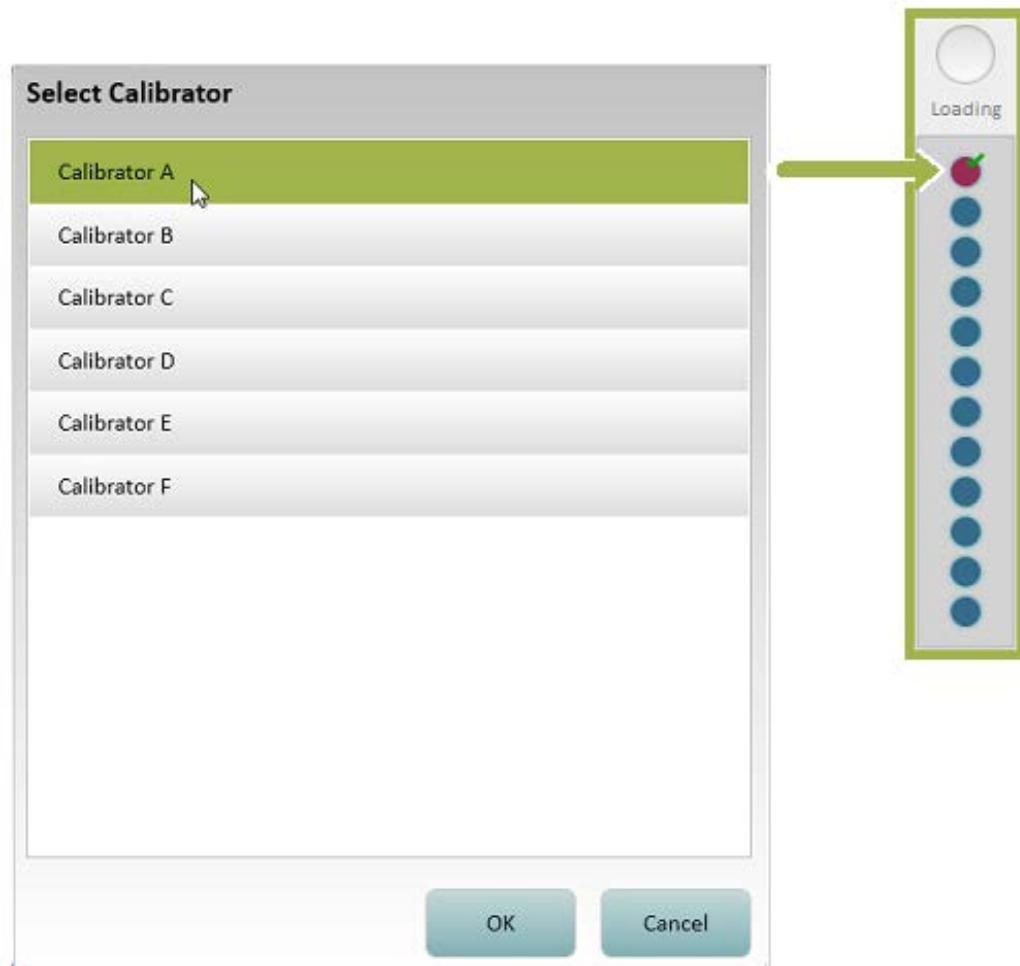

- d Touch **Ascending** or **Descending** under Assign All Calibrators. The software assigns all calibrators that are defined for the assay in the specified order, beginning with the selected tube position. The tube position diagram shows the assignments.

To view a calibrator assignment, touch a tube in the tube rack diagram. The calibrator and concentration appear in the Calibrator box.

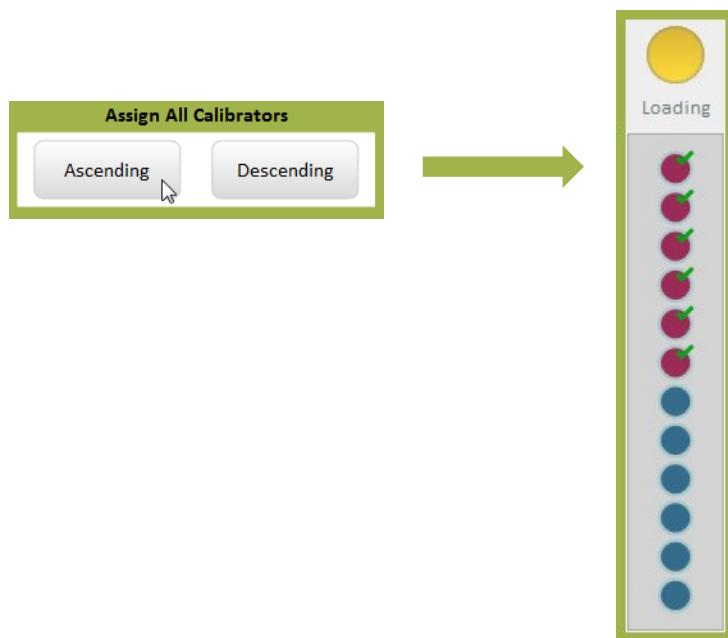

- 6 If you have loaded tubes of samples, assign them now. See “Assigning Samples in Tubes” on page 100. If you are finished, see “Loading Calibrators/Samples into the Instrument” on page 102.

Assigning Samples

To assign samples, you identify the physical location of all samples in a plate or tube rack, and then you specify the number of replicates to run for each sample and dilutions (if the assay definition allows).




Note! You can also import a saved layout of sample (and calibrator) assignments. For information about this option, see “Importing and Exporting a Layout” on page 103.

Assigning Samples in Plates

1 Open the Assign Samples screen:

- If you are in the Assign Calibrators screen, touch **Assign Samples**.
- If you are not in the Assign Calibrators screen, touch the Setup Run tab > touch **Assign Samples**.

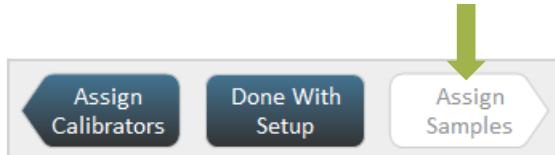

2 Recommended: Place prepared plates of calibrators and samples in plate racks and put the racks into the sample bay, but do not fully insert them.




Note! If you choose to insert the plate racks after you finish assigning calibrators and samples, **you must make sure** that plates containing the samples are loaded in the positions that you identify for them in the next step of this procedure. Placing the racks into the sample bay now makes it easier to correctly identify plate positions.

3 In the plate position diagram, touch a plate that contains samples. (For information on reading the plate position diagram, see page 85.)




- a** If a name appears in the Plate Barcode field, the plate is already named, and you can skip to step 4. This is true if you touch a plate that contains calibrators that you have previously assigned.

- b** If the Plate Barcode field is empty, enter a batch name and a plate name > touch **Enter**. You can use any combination of letters and numbers.

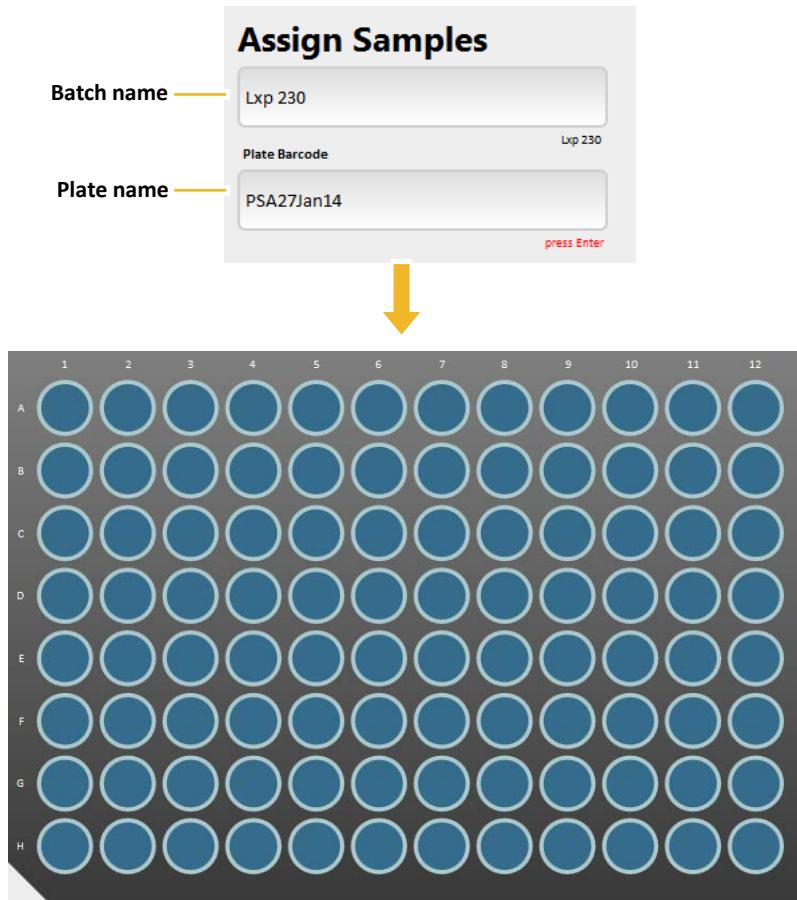


- 4** Select the processing order for the plates. This is a batch-wide setting that applies to all jobs in the run, even if on separate plates.
- Select **Row** to process the samples from left to right.
 - Select **Column** to process the samples from top to bottom.

If you have multiple plates, they are processed in ascending numerical order from 1 to 4.

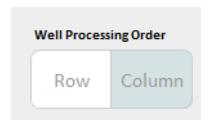


- 5** Assign an ID to the samples in the selected plate (see page 86 for information on autonames):
- Select a range of wells (for information on selecting wells, see page 85).

- b** In the Autoname tool, enter a prefix string, which will apply to every name, and a starting number.
- c** Touch **Row Order** or **Column Order**.
- d** Touch **Generate Now**. The software assigns a name to the selected wells and they turn yellow.

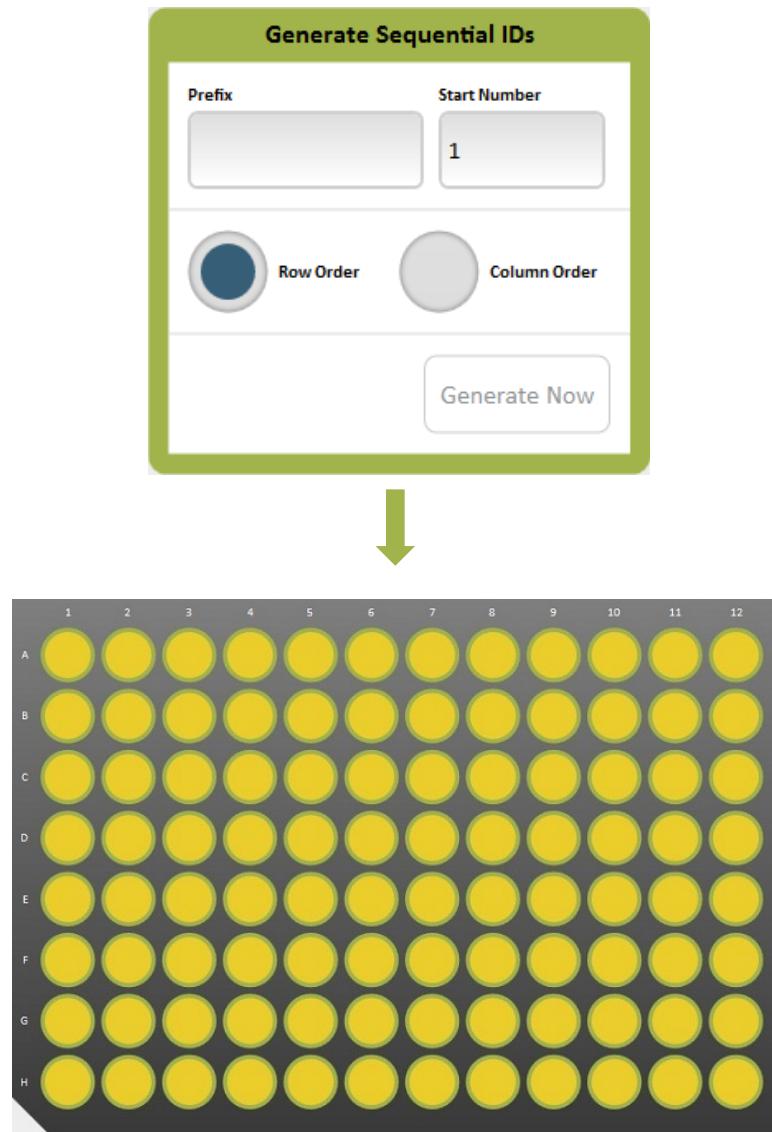


- 6** Specify an assay, replicates, and dilutions for samples:
 - a** Select plate wells.
 - b** Specify an assay (checkmarks appear on the selected wells after you do this):

- ♦ **If you have already loaded reagents:** Set the Assays Loaded only button to **On** so that the screen only displays assay definitions that contain the loaded reagents. Select your assay.

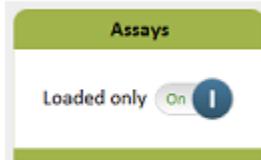

c Specify replicates and dilutions:

- ♦ Select a replicate number from the drop-down list.
- ♦ Touch the dilution that is specified in the assay kit instructions.
- ♦ Touch any other dilution that has a checkmark to remove the checkmark.

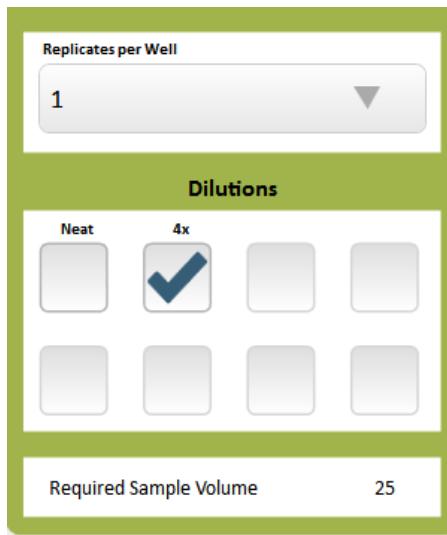


- 7** If you have more than one plate of samples, repeat steps 3 to 5 to assign all the samples.

- 8 If you need to assign calibrators, see “Assigning Calibrators in Plates” on page 88. If you are finished, see “Loading Calibrators/Samples into the Instrument” on page 102.

Assigning Samples in Tubes

- 1 Open the Assign Samples screen:
 - If you are in the Assign Calibrators screen, touch **Assign Samples**.
 - If you are not in the Assign Calibrators screen, touch the Setup Run tab > touch **Assign Samples**.

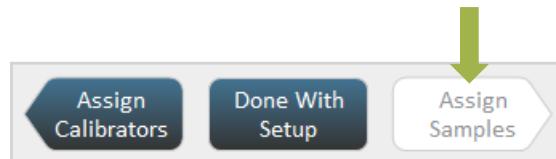


- 2 Recommended: Place prepared tubes of calibrators and samples in tube racks and put the racks into the sample bay, but do not fully insert them.



Note! If you choose to insert the tube racks after you finish assigning calibrators and samples, **you must make sure** that tubes containing the samples are loaded in the positions that you identify for them in the next step of this procedure. Placing the racks into the sample bay now makes it easier to correctly identify tube positions.

- 3 Touch a lane that holds a rack containing sample tubes.


Sample bay lane 1
has been selected

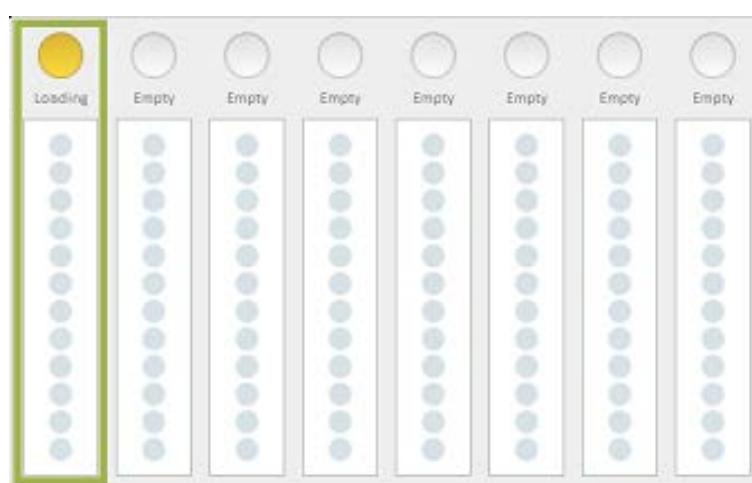

- 4** Assign an ID to the samples in the selected rack (see page 86 for information on autonaming):
- In the tube rack diagram, select a tube position.
-
- Tube 1 has been selected
- Enter a Barcode ID for the tube. This step is not necessary if you already assigned an ID to the samples, as specified in Step 4 above.




- Physically load the rack. See “Loading Calibrators/Samples into the Instrument” on page 102.
 - Repeat these steps to name all tubes in the rack.
- 5** Select an assay for the selected tube:
- If you have already loaded reagents: Set the Assays Loaded only button to On so that the screen only displays assay definitions that contain the loaded reagents. Select your assay.

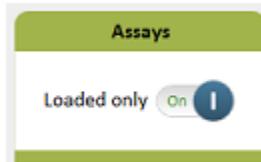

- 6** Specify replicates and dilutions for the selected tube:
- Select a replicate number from the drop-down list.
 - Touch the dilution that is specified in the assay kit instructions.

- c Touch any other dilution that has a checkmark to remove the checkmark.

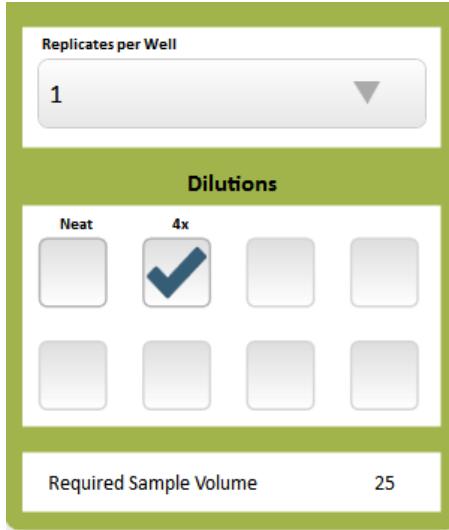

- 7 Repeat steps 3 to 6 to assign all tubes of samples.
- 8 If you need to assign calibrators, assign them now. See “Assigning Calibrators in Tubes” on page 92. If you are finished, see “Loading Calibrators/Samples into the Instrument” on page 102.

Loading Calibrators/Samples into the Instrument

The final step in assigning calibrators and samples is to load them into the instrument so that the software accepts the assignments.

❖ To load assigned calibrators and samples

- 1 Optional: Save the layout. In the Assigned Samples screen, touch **Export Layout** > save the file in the Save As dialog.
- 2 If you have not already done so, place the plate racks or tube racks into the sample bay.



Note! Place the racks in the sample bay according to the assignments that you have made. If you do not put plates or tubes in the positions that you identified for them, the instrument will not be able to locate the calibrators or samples correctly, and your assay results will be undefined.

- 3 Push the plate racks or tube racks into the sample bay until they latch into place.
- 4 Recommended: Touch **List View** to see a list of everything you have loaded. This allows you to check for the correct number of replicates, dilutions, and other consumables.

- 5 Touch **Done With Setup**. The All Resources screen in the System Resources tab opens. See “Starting a Run” on page 103 to initiate the run.



Importing and Exporting a Layout

You can export a finished layout of calibrator and sample assignments to a file and then re-use it the next time you run the assay.

❖ To export a layout

In the Assign Samples screen, touch **Export Layout** > name and save the file. Do not change the specified file extension (.csv).



A rounded rectangular button with a light blue gradient background and white text: "Export Layout".

❖ To import a layout

In the Assign Calibrators screen, touch **Import Layout** > open the file.



A rounded rectangular button with a light blue gradient background and white text: "Import Layout".

Starting a Run

You can begin a run immediately after you assign calibrators and samples if sufficient system resources are available and no mandatory maintenance tasks are pending.

When you finish assigning calibrators and samples, the software calculates the system resources needed to run the assay and posts resource status on the System Resources screen. If enough consumables are available to run the loaded calibrators and samples, the supplies of wash buffer are sufficient, the liquid waste container is not full, and no maintenance tasks are pending, you can start the run.

❖ **To start a run**

- 1 Touch **Done With Setup** in the Assign Calibrators or Assign Samples screen. The System Resources tab opens.
- 2 If the button is active, touch **Start Run**.



Start Run


If the button is inactive, touch **Resource Details** to display a complete readiness status (including status of beads, detectors, SBG, and RGP). You may need to perform some of the following tasks to activate the **Start Run** button again:

- If the Maintenance task icon is red, touch the Maintenance icon and perform any pending mandatory tasks. See Chapter 10, “Maintaining the Simoa HD-X Analyzer,” on page 127.
- If the cuvettes, tips, or discs diagrams are red, add more of those resources. See “Loading Instrument Consumables” on page 38.
- If the system liquid diagram or either wash buffer diagrams are red, refill the container. See “Refilling the System Fluid Containers” on page 48.
- If the liquid waste diagram is red, empty the liquid waste container. See “Emptying the Liquid Waste Containers” on page 55.
- If the solid waste diagram is red, empty the solid waste container. See “Emptying the Solid Waste Container” on page 53.
- If the sealing oil diagram is red, replace the sealing oil bag. See “Changing the Sealing Oil Bag” on page 51.

Understanding the System Resources Screen


Diagram color indicates status:

- **Green:** Adequate amounts
- **Yellow:** Low amounts
- **Red:** Empty, must be refilled, or container is full and must be emptied.

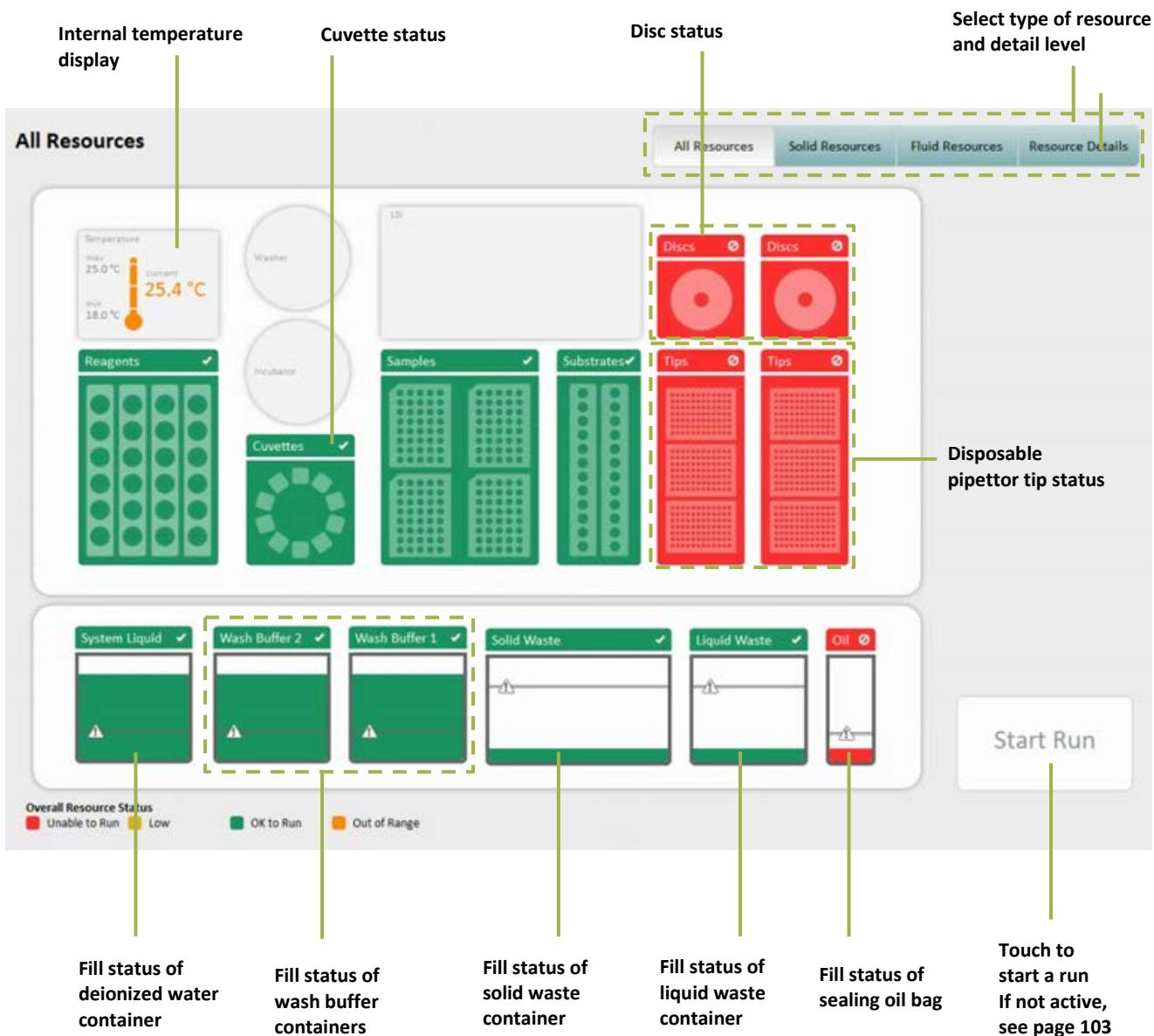

Understanding the Solid Resources Screen

Touch the System Resources tab > touch the cuvette diagram or one of the disposable tip or disc diagrams. The Solid Resources screen opens.

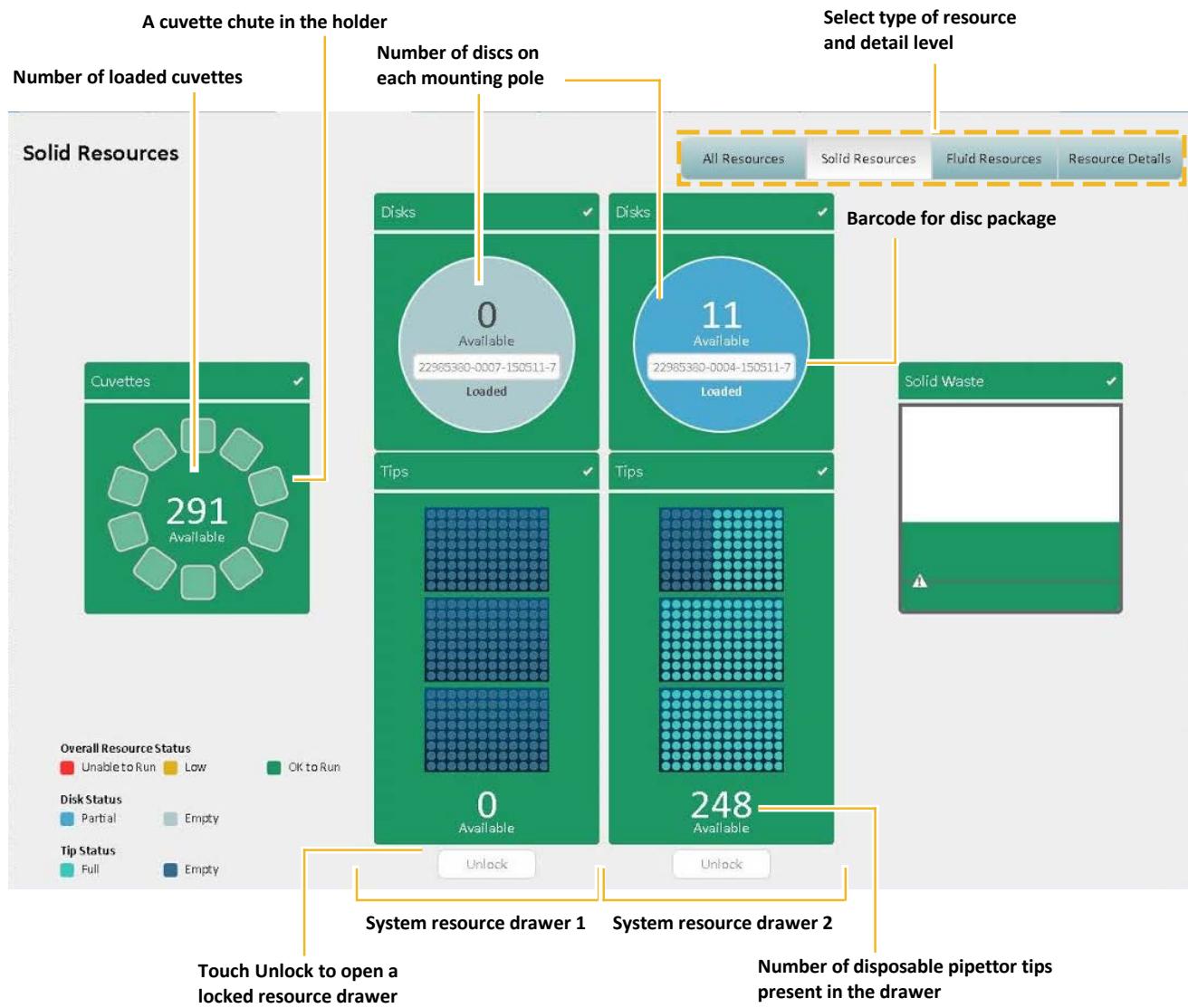

Monitoring a Run


Touch the Current Run tab to observe the status of sample processing during the run. Each loaded plate or tube is graphically represented on this tab. As samples are processed, they change color. See the color key at the bottom of the screen for an explanation of the meaning of each color.

Sample In-Process Run screens for plates and tubes appear below.




CAUTION! Do not attempt to open the system resource drawers; remove or insert a reagent rack, a sample rack, or an RGP rack; or change the software settings once the run has begun. Doing any of these things may cause the instrument to halt the run, with subsequent loss of results. You can refill system liquids and empty one liquid waste container at a time during the run.

Understanding the In-Process Run Screen for Plates

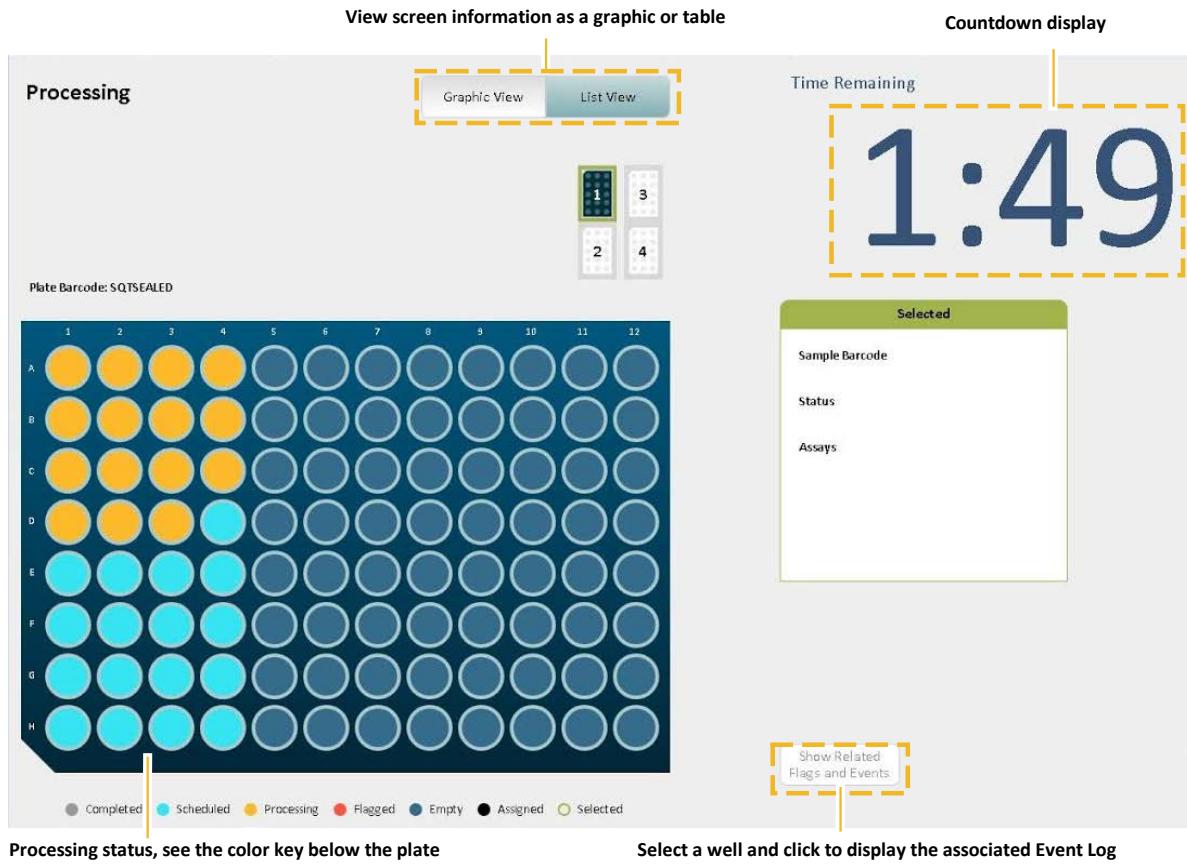

Understanding the In-Process Run Screen for Tubes

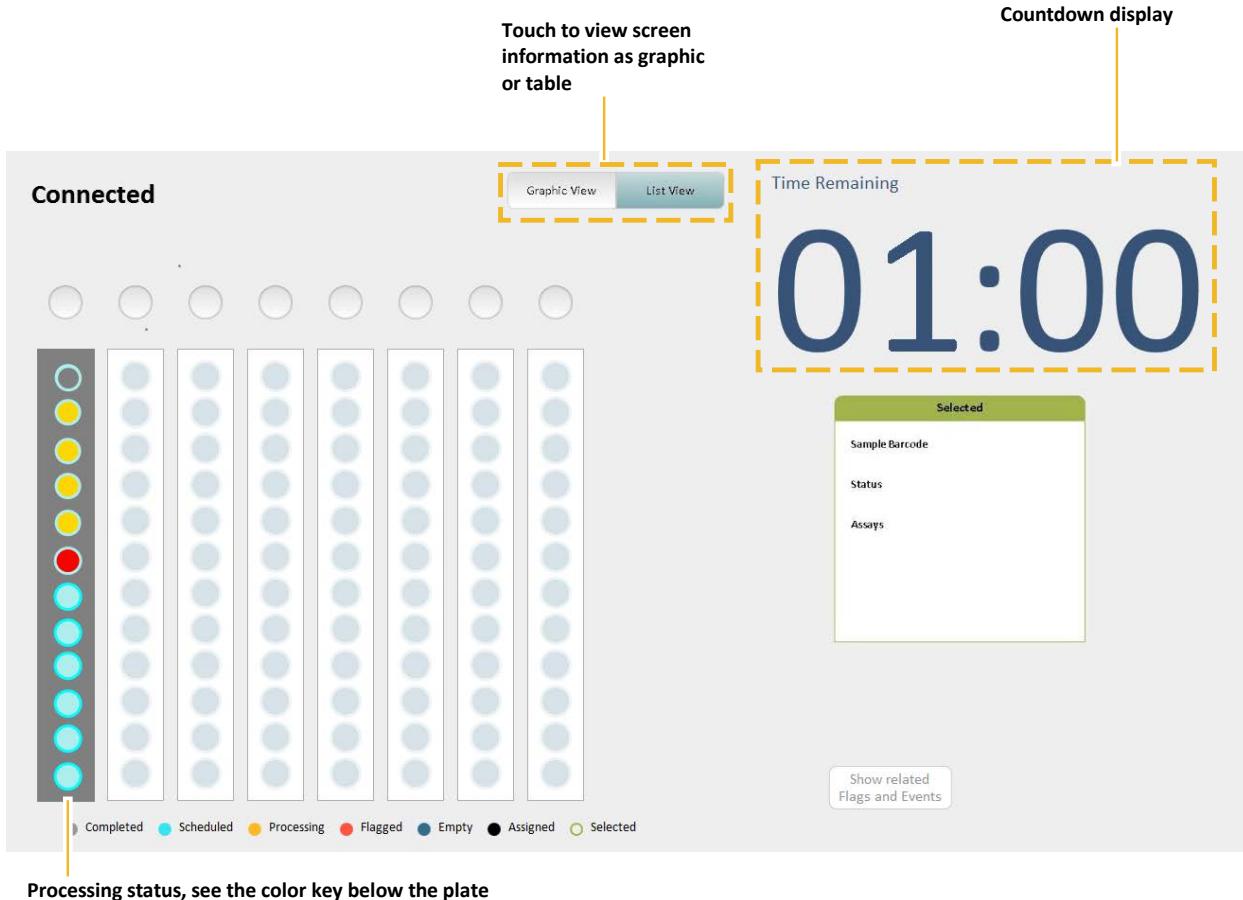

Checking Instrument Status

During the run, periodically check the instrument status indicators at the bottom of the software to determine whether you need to refill any system liquids or empty any system waste containers. The icons for system fluids, consumables, and instrument waste containers change color when any of these components need attention:

- **Green:** Adequate amounts
- **Yellow:** Low amounts
- **Red:** Empty, must be refilled, or container is full and must be emptied

For information about refilling system fluid containers or emptying system waste containers, see Chapter 10, “Maintaining the Simoa HD-X Analyzer,” on page 127.

When the run finishes, the countdown display says 00:00 and the instrument state notice changes to Ready.



Viewing the Event Log

Samples that cannot be processed for some reason are flagged red in the Current Run tab. To determine the reason for the flag, you can view the Event Log.

❖ To view the Event Log

Touch the History & Reports tab > touch **Event Log**. For more information, see Understanding the Event Log Screen on page 122.

8 Analyzing Run Data

Data Analysis Overview.....	111
Working with Run History	111
Exporting Results in CSV Format	114
Exporting Results in PDF or Microsoft Excel Format	115
Analyzing Results in the Data Reduction Tab.....	115

Data Analysis Overview

When the assay run is complete, results appear on the Run History screen. You can perform data analysis on all of the results or a subset.

Results can be exported to an external analysis package or a LIMS system, or you can view and analyze results in the instrument's software. The following options are available:

- Export a CSV file for external analysis
- Export a PDF or Excel file for external analysis
- Analyze the results in the instrument's Data Reduction screen

Working with Run History

The Run History screen provides options that allow you to select the data and columns to include in the analysis, export the data to a CSV file, and move analysis data between multiple instruments.

The Run History view calculates replicate statistics only for samples (including calibrators) that come from the same well or tube (for each group of replicates that originate from the same location).

The screenshot shows the Run History screen with the following numbered elements:

- 1**: Assay Contains: PSA (with X and + buttons)
- 2**: Configure Columns
- 3**: Number of Selected Results: 0 out of 1092
- 4**: Select all Results
- 5**: Deselect all Selected Results
- 6**: Exclude Selected Results from Analysis
- 7**: Include Selected Results into Analysis
- 8**: Show Related Flags and Events
- 9**: Recalculate with Different Curve
- 10**: Export
- 11**: Archive/Restore

The main table displays the following data:

Selected	Instrument SN	SW Version	Sample Barcode	Errors	Assay	Plex	Location	Carrier Barcode	Replicate AEB	Mean AEB	SD AEB	CV AEB	Replicate Conc.	Mean Conc.	SD Conc.	CV
	2710000024	1.4.1502.5001	PSA Calibrator A		PSA	PSA	Lane: 1 - Well: 1	20150518PSA	0.01				0			
	2710000024	1.4.1502.5001	PSA Calibrator A		PSA	PSA	Lane: 1 - Well: 1	20150518PSA	0.008				0			
	2710000024	1.4.1502.5001	PSA Calibrator A		PSA	PSA	Lane: 1 - Well: 1	20150518PSA	0.009				0			
	2710000024	1.4.1502.5001	PSA Calibrator A		PSA	PSA	Lane: 1 - Well: 1	20150518PSA	0.009				0			
	2710000024	1.4.1502.5001	PSA Calibrator A		PSA	PSA	Lane: 1 - Well: 1	20150518PSA	0.008				0			
	2710000024	1.4.1502.5001	PSA Calibrator A		PSA	PSA	Lane: 1 - Well: 1	20150518PSA	0.009				0			
	2710000024	1.4.1502.5001	PSA Calibrator A		PSA	PSA	Lane: 1 - Well: 1	20150518PSA		0.009	0.001	0.105	0	0		
	2710000024	1.4.1502.5001	PSA Calibrator B		PSA	PSA	Lane: 1 - Well: 2	20150518PSA	0.037				0.1			
	2710000024	1.4.1502.5001	PSA Calibrator B		PSA	PSA	Lane: 1 - Well: 2	20150518PSA	0.046				0.1			
	2710000024	1.4.1502.5001	PSA Calibrator B		PSA	PSA	Lane: 1 - Well: 2	20150518PSA	0.039				0.1			

The following table explains how to use the Run History screen.

Element	Definition
1 Add New Filter 	Adds filter criteria that determine what results appear in the table. Multiple filters are allowed. Touch the Add New Filter (+) button to add a filter. Touch the X on the filter button to remove the filter from the table. Loading a large amount of data into the results table can decrease responsiveness. To avoid this, add a new filter before removing unwanted ones.
2 Configure Columns	Opens the Configure Columns window that allows you to configure column display and ordering.


Element	Definition
3 Automatic Replicates Selection	<p>When enabled, the software automatically selects all related replicates when one of the replicate results is selected.</p> <p>Automatic Replicates Selection </p>
4 Select all Results	Selects all currently displayed results, as determined by data filters, if any are in place.
5 Deselect all Selected Results	Clears all selections made on results in the table.
6 Exclude Selected Results from Analysis	<p>Marks selected results for exclusion in analysis. Only results that have been automatically included in analysis can be excluded manually. Isingle is not recalculated when excluding results in the Run History table.</p> <p>Note: Ensure the Automatic Replicates Selection button is turned off > select the replicates to be excluded from analysis, then touch Exclude Selected Results from Analysis button for result.</p>
7 Include Selected Results Into Analysis	Marks selected results for inclusion in analysis. Only results that have been excluded automatically from analysis can be included manually.
8 Show Related Flags and Events	Displays all flags and event messages for selected results.
9 Recalculate with Different Curve	Allows you to recalculate concentration results by manipulating the default calibration curve or to select a new calibration curve.
10 Export	Exports displayed results to a comma separated value (CSV) file. All columns displayed in the table are exported. If you have not selected a column with Configure Columns, it is not exported. Additionally, if you applied data filters to the table, only results matching the filter criteria are exported.
11 Archive/Restore	Provides ability to move result data between two computers running the Simoa HD-X software. Results archived (exported) to an XML file can be restored (imported) into the software running on another computer.

Exporting Results in CSV Format

Export a CSV file from the Run History screen to create a file for external analysis. To export data for troubleshooting analysis, include the Extended Properties column.

❖ To export results in CSV format

- 1 In the History & Reports tab, touch **RunHistory**.
- 2 Clear all filters on the results table.
- 3 Touch the **Add New Filter** button > select a filter from the drop-down list > touch **Next**.
- 4 Specify the filter criteria. Some filters provide drop-down lists and others provide a text entry field that sets the filter.
- 5 Touch **OK**. The results are filtered and the filter name appears above the results table.
- 6 Touch **Configure Columns**.
- 7 Move columns between the Available and Selected lists and change their order with the up and down arrow key buttons. The columns you select identify the data that is exported See *Appendix A – Configure Columns Definitions* for more information about each column.
- 8 Touch **Done**.

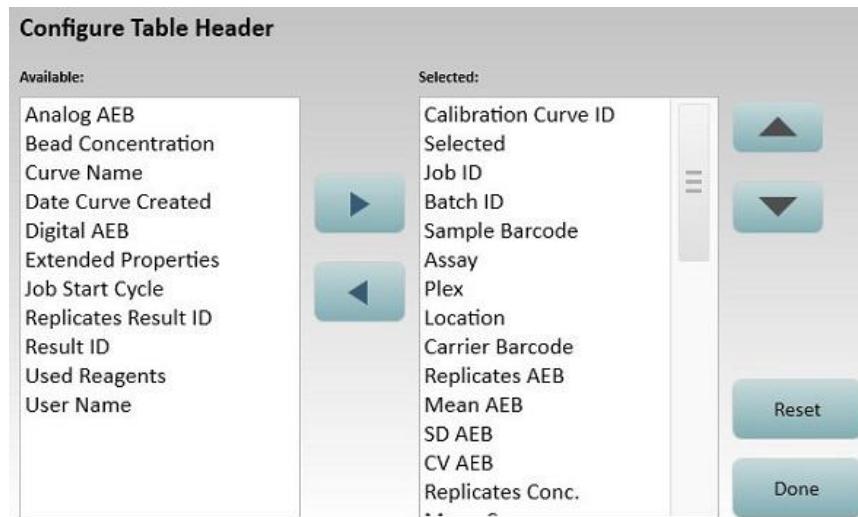

- 9 Touch **Select All**.

- 10 Touch **Export** > specify a file name > touch **Save**.

Exporting Results in PDF or Microsoft Excel Format

Export a file in PDF or Microsoft Excel format from the Reports window to create a file for external analysis.

❖ To export results in PDF or Microsoft Excel format

- 1 In the History & Reports tab, touch **Reports**.
- 2 Touch a report in the displayed list.
- 3 Select individual reports or touch **Select All** for all reports.
- 4 Touch **Preview Report** > touch **Export**.
- 5 Save the data in PDF or Excel format.

Analyzing Results in the Data Reduction Tab

The instrument plots a Calibration Curve representing the selected results on the Data Reduction tab.

Curves are batch-based. All the samples in the batch have the same calibration curve applied to them, regardless of processing order.


The calibration data table contains all calibration data points used in curve fitting. The table includes aggregated replicates data. If calibrators were run as replicates, their averages are used.

Curves are specific to the assay. All samples of the same assay and batch will initially have the same calibration curve, but they can be assigned to a different curve, if calibration data points are included in the batch or Curve Selection Strategy is set to None. Following is a sample Data Reduction tab and a table that explains its components.

❖ To review results in the Data Reduction tab

- 1 In the Data Reduction tab, select the relevant fields in Assay, Plex, Calibration Curve, Curve Fit Algorithm, and Weighting.
- 2 Use the tools described in the following table to analyze the data.

The following screen shows a calibration curve calculated with the default settings. Some default settings, such as fitting algorithm, come from plex definition.

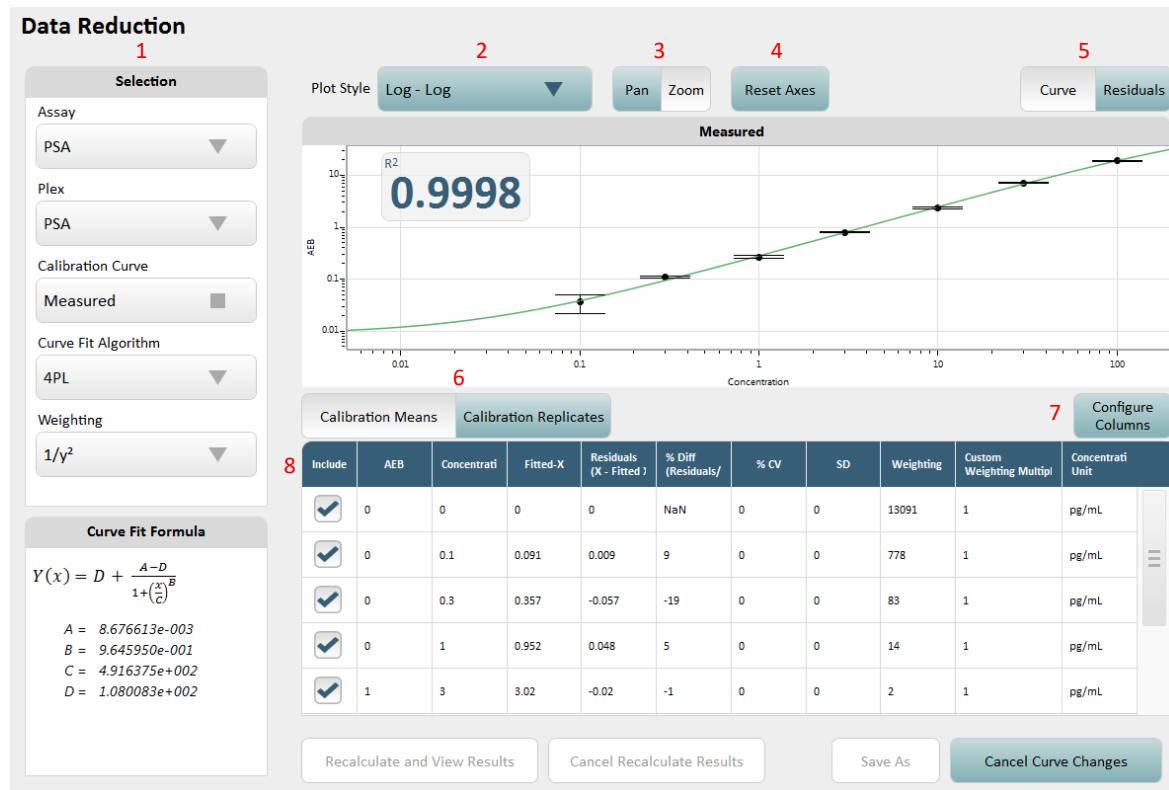

9

10

11

12

The following table explains how to use the Data Reduction tab.

Element	Definition		
1 Selection	Controls which calibration curve is displayed and the fitting algorithm and weighting used to generate the curve. Below are the recommended minimum number of data points needed for the various curve formulas.	Points Needed	Curve Formula
	2		Piecewise Linear, Linear, Exponential, Log-Log
	3		Quadratic
	4		Cubic, 4PL
	5		5PL

Element	Definition
2 Plot Style	Controls how the data is plotted.
3 Pan/Zoom	Pan across the graph or zoom in on a selected area for more detail.
4 Reset Axes	Use this to reset the graph axes to their original values.
5 Curve/Residuals	Use the Curve and Residuals buttons to toggle between the calibration curve and fitting errors.
6 Calibration Means/ Calibration Replicates	Calibration Means shows the replicates' summarized information in the main table. Calibration Replicates shows individual replicate concentrations.
7 Configure Columns	Select the columns to include in the calibration data table.
8 Include	To exclude individual replicate results from analysis, uncheck the Include checkbox in the row you want to exclude.
9 Recalculate and View Results	After selecting a different calibration curve, recalculates results and displays them in the Run History table.
10 Cancel Recalculate Results	Cancel recalculation and return to the History and Reports tab.
11 Save As	Saves the calibration curve. When you modify an existing curve, you must save it before it can be used to recalculate the results.
12 Cancel Curve Changes	Resets any changes made to the existing calibration curve.

9 Using Reports, Run History, and the Event Log

The History and Reports Tab.....	119
Understanding the Reports Screen.....	120
Understanding the Run History Screen.....	121
Understanding the Event Log Screen.....	122
Customizing the Table View in the History & Reports Tab	123
Creating a Report	125
Exporting or Printing a Report	126
Customizing the Logo for a Report	126


The History and Reports Tab

The History and Reports tab provides access to information about your assay runs from the following screens:

- Reports
- Run History
- Event Log

Understanding the Reports Screen

A report is a PDF file of run data that you can print or export to a computer. See “Creating a Report” on page 125 for an explanation of each report.

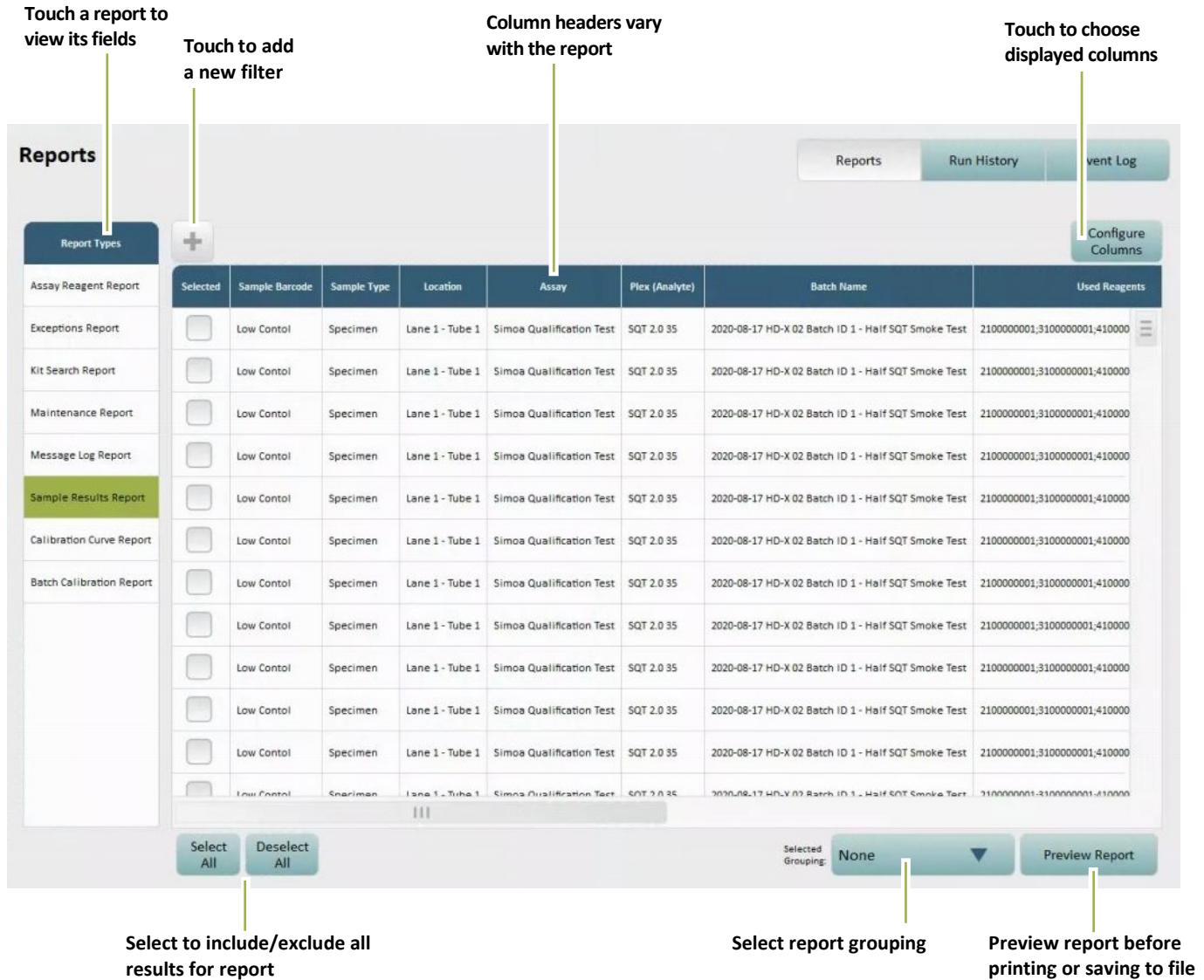

Understanding the Run History Screen

This tab is most frequently used to export the raw data for archival use or inclusion in the database. The Run History view calculates replicate statistics only for samples (including calibrators) that come from the same well or tube (for each group of replicates that originate from the same location).

Run History
Reports
Run History
Event Log

Number of decimal places controlled in System Settings
Choose columns to display

Assay Contains: psa

Filter results

Selected	Sample Barcode	Assay	Plex	Location	Carrier Barcode	Replicates AEB	Mean AEB	SD AEB	CV AEB	Replicates Conc.
<input checked="" type="checkbox"/>	PSA Calibrator A	PSA	PSA Kit Assay plex id=5	Lane: 1 - Well: 1	R1	0.0171243920206295				0
<input checked="" type="checkbox"/>	PSA Calibrator A	PSA	PSA Kit Assay plex id=5	Lane: 1 - Well: 1	R1	0.0200974777815112				0
<input checked="" type="checkbox"/>	PSA Calibrator A	PSA	PSA Kit Assay plex id=5	Lane: 1 - Well: 1	R1	0.0211416546448138				0
<input type="checkbox"/>	PSA Calibrator A	PSA	PSA Kit Assay plex id=5	Lane: 1 - Well: 1	R1		0.0194545081489848	0.00208438413186732	0.107141445874903	
<input type="checkbox"/>	PSA Calibrator B	PSA	PSA Kit Assay plex id=5	Lane: 1 - Well: 2	R1	0.0335304985516174				0.1
<input type="checkbox"/>	PSA Calibrator B	PSA	PSA Kit Assay plex id=5	Lane: 1 - Well: 2	R1	0.0331392657339118				0.1
<input type="checkbox"/>	PSA Calibrator B	PSA	PSA Kit Assay plex id=5	Lane: 1 - Well: 2	R1	0.0280805604070937				0.1

Automatic Replicates Selection (On)
Number of Selected Results: 6 out of 20
 L2DR Result
 Replicate Result
 Flagged Result

Select all Results
Deselect all Selected Results
Exclude Selected Results from Analysis
Include Selected Results into Analysis
Show Related Flags and Events
Recalculate with Different Curve
Export
Archive/Restore

Selects all currently displayed results
Automatically select all replicates when one is selected
Marks results for exclusion from analysis
Marks results for inclusion in analysis
Displays all event messages and flags
Recalculates concentration results with a new or modified curve
Exports displayed results to a CSV file
Archives to or restores batch data from an XML file

Clears selected results


Understanding the Event Log Screen

The Event Log tracks significant events that occur during software operation. Event log information is helpful for troubleshooting problems.

Icons in the Type column specify the severity of the message. These icons also appear at the bottom of every tab in the software. A number in the center of these icons signifies how many messages of that type are listed in the log.

To view the Event Log, touch the **History & Reports** tab > touch **Event Log**.


Using the Configure Columns button, you may customize the columns displayed (page 123). When exporting the Event Log, however, column headings are automatically selected to ensure the best export format (columns fit on a page). You may select either PDF (.pdf) or an Excel file (.xls). See “Exporting or Printing a Report” on page 126.

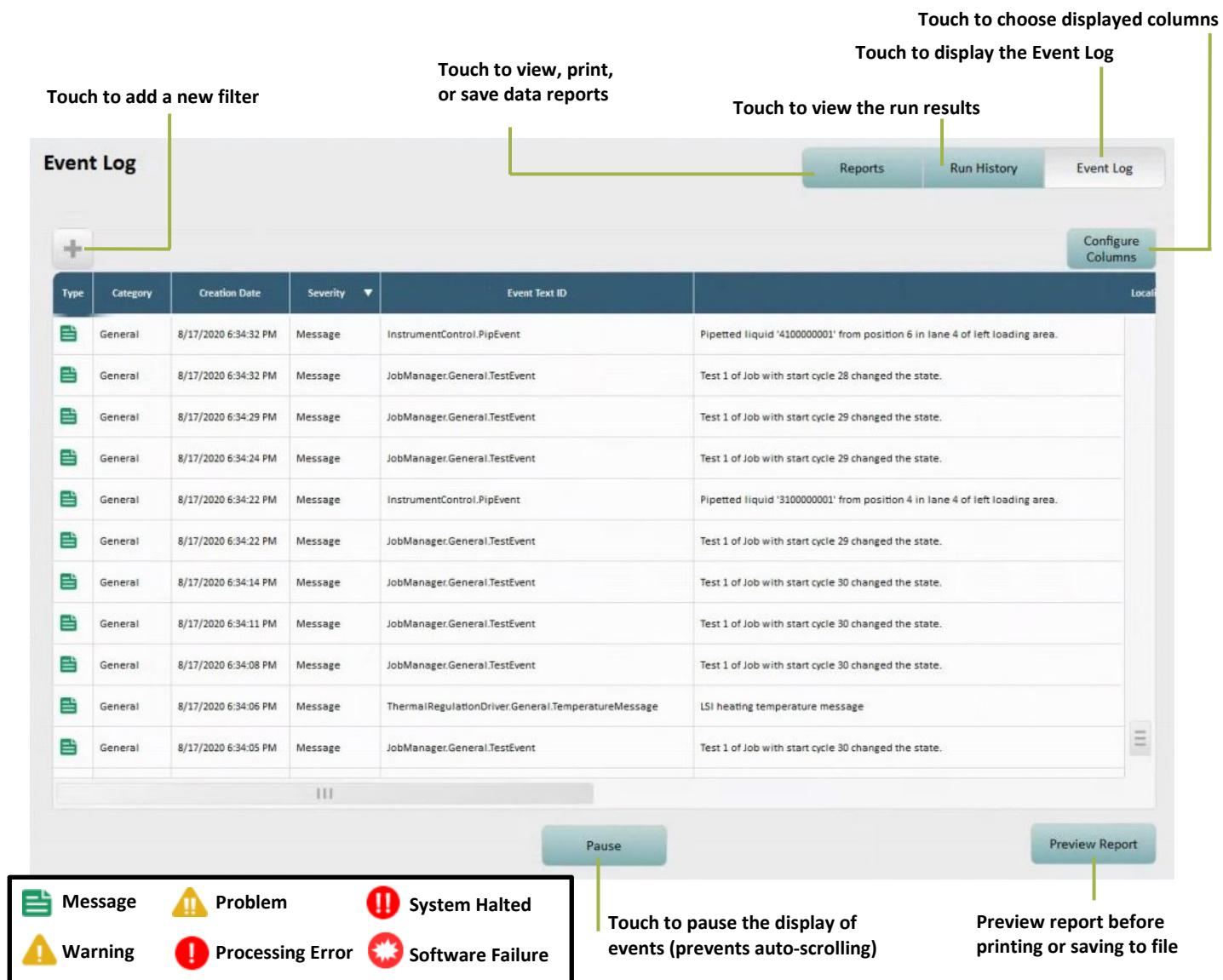

Customizing the Table View in the History & Reports Tab

Changing Columns

To change the columns for any table in the History & Report tab, touch **Configure Columns** > touch the arrows to move column names between the Available and the Selected boxes > touch **Done**.

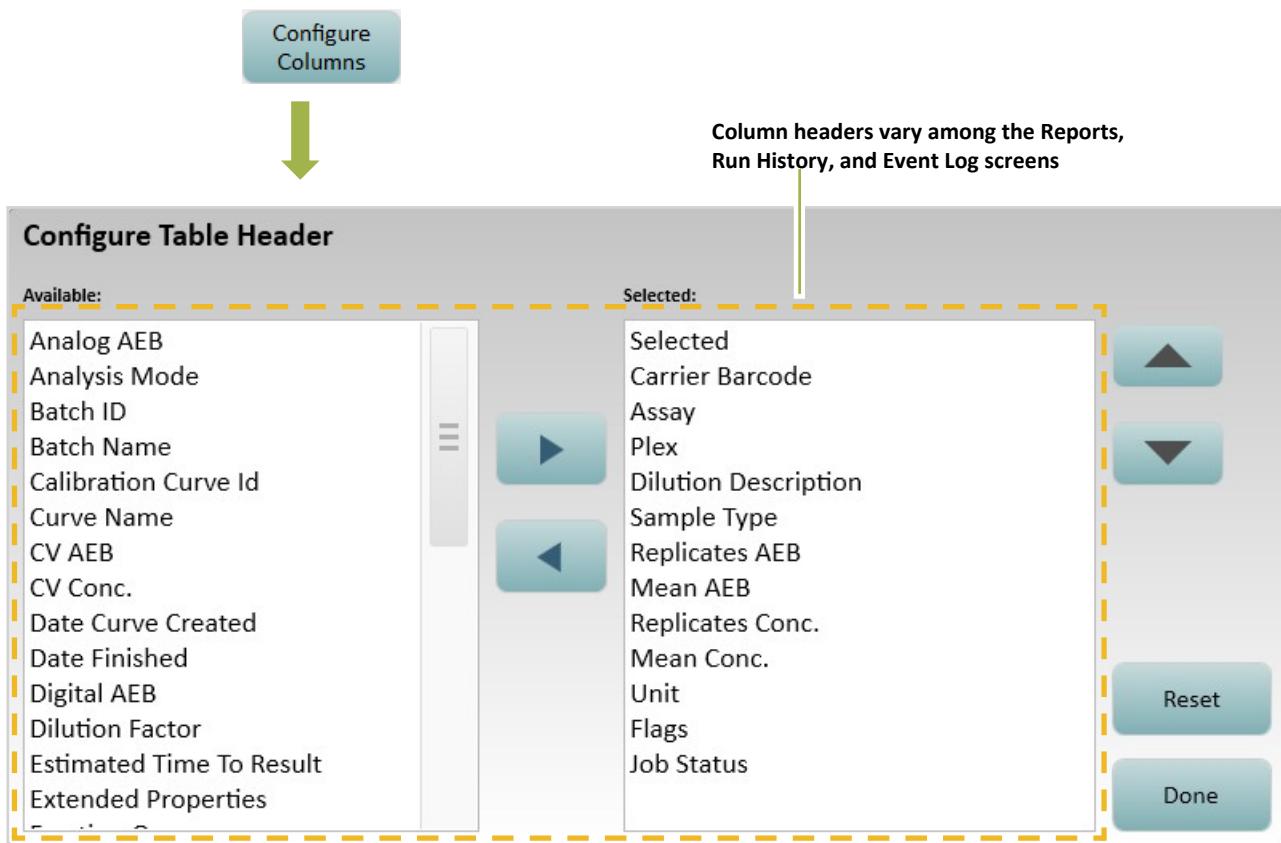

Changing Column Order

To move a column in any table in the History & Reports tab, use the up and down arrow keys, or touch the column header and drag it to the new position. In table view, you can also reorder the columns by dragging.

Sorting Rows

To sort any table in the History & Reports tab in ascending or descending alphabetical order for a specific column, touch the column header.


Adding a Filter

A filter modifies the contents of the Run History, the Event Log, or a report to show only data that matches the filter requirements. For example, you can filter the Event Log to show only warning messages.

You can add multiple filters to any table. Available filters vary with the table type.

❖ To add a filter

- 1 Touch the **Add New Filter** button > select a filter type from the drop-down list > touch **Next**.
- 2 Specify the filter criteria. Some filters provide drop-down lists for this, and others provide a text entry field in which you enter terms that set the filter.
- 3 Touch **OK**. The filter appears above the table.

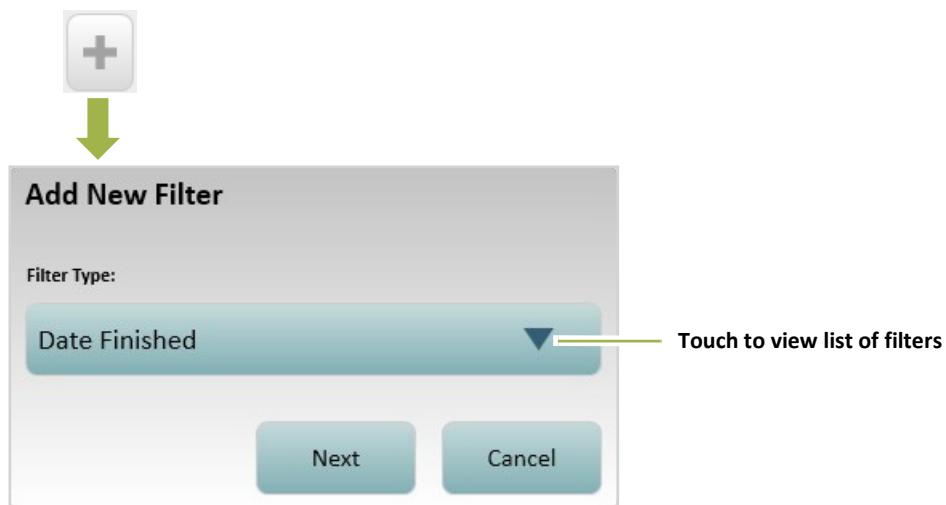

❖ To remove a filter

Touch the large X in the filter.



Creating a Report

A report is a templated summary of run data that can be exported to a user-specified location. The following types of reports are available:


Assay Reagent Report	Reagents used in selected assay runs
Exceptions Report	List of exception (error) messages occurring during the selected runs
Kit Search Report	List of samples associated with a specific kit reagent lot
Maintenance Report	List of completed instrument maintenance tasks
Message Log Report	List of all messages generated by the software during a specified selected time interval
Sample Results Report	Results for selected assay runs
Calibration Curve Report	Calibration curve information for selected run (fitting algorithm, weighting factor, style, etc.)
Batch Calibration Report	Includes calibration curves for all assays and all samples processed for a selected batch. Organized by assay plexes.
Clinical Results Report	Results for samples run with authorized clinical assays only.

To customize reports, change the report table (see “Customizing the Table View in the History & Reports Tab” on page 123.)

❖ To create a report

- 1 Touch the History & Reports tab > touch **Reports**.

- 2 Touch a report type.

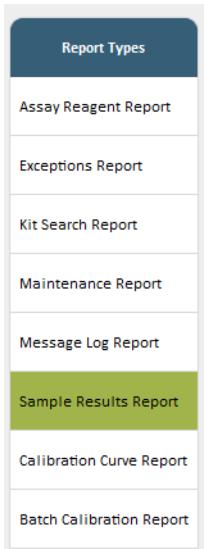

Exporting or Printing a Report

You can preview a report, export it to a file, or print it.

❖ To export or print a report

- 1 Touch **Preview Report** to view the formatted report. After a few moments, the report appears.
- 2 Touch **Export or Print**.
 - **Export:** Opens a Save As window. Enter a filename > select a file type > navigate to a storage location > touch **Save**. You can export a report as a PDF (.pdf) or an Excel file (.xls).
 - **Print:** Opens a Save As window. Enter a filename > navigate to a storage location > touch **Save**. You can export a report as a .xps file for printing on a non-system PC connected to a printer.
- 3 Touch **Done**.

Customizing the Logo for a Report

See “Managing Report Settings” on page 156.

10 Maintaining the Simoa® HD-X Analyzer

Performing Maintenance Tasks	128
Understanding the Maintenance Screen	129
Determining Whether a Maintenance Task Was Performed.....	129
Setting and Viewing Task Details	130
Changing Columns.....	131
Performing the Start of Day Task.....	132
Performing the End of Day Task.....	133
Performing the Idle Fluid Prime Task.....	134
Performing the Database Clean Task.....	135
Cleaning Exterior Instrument Surfaces	136
Cleaning the Touchscreen.....	136
Cleaning the System Resources Drawers.....	136
Cleaning the Reagent and Sample Bays.....	137
Cleaning the Liquid Waste Containers	137
Cleaning the System Fluid Containers and Tube Rinsing	138
Maintaining Idle Instruments.....	143

Performing Maintenance Tasks

Perform instrument maintenance tasks according to the following schedule table. These tasks are listed in the Maintenance tab.



Note! Check the Quanterix® customer portal at <http://portal.quanterix.com/> for updates to monthly maintenance procedures.


Task	Default Interval	Description
Start of Day	After initialization	Prepares instrument systems to start a run. See “Performing the Start of Day Task” on page 132.
End of Day	Daily after the last run of the day	Cleans the system at the end of the day.
Idle Fluid Prime	After 240 idle minutes	Primes the system fluids and resets the idle time counter.
Monthly Fluid Prime	Monthly	Primes the system fluids three times the normal length, approximately 30 minutes.
Replenish Oil	When the sealing oil is empty	Primes sealing oil through the entire line.
Database Clean	Database size limit reached, at least weekly	Cleans up the database. See “Performing the Database Clean Task” on page 135.
Computer Memory Management	Monthly	Removes temporary files and performs disk utility that improves data access speeds.



Note! The quality of system fluids may degrade over time due to limited shelf life, cross-contamination, salt buildup, etc. You must perform the monthly maintenance procedures described in this chapter to maintain optimal performance of the instrument. As a check, you can run the Simoa Qualification Test (SQT) to measure the precision of the system upon completion of the maintenance tasks.

Understanding the Maintenance Screen

The Maintenance screen shows status information about each maintenance task.

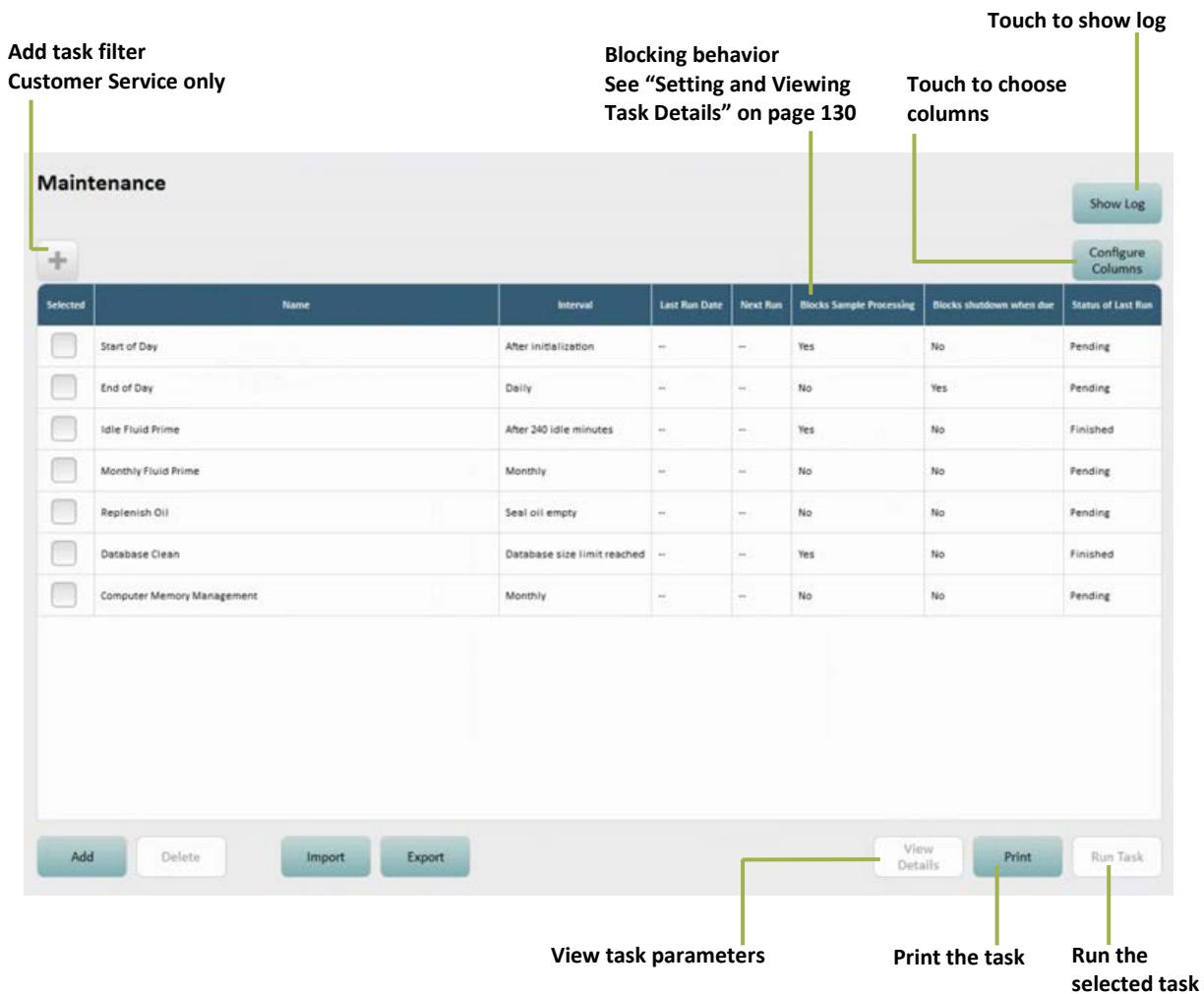

Determining Whether a Maintenance Task Was Performed

The Maintenance screen shows status information about each maintenance task. Two methods are available to do this:


- Check the Last Run column in the Maintenance tab, which reports the date and time when a task was last performed.
- In the Maintenance tab, touch **Show Log** to view a log that lists all completed maintenance tasks.

Setting and Viewing Task Details


Use the Task Details screen to set and display additional task settings.

❖ To view task details


Touch the checkbox for a task > touch **View Details**.

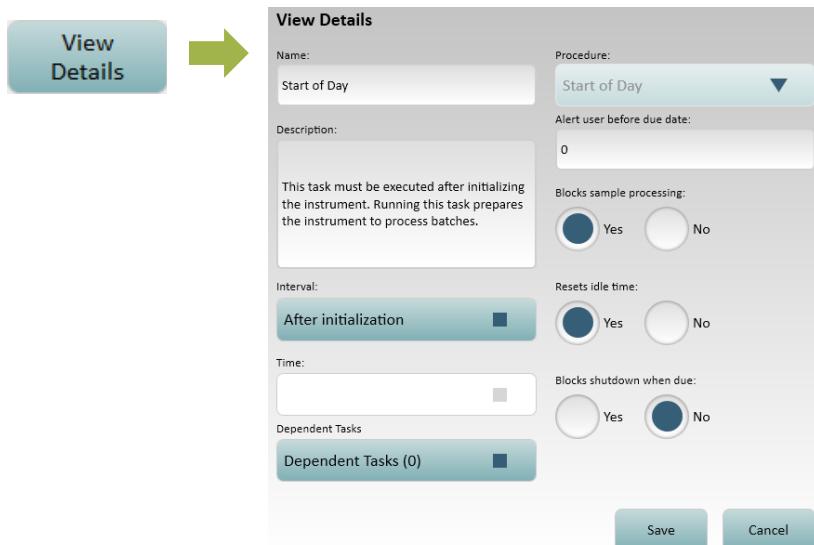

- **Interval** – The frequency of the task.
- **Time** – Time of day the task is performed.
- **Dependent Tasks** – Run specified tasks immediately after completion of the current task.

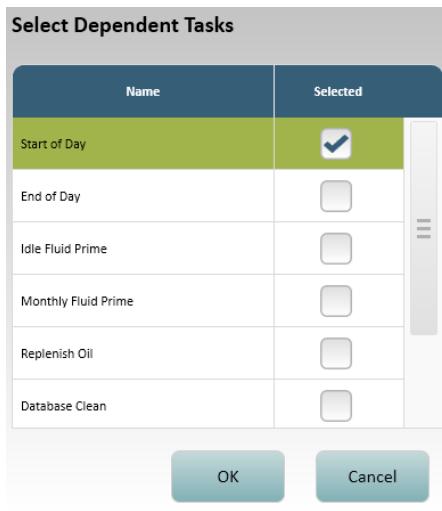


- **Alert user before due date** – Number of days prior to the due date that the software provides a notification prompting you to complete the task. This is only available for the Monthly Fluid Prime and Computer Memory Management tasks.
 - **Blocks sample processing** – Prevents/allows running assays when the task is not complete.
 - **Resets idle time** – Resets/does not reset counter after task is performed.
 - **Blocks shutdown when due** – When enabled, the following prompt appears when you shut down the software.

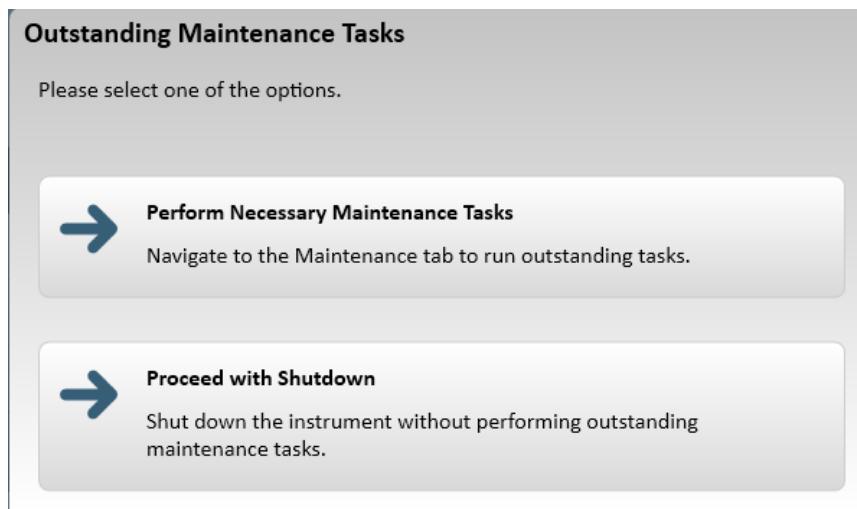

Changing Columns

To change the table columns, touch **Configure Columns** > touch the arrows to move column names between the Available and the Selected boxes.

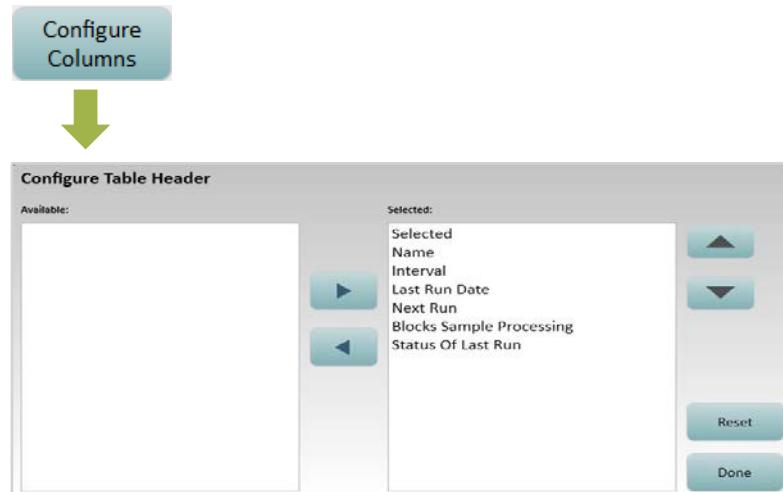

Performing the Start of Day Task


This task sets up the instrument for a new assay run as follows:

- Primes oil line.
- Flushes the pipettor with deionized water, cleans it with system wash buffer 1, and flushes it again with deionized water.
- Primes tubes with wash buffer fluid, then discards the cuvettes.
- Loads 10 cuvettes into the Wash/Incubator rings.


After you perform the End of Day Task, you must perform the Start of Day Task before performing a run.

❖ To perform the Start of Day task

- a Fill a clean reagent bottle a little more than half way with system wash buffer 1. Make sure there are no bubbles in the bottle.
- b Touch the Maintenance Tab > touch the Start of Day task checkbox > touch **Run Task**.

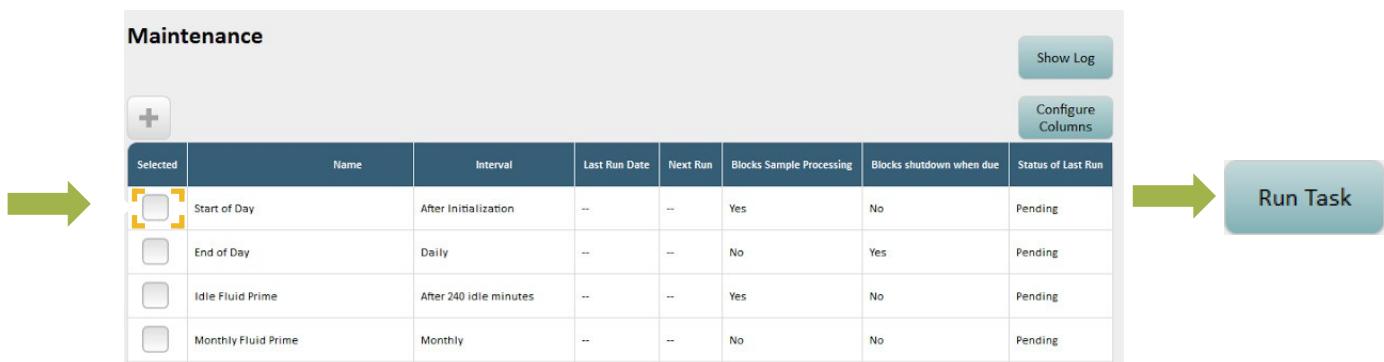

- c When prompted, place the reagent bottle into position 3 on a reagent rack (fourth position from the handle) and put the rack in the rightmost reagent bay position (lane 4).

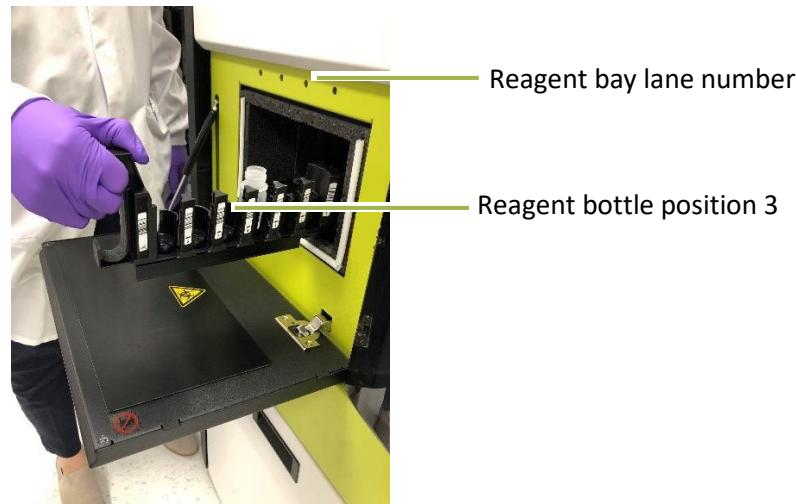

- d After loading the reagent rack into the instrument, select “next step” to continue with the Start of Day task.
- e Allow approximately 20 minutes for completion. When the task is complete, touch **Close**, remove the reagent rack and dispose of the reagent bottle.

Performing the End of Day Task

The End of Day task flushes the pipettor that aspirates reagents and all pipettor and system fluid tubings to prevent the buildup of salts, which can clog the pipettor and instrument tubings.

The task performs these actions:

- Flushes the pipettor with deionized water, cleans it with system wash buffer 1, and flushes it again with deionized water.
- Primes the tubing for system wash buffer 1 and 2 with deionized water.

After you perform the End of Day Task, you must perform the Start of Day Task before performing a run. Quanterix recommends setting a dependent task to automate this process, as described in Setting and Viewing Task Details.

❖ To perform the End of Day task

- 1 Fill a clean reagent bottle a little more than half way with system wash buffer 1. Make sure that there are no bubbles in the bottle.
- 2 Touch the Maintenance tab > touch the End of Day task checkbox > touch **Run Task**.

A screenshot of a mobile application interface showing a list of maintenance tasks. The tasks are listed in a table with columns: Selected, Name, Interval, Last Run Date, Next Run, Blocks Sample Processing, Blocks shutdown when due, and Status of Last Run. The 'End of Day' task is highlighted with a yellow border around its checkbox and name. A green arrow points to the 'Selected' column of the 'End of Day' row. Another green arrow points to a blue button labeled 'Run Task' located to the right of the table.

Selected	Name	Interval	Last Run Date	Next Run	Blocks Sample Processing	Blocks shutdown when due	Status of Last Run
<input type="checkbox"/>	Start of Day	After Initialization	--	--	Yes	No	Pending
<input checked="" type="checkbox"/>	End of Day	Daily	--	--	No	Yes	Pending
<input type="checkbox"/>	Idle Fluid Prime	After 240 idle minutes	--	--	Yes	No	Pending
<input type="checkbox"/>	Monthly Fluid Prime	Monthly	--	--	No	No	Pending
<input type="checkbox"/>	Replenish Oil	Seal Oil Empty	--	--	No	No	Pending
<input type="checkbox"/>	Database Clean	Database size limit reached	--	--	Yes	No	Finished
<input type="checkbox"/>	Temporary Files Clean	Monthly	--	--	No	No	Pending

- 3 When prompted, place the bottle into position 3 on a reagent rack (fourth position from the handle) and put the rack in the rightmost reagent bay position (lane 4). See the photo in the section above.
- 4 Allow approximately 20 minutes for completion. When the task is complete, touch **Close**, and remove the reagent rack and dispose of the reagent bottle.

Performing the Idle Fluid Prime Task

Prime the system fluids after 240 idle minutes. Idle time is the amount of time since the instrument finished a run, or the time since a task that resets the idle counter was run.

- 1 In the Maintenance tab, touch the Idle Fluid Prime checkbox > touch **Run Task**.
- 2 The system fluids are primed and the idle time counter is set to zero.

A screenshot of a mobile application interface showing a list of maintenance tasks. The tasks are listed in a table with columns: Selected, Name, Interval, Last Run Date, Next Run, Blocks Sample Processing, Blocks shutdown when due, and Status of Last Run. The 'Idle Fluid Prime' task is highlighted with a yellow border around its checkbox and name. A green arrow points to the 'Selected' column of the 'Idle Fluid Prime' row. Another green arrow points to a blue button labeled 'Run Task' located to the right of the table.

Selected	Name	Interval	Last Run Date	Next Run	Blocks Sample Processing	Blocks shutdown when due	Status of Last Run
<input type="checkbox"/>	Start of Day	After Initialization	--	--	Yes	No	Pending
<input type="checkbox"/>	End of Day	Daily	--	--	No	Yes	Pending
<input checked="" type="checkbox"/>	Idle Fluid Prime	After 240 idle minutes	--	--	Yes	No	Pending
<input type="checkbox"/>	Monthly Fluid Prime	Monthly	--	--	No	No	Pending
<input type="checkbox"/>	Replenish Oil	Seal Oil Empty	--	--	No	No	Pending
<input type="checkbox"/>	Database Clean	Database size limit reached	--	--	Yes	No	Finished
<input type="checkbox"/>	Temporary Files Clean	Monthly	--	--	No	No	Pending

Performing the Database Clean Task

The Simoa HD-X Analyzer generates a lot of data for each sample—approximately 115 MB per result—mostly consisting of raw camera images. For example, a batch with 96 samples results in approximately 11 GB of data. This data is stored in a database on the instrument computer.

Each week, you must perform the Database Clean task to purge old data so that new data can be stored. The task removes all data for runs older than 30 days from the first time that you perform the task.

If you want to retain data older than 30 days or to protect data against potential disk failure, it is *your responsibility* to back up the database prior to running the Database Clean task (at least monthly).

The Database Clean task permanently deletes data, so it is highly recommended that you back up the database before performing this task, especially if retention of raw sample images is required.

Backup Strategy

Contact Customer Support to discuss best practices regarding database backup and to obtain sample scripts that automate the backup process.

Do not store database backups on the instrument computer, for two reasons:

- Data can be lost if the disk fails.
- Database backups are very large and consume disk space that is needed to operate the instrument.

If the backup strategy you select results in backup files being stored on the instrument hard drive, you must move them off the disk to an external hard drive, network drive, or similar destination.

You can perform a partial backup that contains all of the data in the database except for the raw IPL image files or a complete database backup that includes the raw IPL image files.

❖ To perform a database backup

See the *Simoa HD-X Analyzer Quanterix Customer Support Tool 2.0 User Guide* for information on performing a full backup, or follow the instructions in *Simoa HD-X Analyzer IT Requirements Guide* to make a partial backup.

- ❖ To perform the Database Clean task
-



Note! Be sure to back up the database before you perform this task to avoid permanently losing data.

In the Maintenance tab, touch the Database Clean checkbox > touch **Run Task**. After the Database Clean task completes, exit the HD-X software and reboot the PC.

Cleaning Exterior Instrument Surfaces

Clean exterior instrument surfaces whenever necessary (at least once per month) by wiping with a soft, clean cloth moistened with water or 10% ethanol.

Cleaning the Touchscreen

Use any glass cleaner that does not contain ammonia.

Wet a paper towel with glass cleaner and then gently wipe down the touchscreen. Do not spray glass cleaner directly onto the touchscreen.

Cleaning the System Resources Drawers

Once a month, clean the system resource drawers.



CAUTION! Do not attempt to clean the system resource drawers when an assay is running. If you open a drawer during a run, the instrument halts the run and you lose any accumulated data.

Completely pull out both drawers from the instrument and wipe the metal frame beneath the drawers with water and then with 70% ethanol.

Cleaning the Reagent and Sample Bays

Once a month, clean the inside surfaces of the reagent and sample bays.



WARNING! Turn off the instrument before you perform any interior surface cleaning. The robotic pipettors inside the instrument can injure your hands. Do not attempt to clean the incubator or the washer rings and do not insert your hands into the sample bay or the reagent bay if the instrument is turned on.

❖ To clean the reagent and sample bays

- 1 Turn off the instrument (see “Turning the Instrument On and Off” on page 25).
- 2 Remove all reagent racks, plate racks, tube racks, and RGP racks.
- 3 Wipe the walls, ceiling, and base of each bay first with water and then with 70% ethanol.
- 4 Replace the racks.

Cleaning the Liquid Waste Containers

Once a month, clean liquid waste containers and the platform on which they sit.

During this procedure, discard all fluids and rinsates according to your company’s waste disposal policies.

❖ To clean the liquid waste containers

- 1 Power off the instrument and remove both liquid waste containers by disconnecting the container sensor. If the instrument is not powered off, it will halt when both liquid waste containers are removed.
- 2 Wipe the platform on which the containers sit first with water and then with 70% ethanol.
- 3 Rinse the containers by performing the following procedure three times on each container:
 - a Fill the container at least one-third full of deionized water.
 - b Shake to rinse inner surfaces.
 - c Discard the rinsate.
- 4 Replace the containers and reconnect the container sensors. Align the tab in the sensor with the cutout in the container. When done correctly, you will hear a click when you reinsert each sensor.

Cleaning the System Fluid Containers and Tube Rinsing

Once a month, the fluid resource containers must be cleaned. The tubing from the containers into the system must be primed with water for 30 minutes before they are primed with buffer for another 30 minutes.

During this procedure, discard all fluids and rinsates according to your company's waste disposal policies.

The DI water and wash buffer 1 each have two containers, which are defined as follows:

- **Secondary Container** – This is the removable container with a screw cap where new fluid is added. The wash buffer 1 and DI water secondary containers are identical in size.
- **Primary Container** – This is the pyramidal-shaped container behind the secondary containers. The liquid level sensor and tubing are connected through a screw cap on top.

Wash buffer 2 only has one container, which sits between the DI water and wash buffer 1 secondary containers. It has a screw cap for adding new fluid, a liquid level sensor, and a fluid port that connects directly to the instrument tubing.

The tubing leading to the sensor cap is color coded as follows:

- **Blue** – DI water
- **White** – Wash buffer 1
- **Yellow** – Wash buffer 2


Note! Do not attempt to perform this procedure when an assay is running. The procedure involves disconnecting the fluid containers, which will prevent the instrument from obtaining the system fluids that are required to perform the assay.

❖ To clean the secondary and primary containers

- 1 Close the Simoa HD-X Analyzer software and shut down the instrument.
- 2 Remove the secondary containers and wash buffer 2 container.

Note! Do not reuse the fluids after emptying the containers. Reuse of fluids has been associated with a reduction in data quality.

- **DI water and wash buffer 1:** Lift the container up and then pull out.

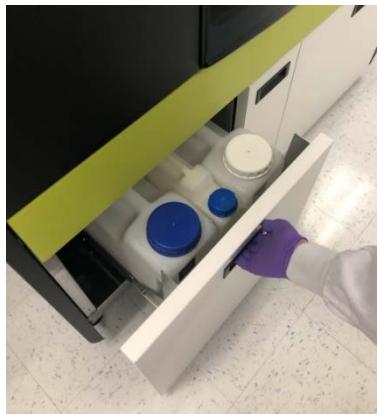

- **Wash buffer 2:** Disconnect the sensor and fluid port and then lift the container up.




3 Empty and dispose of all secondary and wash buffer 2 container fluids.



Note! Cross-contamination of the wash buffers may result in poor assay performance. Take care during this procedure to prevent cross-contamination. Do not reuse the fluids after emptying the containers. Reuse of wash buffers can result in a reduction of data quality.

4 Remove all primary containers. The process is the same for each fluid type.

- Disconnect the sensor by pulling the black connector up.

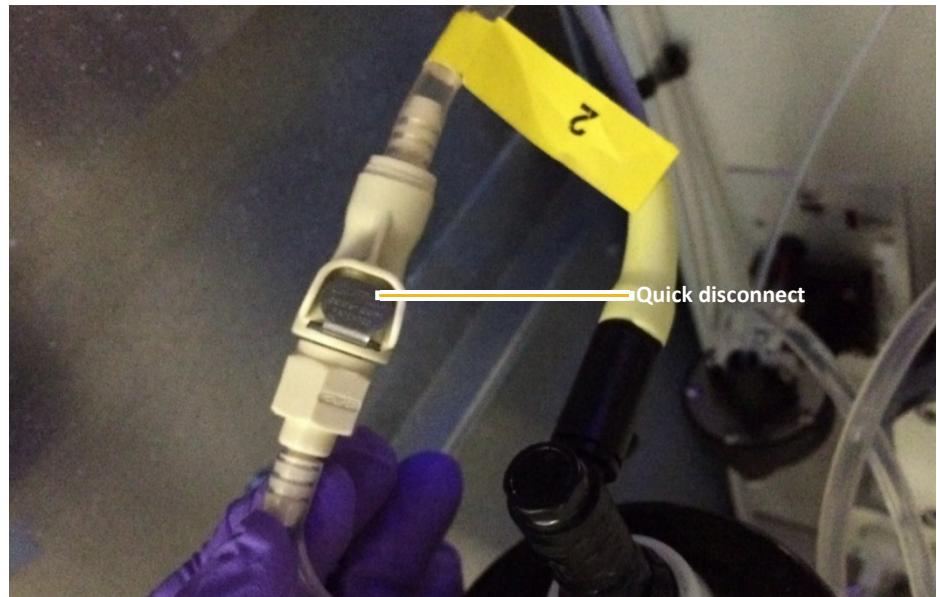

- Press the metal button on the quick disconnect fitting to detach the two ends of the liquid line.
- Remove the primary container by pulling forward and then up.



- Unscrew the black cap and place it on a clean lint-free paper towel.




- 5 Empty and dispose of all primary container fluids.
- 6 Wash all fluid containers (primary, secondary, and wash buffer 2) by triple-rinsing them with DI water. Fill the containers to approximately 1/3 of their total capacity with DI water and shake to rinse all surfaces.
- 7 Clean the liquid level sensors by running under DI water.

❖ To prime the system with DI water

- 1 Fill the containers with DI water as follows. Do not fill the secondary wash buffer 1 container.
 - **DI primary** – Full
 - **Wash buffer 1 primary** – Full
 - **Wash buffer 2** – Half full
 - **DI secondary** – Half full

- 2 Replace the sensor caps for the DI water and wash buffer 1 primary containers. If they have been separated from the containers, identify the caps by the physical differences in the tubing. DI water uses blue tubing, wash buffer 1 uses white and wash buffer 2 uses yellow.

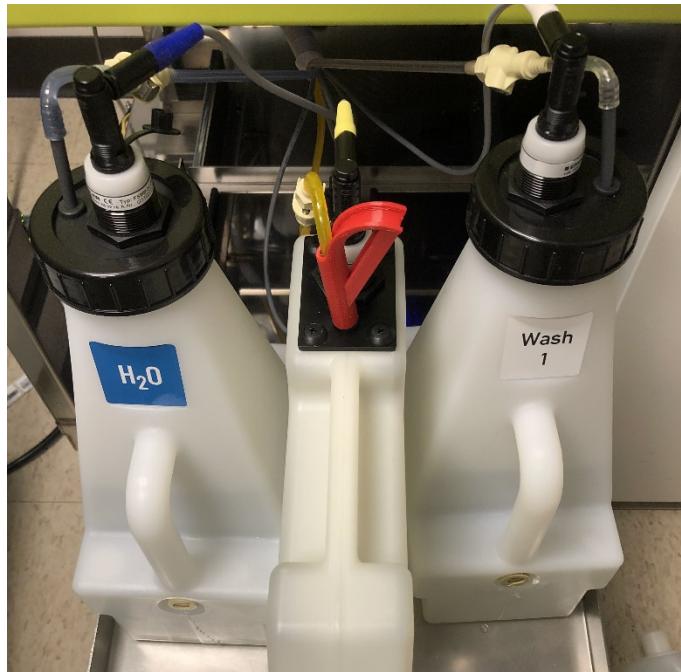

- 3 Replace the primary containers (DI water, wash buffer 1), wash buffer 2 container, and secondary DI container into the instrument. Do not replace the secondary container for wash buffer 1.
- 4 Reconnect the sensor cable and screw the fluid tubing together. Refer to the color code listed previously to ensure proper container connections.
 - Connect the two ends of the quick disconnect fitting. You will hear an audible click sound.
 - Connect the electrical connector into the top of the primary container lid, making sure that the tab on the cable is well-aligned with the slot in the lid. You will hear an audible click sound.
- 5 Turn on the instrument and open the Simoa HD-X software.
- 6 After initialization completes, run the Monthly Fluid Prime task. This step takes approximately 30 minutes to complete.
 - a Touch the Maintenance Tab > touch the Monthly Fluid Prime task checkbox.
 - b Touch **Run Task**. The task takes approximately 30 minutes to complete. The task window has a timer that indicates when the task will complete.
- 7 Close the Simoa HD-X software and shut down the instrument.

❖ To prime the system with wash buffer

- 1 Disconnect and remove the wash buffer 1 primary container and the wash buffer 2 container. Do not remove the DI Water primary container.
- 2 Unscrew the black cap for the wash buffer 1 primary container and place it and the liquid level sensors on a clean lint-free paper towel.
- 3 Empty the remaining water out of the containers.
- 4 Replace the wash buffer 1 primary container while it is still empty. Screw the sensor cap back on and connect the tubing and sensor again. Refer to the earlier instructions for detail.
- 5 Fill the containers with the following volumes:
 - **DI Water Secondary** – at least 5 L
 - **Wash buffer 1 Secondary** – at least 5 L
 - **Wash buffer 2** – at least 2 L
- 6 Replace the secondary containers and wash buffer 2 container into the instrument.
- 7 Turn on the instrument and open the Simoa HD-X software.
- 8 After initialization completes, run the Monthly Fluid Prime task again. This will take another 30 minutes to complete.

Maintaining Idle Instruments

When an instrument has been idle for a period of time, air can be introduced into the fluid lines, potentially impacting data quality and precision. To ensure proper instrument operation, perform maintenance tasks on idle systems as described in the following table.

Idle Period	Procedure
Less than 4 weeks	<p>Perform the following maintenance tasks twice per week for each week of inactivity.</p> <ol style="list-style-type: none"> 1 Power on computer and instrument. 2 Run Start of Day maintenance. 3 Run End of Day maintenance. 4 Power off computer and instrument. <p>If you cannot perform these tasks on schedule, run the Monthly Maintenance task (which performs an extended fluid prime) followed by a Simoa Qualification Test.</p>



Idle Period	Procedure
Greater than 4 weeks	Schedule a visit with Quanterix Service to perform maintenance before and after the scheduled idle period.

11 Configuring Software Settings

Software Settings	145
Changing System Settings	146
Selecting the Sample Loading Mode.....	146
Setting an Automatic Lock Out Time	147
Enabling Password Policies	147
Controlling Tab Switching when Bay Doors are Opened	148
Setting a Password Expiration Time.....	148
Setting Default Decimal Places	149
Specifying the Customer Defined System Name	149
Locking an Account due to Incorrect Login Attempts.....	150
Specifying User Settings	150
Specifying Account Management Settings.....	152
Managing Users.....	153
Managing Roles and Features	155

Software Settings

The Software Settings screen has the following tabs:

- System Settings
 - User Settings
 - Account Management (for users with administrator privileges only)
 - Reporting Settings
- ❖ [To view the software settings](#)

Touch the Settings icon at the top right of the software to open the Software Settings screen.



Changing System Settings

In the System Settings screen, you can change the following settings:

- Sample loading mode
- Automatic lock out time
- Password policies
- Synchronization of loading screen and bay doors
- Default number of decimal places
- Customer-defined system name


It is not necessary to restart the software after changes to the system settings.

Selecting the Sample Loading Mode

Plate mode is selected when you first open the software. To process samples in tubes, select Rack Mode.

❖ To change the sample loading mode

In the System Settings tab, touch **Plate Mode** or **Rack Mode** > touch **Done** in the upper right corner of the tab.

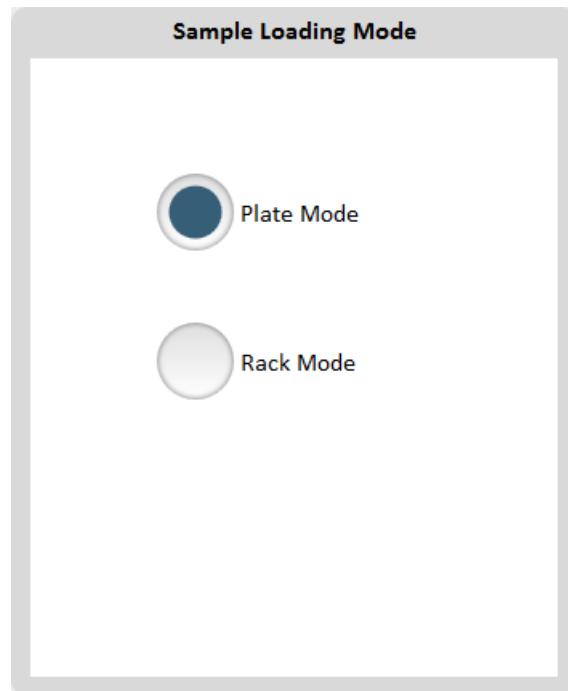


Setting an Automatic Lock Out Time

You can set the software to require users to log in again after a specified time period.

❖ To set a lock out time

In the System Settings tab, touch the **Enable Automatic Lock out** checkbox > type a Lock Out Time > touch **Done** in the upper right corner of the tab.

The Lock out Time value must be greater than 0 and less than or equal to 999.

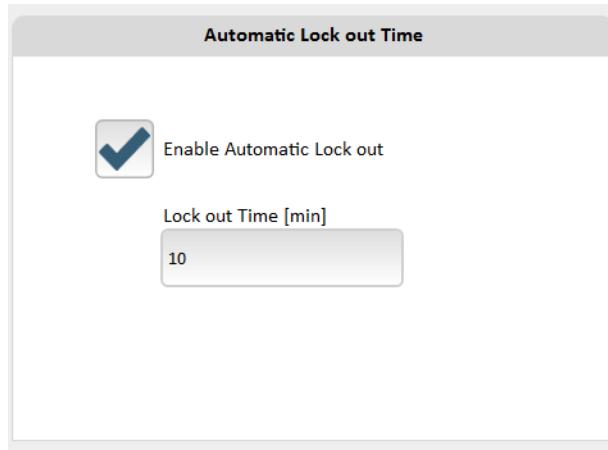


Enabling Password Policies

You can enable or disable the use of password policies. Enabling this setting places additional security requirements on new password values. The security requirements for new passwords are as follows:

- Usage of at least seven characters
- Usage of at least one upper case character
- Usage of at least one lower case character
- Usage of at least one non-alphanumeric character
- Usage of at least one numeric character

❖ To enable password policies

In the System Settings screen, scroll down to the General section > touch the **Enable Password Policies** checkbox > touch **Done** in the upper right corner of the tab.

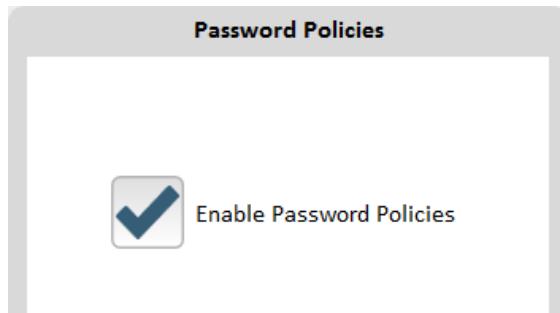

Controlling Tab Switching when Bay Doors are Opened

By default, the software switches to the Load Reagents tab when you open the reagent bay door, and it switches to the Setup Run tab when you open the sample bay door. You can turn on/off tab switching with the **Synchronize Loading Screen and Bay Doors** setting.

❖ To set tab switching behavior


- 1 In the System Settings screen, scroll down to the General section.
- 2 Enable the **Synchronize Loading Screen and Bay Doors** checkbox to switch tabs when you open a bay door, or disable it to keep the current tab open.
- 3 Touch **Done** in the upper right corner of the tab.

Setting a Password Expiration Time

You can enable password expiration and set the amount of time in which you want user passwords to expire.

❖ To set a password expiration time

In the Systems Settings tab, touch the **Enable Password Expiration** checkbox > type the number of days in which you want user passwords to expire > touch Done in the upper right corner of the tab.

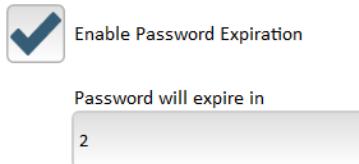


Setting Default Decimal Places

You can specify the number of decimal places that appear for all numbers (except concentrations, which are controlled at the Plex level). This setting controls the display of numbers throughout the application (in result tables, reports, etc.).

❖ To set the default digits to the right of the decimal place

In the System Settings screen, scroll down to the General section > enter an integer from 0 to 15 > touch **Done** in the upper right corner of the tab.

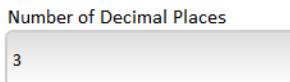

Specifying the Customer Defined System Name


To assist with troubleshooting, you can specify a name that describes the system configuration. For example:

- Version 1.2.3.4
- Simoa 1234
- Lab Machine 12

This name is captured in the Event Log screen and the shellDebug log file when the software is started.

❖ To set the customer defined system name

In the System Settings screen, scroll down to the General section > type a string that describes the system configuration (up to 20 characters) > touch **Done** in the upper right corner of the tab.

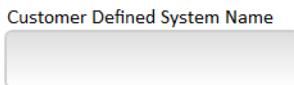


Locking an Account due to Incorrect Login Attempts

You can set the software to lock an account after a specified number of incorrect login attempts.

❖ To lock an account after a configurable number of incorrect login attempts

In the System Settings tab, scroll down to **General** > touch the **Lock Account after Configurable Number of Incorrect Login Attempts** checkbox > type the allowed number of incorrect login attempts > touch **Done** in the upper right corner of the tab.

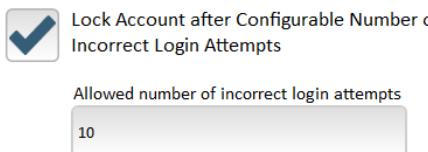

Specifying User Settings

All users can change their own system passwords. If you are a user with administrative privileges, you can:

- change another user's password
- show assays belonging to all users, even if the **Visible only to me** setting (in the Overview screen of the Custom Assay tab) is enabled for the account.

❖ To change a password

In the User Settings tab, type the user name, existing password, and new password > touch **OK**.



The dialog box is titled "Change Password". It contains four input fields: "User Name:" with the value "simoa", "Password:", "New Password:", and "Confirm Password:". At the bottom are two buttons: "OK" and "Cancel".

❖ To display/hide assays from other users




Note! You must be a user with administrative privileges to override the **Visible only to me setting** (in the Overview screen of the Custom Assay tab).

- 1 In the User Settings tab, enable **Show Assays from other Users** to display all assays, or disable it to show only your assays.
- 2 Touch **Done** in the upper right corner of the screen.



Specifying Account Management Settings

Use the Account Management screen to create and manage users and customize roles. The user that you define when you sign into the instrument software for the first time is assigned the AdminRole and appears as a user on this screen.

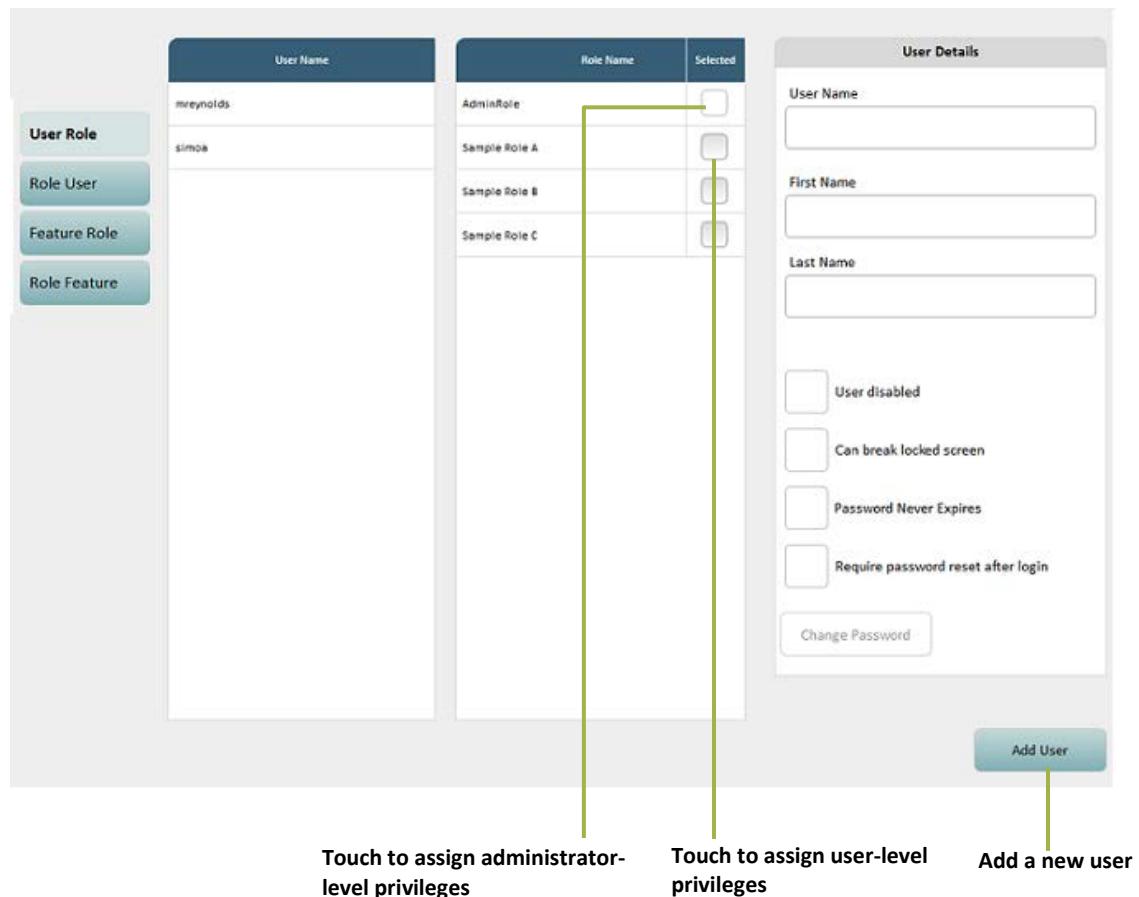

Managing Users

In the Account Management screen, you can add new users and define their access privileges.

❖ To add a new user

- 1 Touch **User Role**.
- 2 Touch **Add User**. Specify the User Name, First Name, Last Name, Password and Confirm Password.
- 3 To require the new user to reset their password after logging in, touch the **Require password reset after login** checkbox.
- 4 Touch **OK**.
- 5 Touch the **Selected** checkbox for the role you want to assign in the Role Name list.
- 6 If you have completed Account Management tasks, touch **Done**.



❖ To change a user's role

- 1 Touch **User Role**.
- 2 Select the user in the User Name list.
- 3 Touch the **Selected** checkbox for the role you want to assign in the Role Name list.
- 4 If you have completed Account Management tasks, touch **Done**.



❖ To change a user's password

- 1 Touch **User Role**.
- 2 Select the user in the User Name list.
- 3 Touch User Settings tab.
- 4 Touch **Change Password**.
- 5 Touch **Done** in the upper right corner of the screen.
- 6 Type the existing password and the new password.
- 7 Touch **OK**.
- 8 If you have completed Account Management tasks, touch **Done**.


❖ To assign multiple users to a single role

- 1 Touch **Role User**.
- 2 Select the role you are assigning in the Role Name list.
- 3 Touch the **Selected** checkbox for each user you want to assign.
- 4 If you have completed Account Management tasks, touch **Done**.

Managing Roles and Features

In the Account Management screen, you can create new roles and define the product features that each role can access. For example, you might want to limit the ability to modify imported assay definitions to administrators only. To do this, you can enable the Hidden option for the Read/Write button on the Custom Assay screen for nonadministrator roles.

The default role, AdminRole, can access all features. You cannot change the status of any feature to Disabled or Hidden for AdminRole.

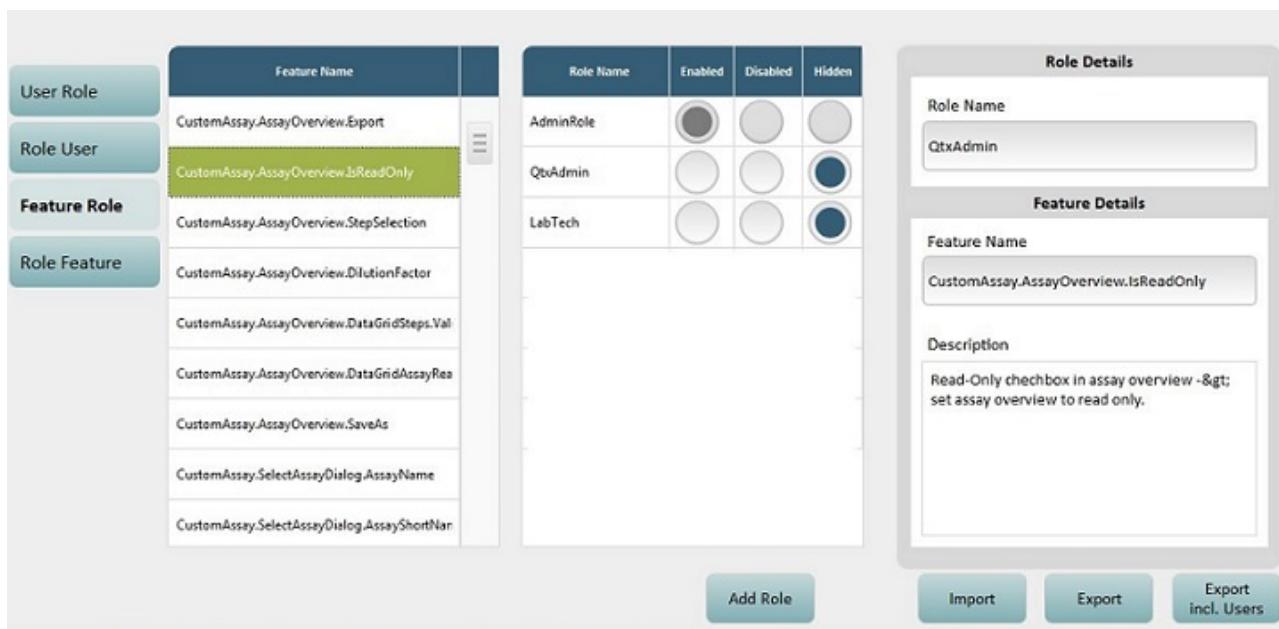

❖ To define features for a role

- 1 Touch **Feature Role**.
- 2 Select the feature you are assigning in the Feature Name list or select multiple features with the Shift key or Control key.
- 3 Touch **Enabled** (role can use the feature), **Disabled** (role cannot use the feature), or **Hidden** (feature does not appear in the software) for each role.
- 4 If you have completed Account Management tasks, touch **Done**.



Note! You can use the **Import**, **Export**, and **Export incl. Users** options to share settings among multiple instruments.

❖ To define multiple features for a role


- 1 Touch **Role Feature**.
- 2 Select the role name you are configuring in the Role Name list or select multiple roles with the Shift key or Control key.
- 3 Touch **Enabled** (role can use the feature), **Disabled** (role cannot use the feature), or **Hidden** (feature does not appear in the software) for each feature.
- 4 If you have completed Account Management tasks, touch **Done**.

❖ To add a new role

- 1 Touch **Feature Role**.
- 2 Touch **Add Role**.
- 3 Enter a **Role Name**.
- 4 Optionally, select a role that is similar to the one you are creating from the dropdown list.
- 5 Touch **OK**.

Managing Report Settings

In the Report Settings tab, you can customize the company name and logo (.jpg, .jpeg, .jpe, .png, or .bmp format) that appear on reports. You can also specify the default output format and path for automatic batch reports, run history table and the Clinical Results Report. There is the option to select All Columns and export the run history report in .csv format.

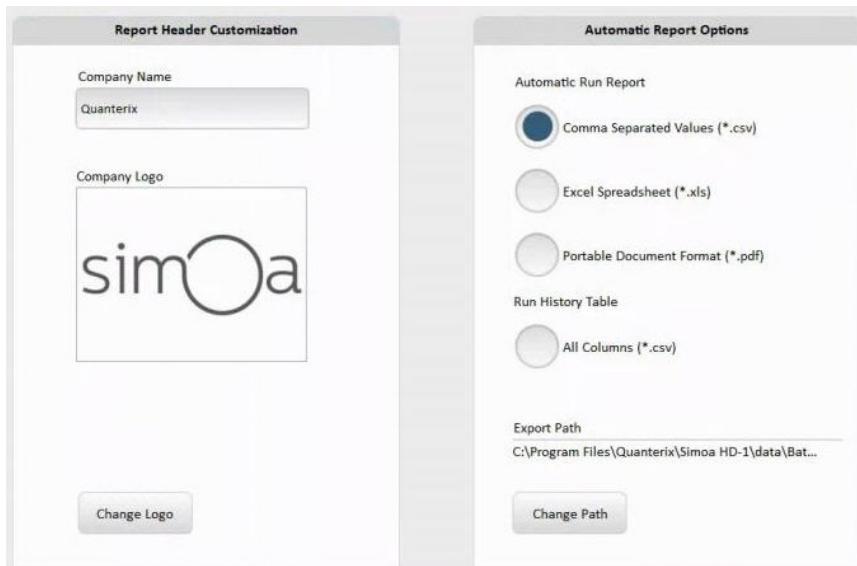

❖ To customize report headers

- 1 In the Report Settings screen, enter the Company Name.
- 2 Touch **Change Logo**.
- 3 Browse to the logo file > touch **Open**.
- 4 Touch **Done** in the upper right corner of the tab.

❖ To customize report export defaults

- 1 In the Report Settings screen, select the Format.
- 2 Touch **Change Path**.
- 3 Browse to the folder you want to use for report output, or touch **New Folder**, specify a new folder name for reports.
- 4 Touch **Select**.
- 5 Touch **Done** in the upper right corner of the tab.

Appendix A – Configure Columns Definitions

The following table describes each column in the Configure Table Header screens. The Configure Column function is accessed on the History & Reports tab for the Assay Reagent, Exceptions, Kit Search, Maintenance, Message Log, Sample Results, Run History and Event Log reports. The Configure Column function for the Maintenance report is located on the Maintenance tab.

Column Name	Comments
Affected Assay	Contains the name of the assay of the related job that is affected by the event (this does not relate to changes made to assays).
Affected Job	Displays which jobs are affected by reported messages/events.
Affected Plex	Contains the name of the plex of the related job or test order result that is affected by the event.
Affected Test Order	Contains the unique ID to identify the test order that is affected by the event.
Analog AEB	AEB value calculated using the analog mode. See Analysis Mode in this table to find out whether digital or analog mode was used to calculate the result AEB value. If the sample was run in replicate, the majority rule is used to determine the mode.
AEB	Average enzymes per bead.
Analysis Mode	Indicates whether the AEB calculation was performed using the digital or analog mode.
Assay	Name of the assay scheduled for the sample.
Assay Revision	The revision of the assay definition used.

Column Name	Comments
Batch ID	Unique system identifier of the sample batch.
Batch Name	Name of the batch specified on the Setup Run screen. If you do not specify a batch name, it defaults to batch creation date and time.
Blocks Sample Processing	Prevents user from starting a run when the task is not complete.
Blocks Shutdown when Due	Enable/disable. When enabled, a prompt appears when you shut down the software to either perform maintenance tasks or proceed with shutdown.
Calibration Curve ID	Unique system identifier of the calibration curve used in sample concentration calculations. This value is blank for calibrators as more than one calibration curve can be generated against calibrator AEB values.
Carrier Barcode	Barcode value of a microtiter plate or tube rack specified in the Setup Run screen.
Category	Categories for the given event. Possible values are: General, Maintenance, Job, L1DR, and L2DR.
Completion Date	Date and time when processing of the sample result completed.
Concentration	-
Creation Date	-
Curve Name	Name of the calibration curve used in sample concentration calculations. This value defaults to Measured. The calibration curve name can be edited in the Data Analysis tab.
CV AEB	Coefficient of variation of sample replicate AEB values for a given plex.

Column Name	Comments
CV Conc.	Coefficient of variation of sample replicate concentration values for a given plex.
Date Curve Created	Date and time when the calibration curve used in the sample concentration calculations was created. This value can be used to distinguish between curves if they all have the default name of Measured.
Dates Reagents Were Used	-
Description	Description of maintenance task performed.
Digital AEB	AEB value calculated using the digital mode. See Analysis Mode in this table to determine if digital or analog mode was used to calculate the result AEB value. If the sample was run in replicate, the majority rule is used to determine the mode.
Dilution Description	Name of the dilution applied to the sample test. The dilution description normally appears as a ratio of sample to diluent, for example, 4x, but it can carry any name, for example, One to Four, etc.
Dilution Factor	Multiplication factor applied to sample results when the concentration is reported. For example, for 4x dilution, the Dilution Factor will be 4.
Entry ID	Numeric value that shows the unique event ID from the database.
Errors	IDs of error flags associated with a sample test. To view additional information in the Event Log, click the error hyperlink.
Estimated Time to Result	Date and time when the system expects to deliver the result value.
Event ID	Represents a unique ID to identify a specific event.

Column Name	Comments
Event ID Hex	Represents a unique ID to identify a specific event - built of hex numbers and dots ".".
Event Text ID	Brief description of event.
Exception Description	Description of error message for a particular run.
Extended Properties	A collection of image analysis diagnostic information. For system diagnostic use only.
First Loading Time	-
Flags	The number of warning or error flags associated with processing of a sample test.
Fraction Monomeric Beads	Fraction of bead content that is monomeric – used to calculate actual bead concentration during bead aggregation assays.
Fraction On	Fraction of beaded wells with enzyme activity for a given plex. Both the digital and analog AEB calculation methods use this value.
Ibead	Average signal growth of beaded wells with enzyme activity. This value is used to calculate the result AEB value in analog mode. In digital mode, this value is not used. Isingle is calculated with Ibead values of jobs meeting the appropriate constraints.
Image Quality Score	Metric that helps the software determine whether a result should be included or excluded from analysis, as determined by image-processing algorithms. Currently, the value defaults to 1 for all sample results. For system diagnostic use only.
Instrument SN	Unique instrument identifier.

Column Name	Comments
Isingle	Estimated average signal growth of beaded wells that contain one enzyme. Certain conditions are required to calculate this value.
Interval	Identifies frequency of the task (e.g., daily, monthly, database size limit reached, etc.).
Job Completion Date	-
Job Creation Date	-
Job ID	Unique system identifier for a set of activities that lead to one result. For example, if a sample is scheduled to run an assay with one specimen in Neat and 4x dilutions in triplicate, the system will generate results for 6 unique jobs.
Job Start Cycle	The internal system clock used to coordinate sample processing activities. For system diagnostic use only.
Job Status	Status of each job.
Last Run Date	Date and time of last maintenance task run.
Localized/English Message	Detailed description of event.
Location	Location where the sample was loaded in the instrument.
Lot Numbers	Components/kits have a lot number that identifies the specific manufactured lot.
Lot Expiration Date	-
Mean AEB	Average of sample replicate AEB values for a given plex.

Column Name	Comments
Mean Conc.	Concentration calculated by averaging available individual replicate concentrations.
Message	Shows the full description of the error.
Name	Name of maintenance task performed (see Page 128).
Next Run	-
Number of Beads	Number of beads identified and used to produce the result for a given plex.
Plex	Name of the assay plex.
Reagent Barcode	Barcode located on reagent bottle identifies the reagent. Instrument uses barcoded information during processing.
Reagent Lots	-
Reagent Name	-
Remaining On-Board Time	Recommended time on-board until reagent expiration.
Replicate AEB	Sample AEB value for a given plex of a given job.
Replicate Conc.	Sample concentration value for a given plex of a given job.
Replicate Result ID	Unique system identifier for a replicate summary result. For example, if a sample was processed in triplicate, the replicate summary result will contain the mean, standard deviation, and coefficient of variation of the replicates (see Mean AEB, Mean Conc., SD AEB, SD Conc., CV AEB, CV Conc. in this table). For system diagnostic use only.
Result ID	Unique system identifier for a single result generated by the system. For system diagnostic use only.

Column Name	Comments
Result Status	Indicates whether the result was included or excluded from result calculations. The value can be one of: Automatically Included, Automatically Excluded, Manually Included, and Manually Excluded.
Sample Barcode	Sample identifier assigned by the user in the Sample Assignment screen. The value is entered manually or by scanning a barcode.
Sample Type	Specimen or calibrator sample type.
Selected	A column of check boxes for selecting records.
Severity	Defined as either a message, warning, problem, processing error, system halted or software failure.
SD AEB	Standard deviation of sample replicate AEB values for a given plex.
SD Conc.	Standard deviation of sample replicate concentration values for a given plex.
Source	Identifies location from which event originates such as PC, firmware, etc.
Status of Last Run	Status will be displayed as pending or finished. If status is pending, it may need to be executed prior to run.
SW Version	Software version. Provide to Quanterix Customer Support during troubleshooting.
Test Order ID	Unique system identifier for a grouping of sample jobs. For system diagnostic use only.
Type	An icon that specifies the severity of the event. The number in center of icon signifies how many messages of that type are listed in the log. See Severity.

Column Name	Comments
Unit	Unit of sample concentration.
Used Disc Lot Numbers	-
Used Reagents	A list of reagent barcodes used to produce the result.
User Name	The user who initiated the sample test that produced this result.

Appendix B – Instrument Specifications

The instrument must be installed by qualified service personnel only. Please do not move the instrument once it is installed.

Physical Specifications

Weight (without fluids)	294 kg (648 lb)
Dimensions	Width: 141.4 cm (55.66 in) Depth: 78.7 cm (31 in) Depth with monitor: 89.7 cm (35.3 in) Height: 161.2 cm (63.5 in)

Installation & Service Dimension Requirements

Doorway width for installation	80.0 cm (31.5 in)
Service clearance at instrument rear	60.0 cm (23.6 in) If clearance at the rear is not sufficient, enough front clearance must be available to move the instrument away from the wall for service.
Service clearance at instrument sides	60.0 cm (23.6 in)
Operator clearance at instrument front	100.0 cm (39.38 in)
Ventilation clearance at instrument rear	Minimum of 20.5 cm (8 in)
Earthquake safety	For sites that require the system to be secured from tipping due to an earthquake, adhesive straps must be affixed to the <u>rear panel</u> of the instrument only.

Power Requirements

Instrument

Instrument	Electrical receptacles (2 each): <ul style="list-style-type: none"> Europe: DIN49441 North America: 5 – 15R
Voltage	100 V – 240 V ± 10%
Frequency	50 – 60 Hz
Input current	7.5 A @100 VAC, 3.2 A @240 VAC
Fuse	Thermal overload protection

Computer and UPS

Computer	Input: 100 – 240 VAC ~50 – 60 Hz, 300W Ethernet (RJ45)
Recommended UPS	APC Smart-UPS SRT 2200VA 120V recommended

Environmental Operating Conditions

Operating temperature	18°C to 25°C (64.4°F to 77°F)
Operating humidity	30% to 80%, noncondensing
BTU output during operating conditions (non-standby)	3288 BTU/Hr

Noise Emission

Noise emission	72 dB(A), distance 1 m (39.4 in)
----------------	----------------------------------

Lasers

The Simoa HD-X Analyzer is a class 2 laser product. It incorporates the following laser sources.

Laser Barcode Scanner

Maximal output radiation	1.3 mW
Maximal pulse duration	110 µs
Wavelength	650 nm
Reading distance	3.0 – 30.0 cm

Laser Light Barrier

Maximal output radiation	1.0 mW
Wavelength	650 nm

Packaging

See delivery note enclosed with the instrument for the contents, weight, and dimensions.

Number of crates	3
Crate 1	
Contents	Instrument
Dimensions (WxDxH)	160.0 cm x 92.0 cm x 178.0 cm (63 in x 36.2 in x 70.1 in)
Weight	378 kg (833 lb)

Crate 2	
Contents	Instrument covers
Dimensions (WxDxH)	157.0 cm x 80.0 cm x 97.0 cm (61.8 in x 31.5 in x 38.2 in) (Box 3 on top of crate 2 bundled)
Weight	130 kg (287 lb)
Crate 3	
Contents	Accessories
Dimensions (WxDxH)	120.0 cm x 80.0 cm x 50.0 cm (47.2 in x 31.5 in x 19.7 in) (packaged in card board box and can be handled separately)
Weight	39 kg (86 lb)

Environmental Conditions

The following table shows the range of conditions needed to run the system safely.

Environmental conditions	The system is made for indoor use.
Temperature	Operating: 18°C to 25°C (64.4°F to 77°F) Storage: 1°C to 45°C (33.8°F to 113°F) Transport: -20°C to 60°C (-4°F to 140°F)
Humidity	Operating: 30% – 80% noncondensing Storage: 5% – 80% noncondensing Transport: 20% – 80% noncondensing
Pollution degree	Degree 2
Overvoltage category	Class 2

Limit class	Class A (For industrial use. Domestic use restricted.)
Sunlight	No direct sunlight. May mislead optical sensors and affect performance.
Altitude	Up to 2000 m (1.24 mi) above mean sea level.
Dust	No excessive dust

Instrument Installation Dimensions

The following diagram shows a top-down view of the instrument.

