

Before you get started:

- Don't forget to activate your kit
- Make sure you wash your hands!

Your kit includes:



Swab



Collection tube



Barcode card



Biohazard bag



Return shipping mailer

Checklist before returning your kit:

- ☐ Did you activate your kit?
- ☐ Did you collect a sample from both nostrils?

Visit color.com/sample-return to learn more.

This self-sample collection kit has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA for use by designated laboratories

This self-sample collection kit has been authorized only for the self-collection and maintenance of nasal swab specimens as an aid in detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens; and,

The emergency use of this self-sample collection kit in combination with the authorized test is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of medical devices during the COVID-19 outbreak under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

color

COVID-19 Self-Swab Collection Kit

At-home Instructions

color

For In Vitro Diagnostic Use, For Prescription Use Only, For Use Under Emergency Use Authorization Only

1

Activate your kit online at color.com/covid/activate.

You need your unique barcode, which you'll find on your collection tube.



Warning: Your sample **must** be dropped off on the same day you activate and collect it.

2

Before collecting your sample, wash your hands thoroughly for 20 seconds and dry them well.


3



Remove the cap of the collection tube.

Keep the cap and collection tube nearby.

4

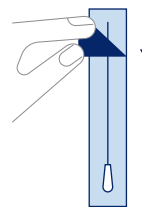


Open the swab package.

Starting from the handle end of the swab package, peel open the paper backing. Pull swab out of its packaging by the handle.

Note: Do not touch the soft tip with your hands or lay it down on any surface.

5

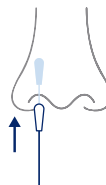

Rotate swab tip in first nostril, 4 times.

Insert the soft tip into one nostril, until it is no longer visible. Using medium pressure against the inside of your nostril the whole time, rub the swab in a large circular motion at least 4 times.

Note: Swab should not go further than ½ inch.

6



Rotate swab tip in other nostril, 4 times.

Repeat the previous step in the second nostril, using the same end of the swab.

Note: You will use the same swab for both nostrils.


7



Put swab into the collection tube.

The soft tip of the swab that went into your nose should go into the tube first.

8

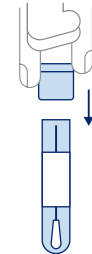

Press the cap of the collection tube firmly closed.

You're almost done! Make sure the top is closed tightly.

Note: After closing the collection tube, wash your hands thoroughly for 20 seconds and dry them well.

9



Put collection tube into the biohazard bag.

Close the bag with the adhesive seal.

10



Put biohazard bag into the return shipping mailer.

Remove the adhesive strip and seal the mailer closed. This packaging is specially designed for carrying samples suspected of carrying infectious substances. Do NOT replace it.

Note: Your return shipping has already been paid.

11

Return kit.

Your sample must be returned on the same day you collect it.

Tip: If you aren't able to return your sample on the same day, contact mycovidtest@color.com.

Questions?

Visit support.color.com or contact us at mycovidtest@color.com or (844) 352-6567.

Before you get started:

- Your sample must be dropped off on the same day you collect it.
- **Make sure there is a collection today.**
Refer to your program details if you're not sure where to return your sample.
- **Don't miss the last pickup.**
Refer to your program details for pickup times.
- Make sure you wash your hands!

Checklist before returning your kit:

- ☐ Did you activate your kit?
- ☐ Did you collect a sample from both nostrils?

Visit color.com/sample-return to learn more.

Questions?

Visit support.color.com or contact us at mycovidtest@color.com or (844) 352-6567.

This self-sample collection kit has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA for use by designated laboratories

This self-sample collection kit has been authorized only for the self-collection and maintenance of nasal swab specimens as an aid in detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens; and,

The emergency use of this self-sample collection kit in combination with the authorized test is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of medical devices during the COVID-19 outbreak under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

color

color

COVID-19 Self-Swab Collection Kit

On-site Instructions

For In Vitro Diagnostic Use, For Prescription Use Only, For Use Under Emergency Use Authorization Only


Your kit includes:




Swab



Collection tube



Barcode card

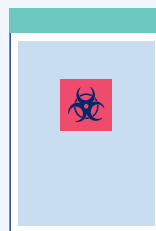

Biohazard bag

1

Activate your kit online at color.com/covid/activate.

You need your unique barcode, which you'll find on your collection tube.



Warning: Your sample will not be tested if the kit is not activated.


2

Before collecting your sample, wash your hands thoroughly for 20 seconds and dry them well.

3

Remove the cap of the collection tube.

Keep the cap and collection tube nearby.

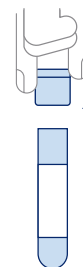


4

Open the swab package.

Starting from the handle end of the swab package, peel open the paper backing. Pull swab out of its packaging by the handle.

Note: Do not touch the soft tip with your hands or lay it down on any surface.

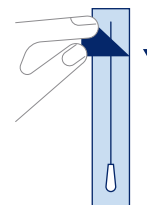


5

Rotate swab tip in first nostril, 4 times.

Insert the soft tip into one nostril, until it is no longer visible. Using medium pressure against the inside of your nostril the whole time, rub the swab in a large circular motion at least 4 times

Note: Swab should not go further than ½ inch.

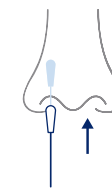

6

Rotate swab tip in other nostril, 4 times.

Repeat the previous step in the second nostril, using the same end of the swab.

Note: You will use the same swab for both nostrils.



7

Put swab into the collection tube.

The soft tip of the swab that went into your nose should go into the tube first.




8

Press the cap of the collection tube firmly closed.

You're almost done! Make sure the top is closed tightly.

Note: After closing the collection tube, wash your hands thoroughly for 20 seconds and dry them well.

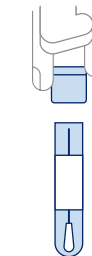

9

Put collection tube into the biohazard bag.

Seal the bag closed.



10

Bring the biohazard bag to your designated test collection point.

Refer to your program details if you're not sure where to return your sample.

Note: Don't miss the last pickup time.