

AMBRY COVID-19 RT-PCR TEST

Saliva Self-Collection & Shipping Instructions

Please follow the directions below to collect your COVID-19 testing sample, which was recommended by a physician for you to complete




IMPORTANT: Read all instructions thoroughly before you start to ensure proper sample collection for testing.

Do NOT eat, drink or chew gum for 30 minutes prior to collecting your sample.

Do NOT remove the plastic film from the funnel lid until you have collected the sample.


You can also watch a video on sample collection steps at: ambrygen.com/salivacollection

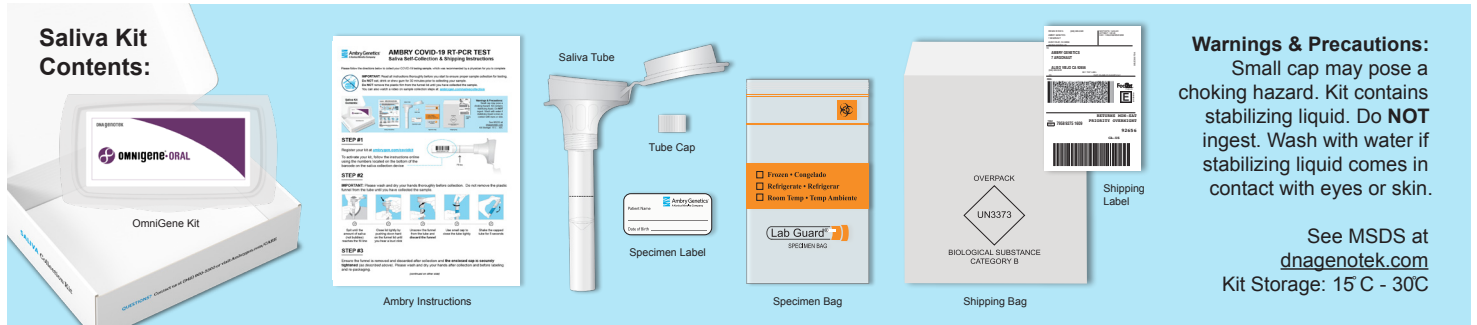

STEP #1


Register your kit at ambrygen.com/covidkit

To activate your kit, follow the instructions online using the numbers located on the bottom of the barcode on the saliva collection device

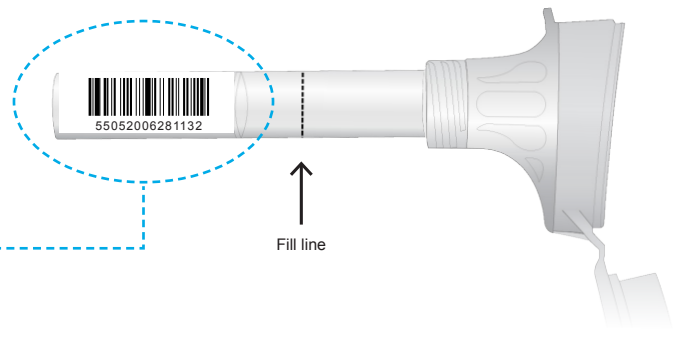


STEP #2


IMPORTANT: Please wash and dry your hands thoroughly before collection. Do not remove the plastic funnel from the tube until you have collected the sample.

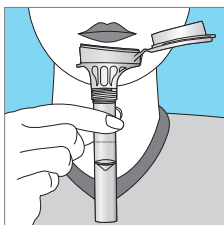


Spit until the amount of saliva (not bubbles) reaches the fill line

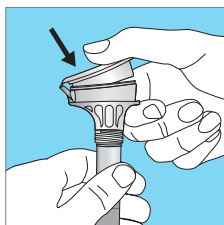


Close lid tightly by pushing down hard on the funnel lid until you hear a loud click

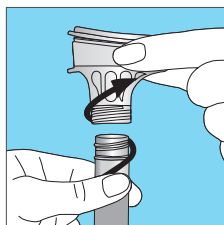

Unscrew the funnel from the tube and **discard the funnel**

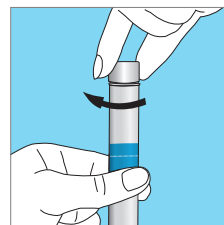

Use small cap to close the tube tightly

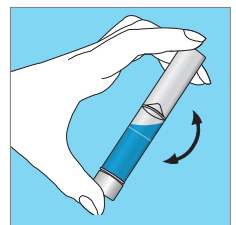

Shake the capped tube for 5 seconds


STEP #3

Ensure the funnel is removed and discarded after collection and **the enclosed cap is securely tightened (as described above)**. Please wash and dry your hands after collection and before labeling and re-packaging.

(continued on other side)

STEP #4


Write your first and last name, and date of birth
(MM/DD/YYYY) on the small Ambry sticker where indicated

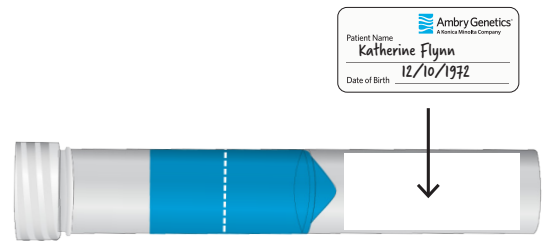

Affix the naming label on back of the tube as illustrated on the white rectangle on back of tube. Ensure the naming sticker does not cover the barcode on the tube.

NOTE: Testing may be delayed or not performed if the tube is received unlabeled, leaking, or damaged

STEP #5

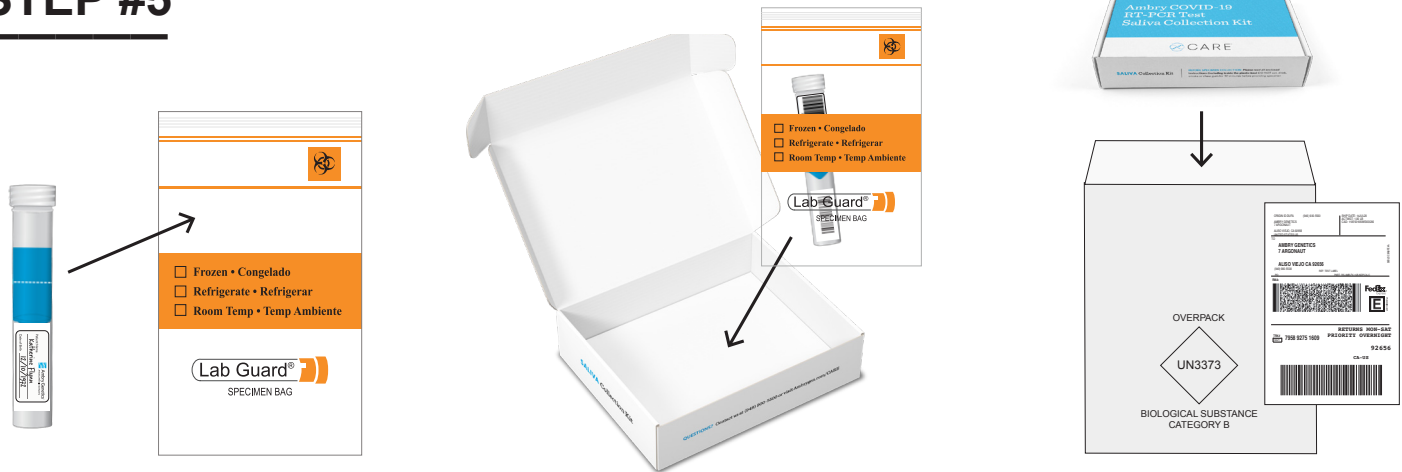

- ✓ Place the closed saliva tube in the biohazard bag; Do not remove the absorbent material from the biohazard bag.
- ✓ Seal the bag and place inside the cardboard box where the saliva collection kit was originally located.
- ✓ Place the cardboard box in the shipping bag that is already affixed with a FedEx return label; then seal the package.

STEP #6

Please drop off the package before the last pick up of the day, Monday - Friday at a FedEx dropbox. For locations, search by zip code at: fedex.com/en-us/shipping/dropbox.html

FOR EMERGENCY USE AUTHORIZATION ONLY
FOR PRESCRIPTION USE ONLY
FOR IN VITRO DIAGNOSTIC USE ONLY

FOR USE BY 18 YEARS OF AGE OR OLDER

This product has not been FDA cleared or approved but has been authorized for emergency use by FDA under an EUA for use by authorized laboratories. This product has been authorized only for the detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens.

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostic tests for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetics Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

For questions on how to send your sample to Ambry Genetics,
please contact our Customer Service team at care4covid@ambrygen.com